Mega» Mega Micromotor System

Easy to Use - Flexible - Safe - Powerful



Autoclavable Brushless Micromotor Guaranteed up to 1000 cycles



Electronic High Precision Calibration

Ensures safe treatment trought accurate torque correction








Implant Plus



Wide light area Light without shadows Power: 31000Lux

«DeLuxe»Multifunction Foot Control

Covered entirely with syntetic rubber: resistant to any leaking fluid. The reostatic effect is ensured by means of a specialized «pressure sensor». There are NO mechanical or moving parts.

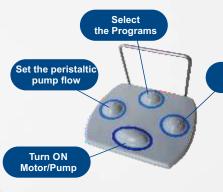

Select FWD/REV



Advanced Irrigation Pump

Easy set up of irragation tubes -Very ergonomic opening and closing system - Better control of irrigation flow with very constant flow rate - Stable and quiet operations - 3 Flow Rates available: MIN-MED-MAX



INCLUDED



I mplant «Plus»



Our powerful Implant Motor with integrated coolant pump is easy and simple to use. Reliable and convenient Autoclavable Brushless motor

Our **OMEGA IMPLANT plus** is a complete system, specifically designed dwith all the power and functionality needed for intelligent and precise implantology and surgery: **powerful performance and precise functionality.**See for yourself!

Complete, powerful and precise!

Our **OMEGA IMPLANT plus** is a compact tabletop unit for implantology that combines the power of a **BRUSHLESS motor** with a user friendly. All necessary treatment stages of implantology can be perormed safely and precisely using just one 20:1contra-angle.

OMEGA IMPLANT plus is equipped with **Contra-Angle calibration**: the precision of the motor is adapted to the contra-angle connected.

OMEGA IMPLANT plus also delivers the highest standards of ergonomics and flexibility. The motor torque is 80 Ncm and the motor speed range is 400 – 40,000 rpm (20 - 2000 with 20:1 contra-angle).

P «Programs»

Scrolls through the programs Choice of 5 programs settings Changed program automatically stored even after switching off -

Speed


Sets the rpm. Motor speeds from 20 to 2000 rpm with 20:1 contraangle

Torque

Changes the torque according to the reduction ratio

Flow

3 different flow rates are offered (Min-Med-Max (Max.= 100ml/min)

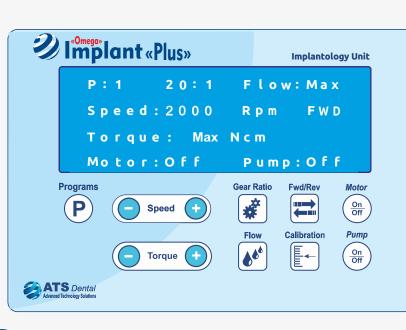

Gear Ratio

Sets the reduction of the contraangle Options: 1:1 – 16:1 – 20:1 – 24:1 – 32:1 – 64:1 – 80:1

Calibration

Compensates for loss of torque caused by the contra-angle

FWD/REV

Switches between forward and reverse drive Acoustic signal when reverse drive is activated

Motor on/off

Activate/Deactivate the Motor

Pump on/off

Activate/Deactivate the Pump

LED Technology



Button LED ON/OFF

LED Light

- wide light area
- Light without shadows
- -Power: 31000Lux

Model	Code	Description
OMEGA IMPLANT Plus	AOMG300DL	Implantology Unit: - Electronic Control Console - Autoclavable Brushless Micromotor - DeLuxe» Foot Control - Solution Holder - Motor Holder - 2 Tube



