

Ball Head Technology®

Pioneers in Angulation



2020 PRODUCT CATALOGUE
CONVENTIONAL PROSTHESIS



30° OF FREEDOM



The evolution towards customization in dental prosthetic design increasingly requires the use of angled abutments to provide versatility, simplicity and ease of use. Our Premium and Universal Solutions are pioneers in covering these needs.



Our main goal is achieving excellence to generate safety, reliability and comfort when using our products, both for the dental lab and clinical professionals, responsible for the patient's treatment.



Since our first patent in 2008 we have developed the most efficient and innovative angled solutions based on the vision and experience of the most qualified restoration and oral surgery professionals.

INDEX

1.0 CONVENTIONAL PROSTHESIS	5/88
1.1 DIRECT CASTING	5/35
1.2 OVERCASTING TO CrCo BASE	36/60
1.3 CEMENTING TO Ti BASE	61/88
 2.0 BHS30 SCREWDRIVER	 89
 3.0 PROTOCOLS	 90/97
 4.0 COMPATIBILITY CHART	 98

Av. de Les Garrigues 17
25001 Lleida · Spain
(+34) 973 831 391
info@ballheadtechnology.com
www.ballheadtechnology.com

All names and trademarks used in this catalogue have been registered and are property of their respective companies.
Images published are merely intended for commercial purposes and have no contractual nature.
The information contained in this catalogue is aimed exclusively to professionals in the dental sector.
Ball Head Technology® reserves the right to make changes to this document without prior notice.
The guarantee conditions of Ball Head Technology® will apply as long as the instructions for each reference are respected.



1.1 DIRECT CASTING



SIMPLICITY EQUALS EFFICIENCY

Engaging & non-engaging connections for all major implant systems.

Continuous surface for an easy wax-up & casting.

Optimized screw channel design and emergence profile.

10, 20 & 30° angle correction.

Soft tissue care.

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE								
I SERIES INTERNAL HEX								
	NARROW	Ø 3 - 3,75	30°	BHS07230	BHS07231	BHS02041	BHS02042	EM045130
			20°	BHS07232	BHS07233	BHS02021		
			10°	BHS07234	BHS07235			
	STANDARD	Ø 3,5 - 6,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
I 10C CONUS								
	RP	Ø 4,2	30°	BHS07067	BHS07068	BHS02023	BHS02024	EM048149
			20°	BHS07069	BHS07070			
			10°	BHS07071	BHS07072			
ACE COMPATIBLE								
INFINITY EXTERNAL HEX								
	NP	Ø 3,3 / 3,75	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013	BHS02001		
			10°	BHS07014	BHS07015			
	RP	Ø 4,0	30°	BHS07001	BHS07002	BHS02001	BHS02002	
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
INFINITY TRI-CAM								
	NP	Ø 3,5	30°	BHS07091	BHS07092	BHS02027	BHS02028	EM046291
			20°	BHS07093	BHS07094	BHS02009		
			10°	BHS07095	BHS07096			
	RP	Ø 4,3	30°	BHS07016	BHS07017	BHS02009	BHS02010	
			20°	BHS07031	BHS07032			
			10°	BHS07046	BHS07047			
	WP	Ø 5,0	30°	BHS07115	BHS07116	BHS02009	BHS02010	
			20°	BHS07117	BHS07118			
			10°	BHS07119	BHS07120			
INFINITY INTERNAL HEX								
	NP	Ø 3,7 / 4,1	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076	BHS02021		
			10°	BHS07077	BHS07078			
	RP	Ø 4,7 / 5,1	30°	BHS07024	BHS07025	BHS02021	BHS02022	
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			
ADIN COMPATIBLE								
SWELL								
		Ø 3,5 - 6	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			


1.1 DIRECT CASTING

	PLATFORM	Ø	△	CASTABLE	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOUAREG S/OS								
		Ø 3,5 - 6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
TOUAREG CLOSEFIT								
	RP	Ø 3,5	30° 20° 10°	BHS07079 BHS07081 BHS07083	BHS07080 BHS07082 BHS07084	BHS02025	BHS02026	EM048148
	WP	Ø 4,3 / 5	30° 20° 10°	BHS07067 BHS07069 BHS07071	BHS07068 BHS07070 BHS07072	BHS02023	BHS02024	EM048149
ALPHABIO COMPATIBLE								
EXTERNAL ALPHABIO								
	NP	Ø 3,5	30° 20° 10°	BHS07010 BHS07012 BHS07014	BHS07011 BHS07013 BHS07015	BHS02005	BHS02006	EM045222
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,1	30° 20° 10°	BHS07218 BHS07220 BHS07222	BHS07219 BHS07221 BHS07223	BHS02039	BHS02040	EM045070
INTERNAL ALPHABIO								
		Ø 3,3 - 6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
ALPHA DENT COMPATIBLE								
ACTIVE								
		Ø 3,75 / 4,2 / 5,0	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
ACTIVE PLUS								
		Ø 4,2 / 5,0	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
CLASSIC								
		Ø 3,3 / 3,75 / 4,2 / 5,0	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
ANCLADEN COMPATIBLE									
ANCLALOCK									
		Ø 3,75 - 5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130	
ANTHOGRYR COMPATIBLE									
ANTHOFIT HE									
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098	
			20°	BHS07004	BHS07005				
			10°	BHS07007	BHS07008				
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082	
			20°	BHS07063	BHS07064				
			10°	BHS07065	BHS07066				
OSSFIT									
	RP	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144	
			20°	BHS07035	BHS07036				
			10°	BHS07050	BHS07051				
	SYN_RP	Ø 4,8	30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375	
			20°	BHS07037	BHS07038				
			10°	BHS07052	BHS07053				
ARDS COMPATIBLE									
SMART									
		Ø 3,75 - 4,5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130	
	CLASSIC								
			Ø 3,3 - 6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
PREMIUM									
			Ø 3,3 - 6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130


1.1 DIRECT CASTING

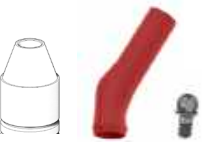	PLATFORM	Ø	△	CASTABLE	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ASTRA COMPATIBLE								
AQUA								
	RP_3.5 / 4.0		30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371
			20°	BHS07099	BHS07100			
			10°	BHS07101	BHS07102			
LILAC								
	WP_4.5 / 5.0		30°	BHS07130	BHS07131	BHS02009	BHS02010	EM046370
			20°	BHS07132	BHS07133			
			10°	BHS07134	BHS07135			
CONO 20°								
	Ø 3,5		30°		BHS07103	BHS02003	BHS02004	EM045638
			20°		BHS07104			
			10°		BHS07105			
CONO 45°								
	Ø 3,5		30°		BHS07106	BHS02003	BHS02004	EM047786
			20°		BHS07107			
			10°		BHS07108			
EV								
	Ø 3,6		30°	BHS07269	BHS07270	BHS02053	BHS02054	
			20°	BHS07271	BHS07272			
			10°	BHS07273	BHS07274			
	Ø 4,2		30°	BHS07263	BHS07264	BHS02055	BHS02056	
			20°	BHS07265	BHS07266			
			10°	BHS07267	BHS07268			
AVINENT COMPATIBLE								
EH								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02031	BHS02032	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
IH		Ø 3,5 / 4,1	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
MULTIUNIT		RP Ø 4,8	30° 20° 10°		BHS07003 BHS07006 BHS07009	BHS02003	BHS02004	EM045169
OCEAN CONICAL (C.C)		Ø 3,5 (Implant Ø 3,5/ Ø 4,0)	30° 20° 10°	BHS07079 BHS07081 BHS07083	BHS07080 BHS07082 BHS07084	BHS02025	BHS02026	EM048148
		Ø 4,1 (Implant Ø 4,5/ Ø 5,0)	30° 20° 10°	BHS07067 BHS07069 BHS07071	BHS07068 BHS07070 BHS07072	BHS02023	BHS02024	EM048149
BEGO COMPATIBLE								
EH		NP_Ø 2,7 / 3,1	30° 20° 10°	BHS07148 BHS07150 BHS07152	BHS07149 BHS07151 BHS07153	BHS02013	BHS02014	EM047957
IH		Ø 3,75 / Ø 4,1	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
BIOGENESIS COMPATIBLE								
3ICON		NP Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02001	BHS02002	EM045468
		RP Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
		WP Ø 5,0	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082

1.1 DIRECT CASTING

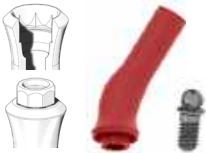	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ITICON								
	RP	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM045144
	RP	Ø 4,8	30° 20° 10°	BHS07022 BHS07037 BHS07052	BHS07023 BHS07038 BHS07053	BHS02013	BHS02014	EM046375
ATICON								
	RP	Ø 3,5 / 4,0	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371
	WP	Ø 4.5 / 5.0	30° 20° 10°	BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135	BHS02009	BHS02010	EM046370
KISSES								
	RP	Ø 3,5 / 4,0	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371
	WP	Ø 4.5 / 5.0 / 5,5	30° 20° 10°	BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135	BHS02009	BHS02010	EM046370
CONO 20°								
	RP / WP		30° 20° 10°		BHS07103 BHS07104 BHS07105	BHS02003	BHS02004	EM045638
CONO 45°								
	RP / WP		30° 20° 10°		BHS07106 BHS07107 BHS07108	BHS02003	BHS02004	EM047786
BIOHORIZONS COMPATIBLE								
EXTERNAL								
	NP	Ø 3,7	30° 20° 10°	BHS07236 BHS07238 BHS07240	BHS07237 BHS07239 BHS07241	BHS02031	BHS02032	
INTERNAL								
	NP	Ø 3,5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
	RP	Ø 4,5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218
	WP	Ø 5,7	30° 20° 10°	BHS07166 BHS07168 BHS07170	BHS07167 BHS07169 BHS07171	BHS02021	BHS02022	

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BIOMET 3i COMPATIBLE								
EH OSSEOTITE								
	NP	Ø 3,5	30°	BHS07194	BHS07195	BHS02001	BHS02002	EM045468
			20°	BHS07196	BHS07197			
			10°	BHS07198	BHS07199			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
CERTAIN								
	NP	Ø 3,4	30°	BHS07029	BHS07030	BHS02007	BHS02008	EM045166
			20°	BHS07044	BHS07045			
			10°	BHS07059	BHS07060			
	RP	Ø 4,1	30°	BHS07018	BHS07019	BHS02007	BHS02008	EM045165
			20°	BHS07033	BHS07034			
			10°	BHS07048	BHS07049			
	WP	Ø 5,0	30°	BHS07121	BHS07122	BHS02007	BHS02008	EM045161
			20°	BHS07123	BHS07124			
			10°	BHS07125	BHS07126			
BIONER COMPATIBLE								
IKELT / BIKELT								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
HIKELT								
	NP	Ø 3,8	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	RP	Ø 4,7	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			
TOP DM								
		Ø 3,5 / 4,0 / 5,0	30°	BHS07079	BHS07080	BHS02025	BHS02026	EM048148
			20°	BHS07081	BHS07082			
			10°	BHS07083	BHS07084			
MULTIUNIT								
	RP	Ø 4,8	30°		BHS07003	BHS02003	BHS02004	EM045169
			20°		BHS07006			
			10°		BHS07009			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BTI COMPATIBLE								
EH								
	NP Tiny	Ø 3,5	30° 20° 10°	BHS07148 BHS07150 BHS07152	BHS07149 BHS07151 BHS07153	BHS02013	BHS02014	EM047957
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,5	30° 20° 10°	BHS07154 BHS07156 BHS07158	BHS07155 BHS07157 BHS07159	BHS02013	BHS02014	EM047967
CUATRILOBULAR								
	RP	Ø 4,1	30° 20° 10°	BHS07160 BHS07162 BHS07164	BHS07161 BHS07163 BHS07165	BHS02021	BHS02022	EM047743
BTK COMPATIBLE								
8NECK INT								
	SR	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM045144
EVO DL EXT								
	NP	Ø 3,4/ Ø3,7	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02035	BHS02036	
	RP	Ø 4,1/ Ø4,6	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,0/ Ø5,6	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082
EVO DL INT								
		Ø3,5 / Ø3,7 / Ø4,6 / Ø5,0 / Ø5,6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
KLASSIC EXT								
	NP	Ø 3,4	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02035	BHS02036	
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,0	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLASSIC INT								
		Ø3,5 / Ø4,5 / Ø5,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
KONIC EXT								
	NP	Ø 3,4	30°	BHS07194	BHS07195	BHS02035	BHS02036	
			20°	BHS07196	BHS07197			
			10°	BHS07198	BHS07199			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
KONIC INT								
		Ø3,5 / Ø5,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
SAFE TL INT								
	SR	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
B&W COMPATIBLE								
HEXAGONO EXT								
	NP	Ø 3,3	30°	BHS07136	BHS07137	BHS02035	BHS02036	EM048022
			20°	BHS07138	BHS07139			
			10°	BHS07140	BHS07141			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
HEXAGONO INT								
	NP	Ø 3,3	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	RP	Ø 4,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	WP	Ø 5,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA		
CAMLOG COMPATIBLE										
CAMLOG										
		Ø 3,3	30°	BHS07275	BHS07276	BHS02049	BHS02050			
			20°	BHS07277	BHS07278					
			10°	BHS07279	BHS07280					
			Ø 3,8	30°	BHS07281	BHS07282	BHS02049		BHS02050	
				20°	BHS07283	BHS07284				
				10°	BHS07285	BHS07286				
			Ø 4,3	30°	BHS07287	BHS07288	BHS02049		BHS02050	
				20°	BHS07289	BHS07290				
				10°	BHS07291	BHS07292				
			Ø 5,0	30°	BHS07293	BHS07294	BHS02051		BHS02052	
				20°	BHS07295	BHS07296				
				10°	BHS07297	BHS07298				
CONEXÃO SIST. PROTESE COMPATIBLE										
MORSE INDEXADO										
		Ø 3,5 / 3,75 / 4,0 / 5,0	30°	BHS07079	BHS07080	BHS02025	BHS02026	EM048148		
			20°	BHS07081	BHS07082					
			10°	BHS07083	BHS07084					
			HEXAGONO EXT							
		Ø 3,3	30°	BHS07136	BHS07137	BHS02035	BHS02036	EM048022		
			20°	BHS07138	BHS07139					
			10°	BHS07140	BHS07141					
		Ø 3,75 / 4,0	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098		
			20°	BHS07004	BHS07005					
			10°	BHS07007	BHS07008					
		Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082		
			20°	BHS07063	BHS07064					
			10°	BHS07065	BHS07066					
		HEXAGONO INT								
				Ø 3,75 / 4,0 / 5,0	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
					20°	BHS07039	BHS07040			
10°	BHS07054				BHS07055					

1.1 DIRECT CASTING

	PLATFORM	Ø	△	CASTABLE	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CORTEX COMPATIBLE								
DYNAMIX								
		Ø 3,3 / 3,8 / 4,2 / 5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
CLASSIX								
		Ø 3,3 / 3,8 / 4,2 / 5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
SATURN								
		Ø 3,8 / 4,2	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
COWELLMEDI COMPATIBLE								
EXTERNAL								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
INTERNAL								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
SUBMERGED (INT)								
		Ø 3,5 / 4 / 4,5 / 5 / 6	30°	BHS07142	BHS07143	BHS02009	BHS02010	
			20°	BHS07144	BHS07145			
			10°	BHS07146	BHS07147			

1.1 DIRECT CASTING

	PLATFORM	Ø	△	CASTABLE	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DENTALTECH COMPATIBLE								
IMPLAGIC GII								
	RP	Ø 4,5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218
IMPLASSIC HX2								
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
DENTIS COMPATIBLE								
S-CLEAN								
	RP	Ø 3,4	30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010	
I-CLEAN								
	RP	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM045144
DENTIUM COMPATIBLE								
SIMPLE LINE II								
	RN	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM045144
IMPLANTUM								
		Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010	
SUPER LINE								
		Ø 3,6 / 4 / 4,5 / 5 / 6	30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010	

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DENTSPLY SIRONA COMPATIBLE								
XIVES								
	D 3.4	30°		BHS07308	BHS07309	BHS02047	BHS02048	
			20°	BHS07310	BHS07311			
			10°	BHS07312	BHS07313			
	D 3.8	30°		BHS07314	BHS07315	BHS02047	BHS02048	
			20°	BHS07316	BHS07317			
			10°	BHS07318	BHS07319			
	D 4.5	30°		BHS07320	BHS07321	BHS02047	BHS02048	
			20°	BHS07322	BHS07323			
			10°	BHS07324	BHS07325			
DIO IMPLANTS COMPATIBLE								
EXTERNAL								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
UF II								
	Ø 3,8 / 4 / 4,5 / 5 / 5,5	30°		BHS07142	BHS07143	BHS02009	BHS02010	
			20°	BHS07144	BHS07145			
			10°	BHS07146	BHS07147			
EDISON MEDICAL COMPATIBLE								
LAMINA (IH)								
	Ø 3,3/3,75/4,2/5,0/6,0	30°		BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
FORTIS (IH)								
	Ø 3,3/3,75/4,2/5,0/6,0	30°		BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			

1.1 DIRECT CASTING

	PLATFORM	Ø	△	CASTABLE		○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ELITE MEDICA COMPATIBLE									
EXTERNAL CONNECT									
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222	
			20°	BHS07012	BHS07013				
			10°	BHS07014	BHS07015				
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098	
			20°	BHS07004	BHS07005				
			10°	BHS07007	BHS07008				
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070	
			20°	BHS07220	BHS07221				
			10°	BHS07222	BHS07223				
EUROTEKNIKA COMPATIBLE									
UNEVA									
		Ø 3,6 / 4,1 / 4,8	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098	
			20°	BHS07004	BHS07005				
			10°	BHS07007	BHS07008				
NATURACTICS									
		Ø 3,5 / 4 / 4,5 / 5	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371	
			20°	BHS07099	BHS07100				
			10°	BHS07101	BHS07102				
NATURALL+									
		Ø 3,5 / 4 / 4,5 / 5	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371	
			20°	BHS07099	BHS07100				
			10°	BHS07101	BHS07102				
NATEA									
	NP	Ø 3,6 / 4,1	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371	
			20°	BHS07099	BHS07100				
			10°	BHS07101	BHS07102				
	RP	Ø 4,8 / 6	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371	
			20°	BHS07099	BHS07100				
			10°	BHS07101	BHS07102				
AESTHETICA									
		Ø 3,6 / 4,1 / 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144	
			20°	BHS07035	BHS07036				
			10°	BHS07050	BHS07051				

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
GALIMPLANT COMPATIBLE								
IPXE - EXTERNA								
	RP	Ø 3,5 / 4 / 5	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
ICI - INTERNA								
	RP	Ø 3,2 / 3,5 / 4 / 5	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371
			20°	BHS07099	BHS07100			
			10°	BHS07101	BHS07102			
IPX - INTERNA								
	RP	Ø 3,5 / 4 / 4,5 / 5	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371
			20°	BHS07099	BHS07100			
			10°	BHS07101	BHS07102			
GMI COMPATIBLE								
FRONTIER (IH)								
	RP Standard	Ø 3,3	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	WP Wide	Ø 4,3	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			
FRONTIER multi-aesthetic								
	RP	Ø 4,8	30°		BHS07003	BHS02003	BHS02004	EM045169
			20°		BHS07006			
			10°		BHS07009			
	WP	Ø 6,0	30°		BHS07248	BHS02031	BHS02032	
			20°		BHS07249			
			10°		BHS07250			
PHOENIX (EH)								
	Narrow	Ø 3,3	30°	BHS07236	BHS07237	BHS02031	BHS02032	
			20°	BHS07238	BHS07239			
			10°	BHS07240	BHS07241			
	Standard	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	Wide (Hex. 2,7)	Ø 5,1	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHOENIX (EH)								
	Wide (Hex. 3,4)	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
	Extra-Wide	Ø 6,0	30°	BHS07242	BHS07243	BHS02039	BHS02040	
			20°	BHS07244	BHS07245			
			10°	BHS07246	BHS07247			
PHOENIX multi-aesthetic								
	RP	Ø 4,8	30°		BHS07003	BHS02003	BHS02004	EM045169
			20°		BHS07006			
			10°		BHS07009			
	WP	Ø 6,0	30°		BHS07248	BHS02031	BHS02032	
			20°		BHS07249			
			10°		BHS07250			
GT MEDICAL COMPATIBLE								
BEST FIT HXE								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
BEST FIT HXI								
	WP	Ø 4,5 / 5,0	30°	BHS07130	BHS07131	BHS02009	BHS02010	EM046370
			20°	BHS07132	BHS07133			
			10°	BHS07134	BHS07135			
BEST FIT 002								
	NP	Ø 3,5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
BEST FIT OCT								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
	SYN_RN	Ø 4,8	30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375
			20°	BHS07037	BHS07038			
			10°	BHS07052	BHS07053			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
HI TECH COMPATIBLE									
SELF THREAD									
	SLT	Ø 3,3 / 3,75	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130	
	GTI	Ø 4,2 / 5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218	
	LGI PLUS								
			Ø 3,5 / 4,3 / 5 / 6	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
		IMPLANT DIRECT COMPATIBLE							
		INTERACTIVE							
	NP	Ø 3,5	30° 20° 10°	BHS07079 BHS07081 BHS07083	BHS07080 BHS07082 BHS07084	BHS02025	BHS02026	EM048148	
	RP	Ø 4,3	30° 20° 10°	BHS07067 BHS07069 BHS07071	BHS07068 BHS07070 BHS07072	BHS02023	BHS02024	EM048149	
	LEGACY								
		NP	Ø 3,5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
		RP	Ø 4,5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218
		WP	Ø 5,7	30° 20° 10°	BHS07166 BHS07168 BHS07170	BHS07167 BHS07169 BHS07171	BHS02021	BHS02022	
REPLUS / REPLANT									
		NP	Ø 3,5	30° 20° 10°	BHS07091 BHS07093 BHS07095	BHS07092 BHS07094 BHS07096	BHS02027	BHS02028	EM046291
		RP	Ø 4,3	30° 20° 10°	BHS07016 BHS07031 BHS07046	BHS07017 BHS07032 BHS07047	BHS02009	BHS02010	EM046292
		WP	Ø 5,0	30° 20° 10°	BHS07115 BHS07117 BHS07119	BHS07116 BHS07118 BHS07120	BHS02009	BHS02010	EM046293

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SWISHPLUS / SWISHPLANT								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
	SYN_RN	Ø 4,8	30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375
			20°	BHS07037	BHS07038			
			10°	BHS07052	BHS07053			
JDENTALCARE COMPATIBLE								
JDEVOLUTION								
		Ø 3,7 / 4,3 / 5 / 6	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
KEYSTONE DENTAL COMPATIBLE								
RESTORE								
	NP	Ø 3,3	30°	BHS07136	BHS07137	BHS02035	BHS02036	EM048022
			20°	BHS07138	BHS07139			
			10°	BHS07140	BHS07141			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
OCT-XT								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
	SYN_RN	Ø 4,8	30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375
			20°	BHS07037	BHS07038			
			10°	BHS07052	BHS07053			
RENOVA								
	NP	Ø 3,75	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	RP	Ø 4,5 / 4,75	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLOCKNER COMPATIBLE								
SK2 (EH)								
	RP	Ø 4,2	30°	BHS07109	BHS07110	BHS02031	BHS02032	EM046747
			20°	BHS07111	BHS07112			
			10°	BHS07113	BHS07114			
KL (EH)								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
LEADER COMPATIBLE								
TIXOS INTERNAL HEX								
	NP	Ø 3,5 / 4,0	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
TIXOS EXTERNAL HEX								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
MEDENTIS COMPATIBLE								
ICX								
		Ø 3,45 / 3,7 - 4,1 / 4,8	30°	BHS07257	BHS07258	BHS02015	BHS02016	
			20°	BHS07259	BHS07260			
			10°	BHS07261	BHS07262			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
MEGAGEN COMPATIBLE									
ANYONE INTERNAL									
	NP		30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010		
	RP		30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010		
	WP		30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010		
	ANYONE EXTERNAL								
	NP		30° 20° 10°	BHS07010 BHS07012 BHS07014	BHS07011 BHS07013 BHS07015	BHS02005	BHS02006	EM045222	
	RP		30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098	
	WP		30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082	
	OCTA Abutment								
		SYN_RN	Ø4,8	30° 20° 10°	BHS07022 BHS07037 BHS07052	BHS07023 BHS07038 BHS07053	BHS02013	BHS02014	EM046375
MICRODENT COMPATIBLE									
EKTOS (IH)									
	NP	Ø 3,5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130	
	RP	Ø 4,5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218	
	WP	Ø 5,7	30° 20° 10°	BHS07166 BHS07168 BHS07170	BHS07167 BHS07169 BHS07171	BHS02021	BHS02022		
	MINI (EH)								
	NP	Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02031	BHS02032		
	SUV (EH)								
		SUV	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
		SUV	Ø 5,1	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SYSTEM (EH)								
	MST	Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02031	BHS02032	
	MS	Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02031	BHS02032	
	M / MV	Ø 4,2	30° 20° 10°	BHS07332 BHS07334 BHS07336	BHS07333 BHS07335 BHS07337	BHS02013	BHS02014	EM048014
	MK	Ø 5,6	30° 20° 10°	BHS07242 BHS07244 BHS07246	BHS07243 BHS07245 BHS07247	BHS02039	BHS02040	
UNIVERSAL (EH)								
	MST	Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02031	BHS02032	
	MS	Ø 3,5	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02031	BHS02032	
	SU	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	SU	Ø 5,1	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082
	MK	Ø 5,6	30° 20° 10°	BHS07242 BHS07244 BHS07246	BHS07243 BHS07245 BHS07247	BHS02039	BHS02040	
MIS COMPATIBLE								
LANCE								
	STD	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,0	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082
SEVEN / M4								
	NARROW	Ø 3,30	30° 20° 10°	BHS07230 BHS07232 BHS07234	BHS07231 BHS07233 BHS07235	BHS02041	BHS02042	
	STANDARD	Ø 3,75 / 4,2	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
	WIDE	Ø 5,0 / 6,0	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218
MULTIUNIT MIS								
	RP	Ø 4,8	30° 20° 10°		BHS07127 BHS07128 BHS07129	BHS02037	BHS02038	

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
MOZO GRAU COMPATIBLE								
INHEX								
	RP	Ø 3,75 / 4,25	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371
	WP	Ø 5,0	30° 20° 10°	BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135	BHS02009	BHS02010	EM046370
	CONO20		30° 20° 10°		BHS07103 BHS07104 BHS07105	BHS02003	BHS02004	EM045638
	CONO45		30° 20° 10°		BHS07106 BHS07107 BHS07108	BHS02003	BHS02004	EM047786
OSSEUS								
	NP	Ø 3,4	30° 20° 10°	BHS07194 BHS07196 BHS07198	BHS07195 BHS07197 BHS07199	BHS02001	BHS02002	EM045468
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,0	30° 20° 10°	BHS07218 BHS07220 BHS07222	BHS07219 BHS07221 BHS07223	BHS02039	BHS02040	EM045070
MPI COMPATIBLE								
EXTERNA EH								
	NP	Ø 3,5	30° 20° 10°	BHS07148 BHS07150 BHS07152	BHS07149 BHS07151 BHS07153	BHS02001	BHS02002	EM047957
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098
	WP	Ø 5,0	30° 20° 10°	BHS07061 BHS07063 BHS07065	BHS07062 BHS07064 BHS07066	BHS02001	BHS02002	EM045082
INTERNA CM								
	RP	Ø 4,0	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371
	WP	Ø 5,0	30° 20° 10°	BHS07130 BHS07132 BHS07134	BHS07131 BHS07133 BHS07135	BHS02009	BHS02010	EM046370

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NEOBIOTECH COMPATIBLE								
EB EXTERNAL SYSTEM								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
IS INTERNAL SYSTEM								
		Ø 3,5 / 4 / 4,5 / 5	30°	BHS07142	BHS07143	BHS02009	BHS02010	
			20°	BHS07144	BHS07145			
			10°	BHS07146	BHS07147			
NEODENT COMPATIBLE								
TITAMAX HE SMART								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
MINI PILAR CM								
	RP	Ø 4,8	30°		BHS07003	BHS02003	BHS02004	EM045169
			20°		BHS07006			
			10°		BHS07009			
NOBEL COMPATIBLE								
BRÄNEMARK								
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			


1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
REPLACE									
	NP	Ø 3,5	30° 20° 10°	BHS07091 BHS07093 BHS07095	BHS07092 BHS07094 BHS07096	BHS02027	BHS02028	EM046291	
	RP	Ø 4,3	30° 20° 10°	BHS07016 BHS07031 BHS07046	BHS07017 BHS07032 BHS07047	BHS02009	BHS02010	EM046292	
	WP	Ø 5,0	30° 20° 10°	BHS07115 BHS07117 BHS07119	BHS07116 BHS07118 BHS07120	BHS02009	BHS02010	EM046293	
	ACTIVE								
	NP	Ø 3,5	30° 20° 10°	BHS07079 BHS07081 BHS07083	BHS07080 BHS07082 BHS07084	BHS02025	BHS02026	EM048148	
	RP	Ø4,3 / Ø5,0	30° 20° 10°	BHS07067 BHS07069 BHS07071	BHS07068 BHS07070 BHS07072	BHS02023	BHS02024	EM048149	
	MULTIUNIT								
	RP	Ø 4,8	30° 20° 10°		BHS07003 BHS07006 BHS07009	BHS02003	BHS02004	EM045169	
	NT TRADING COMPATIBLE								
2-CONnect									
	M1,4	Ø 5,0	30° 20° 10°		BHS07305 BHS07306 BHS07307	BHS02063	BHS02064		
	M1,6	Ø 5,0	30° 20° 10°		BHS07305 BHS07306 BHS07307	BHS02065	BHS02066		
	M1,8	Ø 5,0	30° 20° 10°		BHS07305 BHS07306 BHS07307	BHS02067	BHS02068		
	OSSTEM IMPLANT COMPATIBLE								
	OSSTEM TS								
		Ø 4,0 / 4,5 / 5,0 / 6,0	30° 20° 10°	BHS07142 BHS07144 BHS07146	BHS07143 BHS07145 BHS07147	BHS02009	BHS02010		

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OSTEM US								
	MINI	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222
			20°	BHS07012	BHS07013			
			10°	BHS07014	BHS07015			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,1	30°	BHS07218	BHS07219	BHS02039	BHS02040	EM045070
			20°	BHS07220	BHS07221			
			10°	BHS07222	BHS07223			
OSTEOPLUS COMPATIBLE								
ARGENTIUM								
	RP	Ø 4,2	30°	BHS07109	BHS07110	BHS02031	BHS02032	EM046747
			20°	BHS07111	BHS07112			
			10°	BHS07113	BHS07114			
OSTEOPLUS - SHe								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
OSTEOPLUS - SHi								
	RP	Ø 3,5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	WP	Ø 4,5	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			
OXTEIN COMPATIBLE								
HEX EXTERNA								
	MINI	Ø 3,3	30°	BHS07136	BHS07137	BHS02035	BHS02036	EM048022
			20°	BHS07138	BHS07139			
			10°	BHS07140	BHS07141			
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
	WP	Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
HEX INTERNA								
	NP	Ø 3,5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	RP	Ø 4,5	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			

1.1 DIRECT CASTING

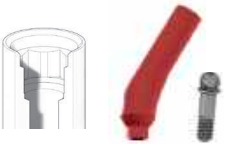	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
OCTA INTERNA									
	RN	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM045144	
	RN	Ø 4,8	30° 20° 10°	BHS07022 BHS07037 BHS07052	BHS07023 BHS07038 BHS07053	BHS02013	BHS02014	EM046375	
	PALTOP COMPATIBLE								
	INTERNAL HEX CONNECTION								
		NARROW	Ø 3,25	30° 20° 10°	BHS07230 BHS07232 BHS07234	BHS07231 BHS07233 BHS07235	BHS02041	BHS02042	
		STANDARD	Ø 3,75 / 4,2 / 5,0	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130
WIDE		6,0	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218	
PCA Implant System (Paltop Conical Active)									
NARROW		Ø 3,25	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371	
STANDARD		Ø 3,75 / 4,2 / 5,0	30° 20° 10°	BHS07097 BHS07099 BHS07101	BHS07098 BHS07100 BHS07102	BHS02025	BHS02026	EM046371	
MULTI UNIT									
			Ø 5,0	30° 20° 10°		BHS07326 BHS07327 BHS07328	BHS02031	BHS02032	
SINGLE UNIT									
			Ø 4,5	30° 20° 10°	BHS07329 BHS07330 BHS07331		BHS02031	BHS02032	
PHIBO COMPATIBLE									
TSH									
	S2	Ø 3,3	30° 20° 10°	BHS07136 BHS07138 BHS07140	BHS07137 BHS07139 BHS07141	BHS02035	BHS02036	EM048022	
	S3 / S4	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098	

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
AUREA EVO									
	NP	Ø 3,5	30°	BHS07178	BHS07179	BHS02041	BHS02042		
			20°	BHS07180	BHS07181				
			10°	BHS07182	BHS07183				
	RP	Ø 4,3	30°	BHS07186	BHS07187	BHS02021	BHS02022		
			20°	BHS07188	BHS07189				
			10°	BHS07190	BHS07191				
RADHEX IMPLANTS COMPATIBLE									
PHE / PHEA									
	NP	Ø 3,5	30°	BHS07010	BHS07011	BHS02005	BHS02006	EM045222	
			20°	BHS07012	BHS07013				
			10°	BHS07014	BHS07015				
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098	
			20°	BHS07004	BHS07005				
			10°	BHS07007	BHS07008				
PHI / PHIA									
	NP	Ø 3,5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130	
			20°	BHS07075	BHS07076				
			10°	BHS07077	BHS07078				
	RP	Ø 4,5	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218	
			20°	BHS07039	BHS07040				
			10°	BHS07054	BHS07055				
PCI									
		280	30°	BHS07097	BHS07098	BHS02025	BHS02026	EM046371	
			20°	BHS07099	BHS07100				
			10°	BHS07101	BHS07102				
		350	30°	BHS07130	BHS07131	BHS02009	BHS02010	EM046370	
			20°	BHS07132	BHS07133				
			10°	BHS07134	BHS07135				
SGS COMPATIBLE									
P1 / P7									
	RP	Ø 3,75	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130	
			20°	BHS07075	BHS07076				
			10°	BHS07077	BHS07078				
SIGNO VINCES COMPATIBLE									
DUO / INTEGRA									
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098	
			20°	BHS07004	BHS07005				
			10°	BHS07007	BHS07008				

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DUOCON / CM								
	WP_4.5 / 5.0		30°	BHS07130	BHS07131	BHS02009	BHS02010	EM046370
			20°	BHS07132	BHS07133			
			10°	BHS07134	BHS07135			
STRAUMANN COMPATIBLE								
OCTA INTERN								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
	SYN_RN	Ø 4,8	20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
			30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375
			20°	BHS07037	BHS07038			
			10°	BHS07052	BHS07053			
BONE LEVEL								
	NC	Ø 3,3	30°	BHS07085	BHS07086	BHS02015	BHS02016	EM048048
	RC	Ø 4,1 / 4,8	20°	BHS07087	BHS07088			
			10°	BHS07089	BHS07090			
			30°	BHS07027	BHS07028	BHS02015	BHS02016	EM046904
			20°	BHS07042	BHS07043			
			10°	BHS07057	BHS07058			
BL MULTI BASE								
	RC	Ø 4,8	30°		BHS07026	BHS02017	BHS02018	
		20°		BHS07041				
		10°		BHS07056				
SWEDEN & MARTINA COMPATIBLE								
PREMIUM KOHNO								
	NP	Ø 3,3	30°	BHS07200	BHS07201	BHS02021	BHS02022	
	RP	Ø 3,8	20°	BHS07202	BHS07203			
			10°	BHS07204	BHS07205			
			30°	BHS07206	BHS07207	BHS02021	BHS02022	
			20°	BHS07208	BHS07209			
			10°	BHS07210	BHS07211			
WP	Ø 4,25	30°	BHS07212	BHS07213	BHS02011	BHS02012		
		20°	BHS07214	BHS07215				
		10°	BHS07216	BHS07217				
SYRA (EH)								
		Ø4,10 / 4,25	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
		Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SYRA SHORT (EH)								
		Ø4,10	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
		Ø 5,0	30°	BHS07061	BHS07062	BHS02001	BHS02002	EM045082
			20°	BHS07063	BHS07064			
			10°	BHS07065	BHS07066			
TRI DENTAL IMPLANTS COMPATIBLE								
TRI-VENT								
		Ø 3,75 / 4,1 / 4,7	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
TRINON COMPATIBLE								
TRINON - Q2								
	RP	Ø 4,1	30°	BHS07001	BHS07002	BHS02001	BHS02002	EM045098
			20°	BHS07004	BHS07005			
			10°	BHS07007	BHS07008			
TRINON - QK								
	RN	Ø 4,8	30°	BHS07020	BHS07021	BHS02011	BHS02012	EM045144
			20°	BHS07035	BHS07036			
			10°	BHS07050	BHS07051			
	RN	Ø 4,8	30°	BHS07022	BHS07023	BHS02013	BHS02014	EM046375
			20°	BHS07037	BHS07038			
			10°	BHS07052	BHS07053			
ZIACOM COMPATIBLE								
ZINIC - IH								
	NP	Ø 3,2	30°	BHS07230	BHS07231	BHS02041	BHS02042	
			20°	BHS07232	BHS07233			
			10°	BHS07234	BHS07235			
	RP	Ø 3,5	30°	BHS07073	BHS07074	BHS02021	BHS02022	EM045130
			20°	BHS07075	BHS07076			
			10°	BHS07077	BHS07078			
	WP	Ø 4,5	30°	BHS07024	BHS07025	BHS02021	BHS02022	EM046218
			20°	BHS07039	BHS07040			
			10°	BHS07054	BHS07055			

1.1 DIRECT CASTING

	PLATFORM	Ø	∠	CASTABLE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
OEX - EH									
	NP	Ø 3,3	30° 20° 10°	BHS07136 BHS07138 BHS07140	BHS07137 BHS07139 BHS07141	BHS02035	BHS02036	EM048022	
	RP	Ø 4,1	30° 20° 10°	BHS07001 BHS07004 BHS07007	BHS07002 BHS07005 BHS07008	BHS02001	BHS02002	EM045098	
	OIN - IH (M2)								
	RP	Ø 4,0	30° 20° 10°	BHS07338 BHS07340 BHS07342	BHS07339 BHS07341 BHS07343	BHS02011	BHS02012		
	WP	Ø 5,0	30° 20° 10°	BHS07338 BHS07340 BHS07342	BHS07339 BHS07341 BHS07343	BHS02011	BHS02012		
	XDrive (Multi-unit)								
	RP	Ø 4,8	30° 20° 10°		BHS07003 BHS07006 BHS07009	BHS02003	BHS02004	EM045169	
ZIMMER COMPATIBLE									
INTERNAL HEX									
	NP	Ø 3,5	30° 20° 10°	BHS07073 BHS07075 BHS07077	BHS07074 BHS07076 BHS07078	BHS02021	BHS02022	EM045130	
	RP	Ø 4,5	30° 20° 10°	BHS07024 BHS07039 BHS07054	BHS07025 BHS07040 BHS07055	BHS02021	BHS02022	EM046218	
	WP	Ø 5,7	30° 20° 10°	BHS07166 BHS07168 BHS07170	BHS07167 BHS07169 BHS07171	BHS02021	BHS02022		
	SWISS-PLUS								
		RP	Ø 4,8	30° 20° 10°	BHS07020 BHS07035 BHS07050	BHS07021 BHS07036 BHS07051	BHS02011	BHS02012	EM45144
		RP	Ø 4,8	30° 20° 10°	BHS07022 BHS07037 BHS07052	BHS07023 BHS07038 BHS07053	BHS02013	BHS02014	EM46375

1.2 OVERCASTING TO CrCo BASE



PRECISION AND SAFETY FOR YOUR RESTORATIONS

- Continuous surface for an easy wax-up, casting & overcasting.
- Conical design for multiple restorations & engaging flat wall for single ones.
- Optimal emergence profile. Soft tissue care.
- Optimized access screw channel design.
- 10, 20 & 30° angle correction.
- Optimum height allowing 30° angulation for BHS screwdriver.
- Guide grip for BHS castable axial retention & free rotational movement.
- Retention pins & guide grip that enables axial retention.
- Engaging & non-engaging design for multiple and single restorations.

1.2 OVERCASTING TO CrCo BASE

		PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE									
I SERIES INTERNAL HEX									
	NARROW	Ø 3,0 - 3,75	30°	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	EM045130
			20°	BHS04057					
			10°	BHS04058					
	STANDARD	Ø 3,5 - 6,0	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	
			20°	BHS04057					
			10°	BHS04058					
I 10C CONUS									
	RP	Ø 4,2	30°	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	EM048149
			20°	BHS04105					
			10°	BHS04106					
ACE COMPATIBLE									
INFINITY EXTERNAL HEX									
	RP	Ø 4,0	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
			20°	BHS04005					
			10°	BHS04009					
INFINITY INTERNAL HEX									
	NP	Ø 3,7 / 4,1	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
			20°	BHS04057					
			10°	BHS04058					
	RP	Ø 4,7 / 5,1	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
			20°	BHS04039					
			10°	BHS04053					
ADIN COMPATIBLE									
SWELL									
		Ø 3,5 - 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
			20°	BHS04057					
			10°	BHS04058					
TOUAREG S / OS									
		Ø 3,5 - 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
			20°	BHS04057					
			10°	BHS04058					

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOUAREG CLOSEFIT								
	RP Ø 3,5	30°	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	EM048148
		20°	BHS04069					
		10°	BHS04070					
	WP Ø 4,3 / 5	30°	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	EM048149
		20°	BHS04105					
		10°	BHS04106					
ALPHABIO COMPATIBLE								
EXTERNAL ALPHABIO								
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
INTERNAL ALPHABIO								
	Ø 3,3 - 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
ALPHA DENT COMPATIBLE								
ACTIVE								
	Ø 3,75 / 4,2 / 5,0	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
ACTIVE PLUS								
	Ø 4,2 / 5,0	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
CLASSIC								
	Ø 3,3 / 3,75 / 4,2 / 5,0	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
ANCLADEN COMPATIBLE								
ANCLALOCK								
	Ø 3,75 - 5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					

1.2 OVERCASTING TO CrCo BASE

		PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA							
ANTHOGRYR COMPATIBLE																
ANTHOFIT EH																
   	RP	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098							
			20°	BHS04005												
			10°	BHS04009												
ARDS COMPATIBLE																
SMART																
   		Ø 3,75 - 4,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130							
			20°	BHS04057												
			10°	BHS04058												
CLASSIC																
   		Ø 3,3 - 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130							
			20°	BHS04057												
			10°	BHS04058												
PREMIUM																
   		Ø 3,3 - 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130							
			20°	BHS04057												
			10°	BHS04058												
ASTRA COMPATIBLE																
AQUA																
   	RP_3.5 / 4.0		30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371							
			20°	BHS04087												
			10°	BHS04088												

1.2 OVERCASTING TO CrCo BASE

		PLATFORM Ø		CrCo CASTABLE	 CrCo BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
EV									
		Ø 3,6	30°	BHS04056	BHS03154	BHS03155	BHS02053	BHS02054	
			20°	BHS04057					
			10°	BHS04058					
		Ø 4,2	30°	BHS04056	BHS03158	BHS03159	BHS02055	BHS02056	
			20°	BHS04057					
			10°	BHS04058					
AVINENT COMPATIBLE									
EH									
		RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
			20°	BHS04005					
			10°	BHS04009					
IH									
		Ø 3,5 / 4,1	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
			20°	BHS04057					
			10°	BHS04058					
MULTIUNIT									
		RP Ø 4,8	30°	BHS04002		BHS03008	BHS02003	BHS02004	EM045169
			20°	BHS04006					
			10°	BHS04010					
OCEAN CONICAL (C.C)									
		Ø 3,5 (Implant Ø 3,5 / Ø 4,0)	30°	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	EM048148
			20°	BHS04069					
			10°	BHS04070					
		Ø 4,1 (Implant Ø 4,5 / Ø 5,0)	30°	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	EM048149
			20°	BHS04105					
			10°	BHS04106					
BEGO COMPATIBLE									
IH									
		Ø 3,75 / 4,1	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
			20°	BHS04057					
			10°	BHS04058					

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BIOGENESIS COMPATIBLE								
3ICON								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
ATICON								
	RP Ø 3,5 / 4,0	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
KISSES								
	RP Ø 3,5 / 4,0	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
BIOHORIZONS COMPATIBLE								
INTERNAL								
	NP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
	RP Ø 4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
		20°	BHS04039					
		10°	BHS04053					
BIOMET 3i COMPATIBLE								
EH OSSEOTITE								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA						
BIONER COMPATIBLE														
IKELT / BIKELT														
	RP	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098					
			20°							BHS04005				
			10°							BHS04009				
HIKELT														
	NP	Ø 3,8	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130					
			20°	BHS04057										
			10°	BHS04058										
	RP	Ø 4,7	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218					
			20°	BHS04039										
			10°	BHS04053										
TOP DM														
	Ø 3,5 / 4,0 / 5,0	30°	30°	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	EM048148					
			20°	BHS04069										
			10°	BHS04070										
MULTIUNIT														
	RP	Ø 4,8	30°	BHS04002		BHS03008	BHS02003	BHS02004	EM045169					
			20°	BHS04006										
			10°	BHS04010										
BTI COMPATIBLE														
EH														
	RP	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098					
			20°	BHS04005										
			10°	BHS04009										
BTK COMPATIBLE														
EVO DL EXT														
	RP	Ø 4,1 / 4,6	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098					
			20°	BHS04005										
			10°	BHS04009										

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
EVO DL INT								
	Ø 3,5 / 3,7 / 4,6 / 5,0 / 5,6	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
KLASSIC EXT								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
KLASSIC INT								
	Ø 3,5 / 4,5 / 5,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
KONIC EXT								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
KONIC INT								
	Ø 3,5 / 5,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
B&W COMPATIBLE								
HEXAGONO EXT								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
HEXAGONO INT								
	NP Ø 3,3	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
	RP Ø 4,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
	WP Ø 5,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130

1.2 OVERCASTING TO CrCo BASE

		PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA		
CAMLOG COMPATIBLE											
CAMLOG											
	Ø3,3	30º	BHS04056	BHS03130	BHS03131	BHS02049	BHS02050				
		20º	BHS04057								
		10º	BHS04058								
	Ø3,8	30º	BHS04056	BHS03134	BHS03135	BHS02049	BHS02050				
		20º	BHS04057								
		10º	BHS04058								
	Ø4,3	30º	BHS04026	BHS03138	BHS03139	BHS02049	BHS02050				
		20º	BHS04039								
		10º	BHS04053								
	Ø5,0	30º	BHS04065	BHS03142	BHS03143	BHS02051	BHS02052				
		20º	BHS04066								
		10º	BHS04067								
CONEXÃO SIST. PROTESE COMPATIBLE											
MORSE INDEXADO											
	Ø 3,5 / 3,75 / 4,0 / 5,0	30º	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	EM048148			
		20º	BHS04069								
		10º	BHS04070								
		HEXAGONO EXT									
		Ø 3,75 / 4,0	30º	BHS04001	BHS03003	BHS03004	BHS02001		BHS02002	EM045098	
	20º		BHS04005								
	10º		BHS04009								
	HEXAGONO INT										
	Ø 3,75 / 4,0 / 5,0		30º	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218		
		20º	BHS04039								
		10º	BHS04053								

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CORTEX COMPATIBLE								
DYNAMIX								
	Ø 3,3 / 3,8 / 4,2 / 5	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
CLASSIX								
	Ø 3,3 / 3,8 / 4,2 / 5	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
SATURN								
	Ø 3,8 / 4,2	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
COWELLMEDI COMPATIBLE								
EXTERNAL								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
SUBMERGED (INT)								
	Ø 3,5 / 4 / 4,5 / 5 / 6	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03102	BHS03103	BHS02009	BHS02010	
DENTALTECH COMPATIBLE								
IMPLAGIC GII								
	RP Ø 4,5	30° 20° 10°	BHS04026 BHS04039 BHS04053	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
IMPLASSIC HX2								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA		
DENTIS COMPATIBLE										
S - CLEAN										
	RP	Ø 3,4	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010		
			20°	BHS04057						
			10°	BHS04058						
DENTIUM COMPATIBLE										
IMPLANTIUM										
		Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010		
			20°	BHS04057						
			10°	BHS04058						
SUPER LINE										
				Ø 3,6 / 4 / 4,5 / 5 / 6	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010
	20°	BHS04057								
	10°	BHS04058								
DENTSPLY SIRONA COMPATIBLE										
XIVE S										
	D 3.4	30°	BHS04056	BHS03118	BHS03119	BHS02047	BHS02048			
		20°	BHS04057							
		10°	BHS04058							
	D 3.8	30°	BHS04056	BHS03122	BHS03123	BHS02047	BHS02048			
		20°	BHS04057							
		10°	BHS04058							
	D 4.5	30°	BHS04056	BHS03126	BHS03127	BHS02047	BHS02048			
		20°	BHS04057							
		10°	BHS04058							

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DIO IMPLANTS COMPATIBLE								
EXTERNAL								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
UF II								
	Ø 3,8 / 4 / 4,5 / 5 / 5,5	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03102	BHS03103	BHS02009	BHS02010	
EDISON MEDICAL COMPATIBLE								
LAMINA (IH)								
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
FORTIS (IH)								
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
ELITE MEDICA COMPATIBLE								
EXTERNAL CONNECT								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
EUROTEKNIKA COMPATIBLE								
UNEVA								
	Ø 3,6 / 4,1 / 4,8	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
NATURACTIS								
	Ø 3,5 / 4 / 4,5 / 5	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
NATURALL +								
	Ø 3,5 / 4 / 4,5 / 5	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
NATEA								
	NP Ø 3,6 / 4,1	30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
		20°	BHS04087					
		10°	BHS04088					
	RP Ø 4,8 / 6	30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
		20°	BHS04087					
		10°	BHS04088					
GALIMPLANT COMPATIBLE								
IPXE- EXTERNA								
	RP Ø 3,5 / 4 / 5	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
ICI - INTERNA								
	RP Ø 3,2 / 3,5 / 4 / 5	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371

1.2 OVERCASTING TO CrCo BASE

			PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA							
IPX - INTERNA																	
	RP	Ø 3,5 / 4 / 4,5 / 5	30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371								
			20°	BHS04087													
			10°	BHS04088													
GMI COMPATIBLE																	
FRONTIER (IH)																	
	RP Standard	Ø 3,3	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130								
			20°	BHS04057													
			10°	BHS04058													
	WP Wide	Ø 4,3	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218								
			20°	BHS04039													
			10°	BHS04053													
FRONTIER multi-aesthetic																	
	RP	Ø 4,8	30°	BHS04002		BHS03008	BHS02003	BHS02004	EM045169								
			20°	BHS04006													
			10°	BHS04010													
PHOENIX (EH)																	
	Standard	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098								
			20°	BHS04005													
			10°	BHS04009													
PHOENIX multi-aesthetic																	
	RP	Ø 4,8	30°	BHS04002		BHS03008	BHS02003	BHS02004	EM045169								
			20°	BHS04006													
			10°	BHS04010													
GT MEDICAL COMPATIBLE																	
BEST FIT HXE																	
	RP	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098								
			20°	BHS04005													
			10°	BHS04009													

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BEST FIT 002								
   	NP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
HI TECH COMPATIBLE								
SELF THREAD								
   	SLT Ø 3,3 / 3,75	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
   	GTI Ø 4,2 / 5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
		20°	BHS04039					
		10°	BHS04053					
LGI PLUS								
   	Ø 3,5 / 4,3 / 5 / 6	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
IMPLANT DIRECT COMPATIBLE								
INTERACTIVE								
   	NP Ø 3,5	30°	BHS04068	BHS03058	BHS03059	BHS02025	BHS02026	EM048148
		20°	BHS04069					
		10°	BHS04070					
   	RP Ø 4,3	30°	BHS04104	BHS03046	BHS03047	BHS02023	BHS02024	EM048149
		20°	BHS04105					
		10°	BHS04106					
LEGACY								
   	NP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
   	RP Ø 4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
		20°	BHS04039					
		10°	BHS04053					

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
JDEENTALCARE COMPATIBLE								
JDEVOLUTION								
	Ø 3,7 / 4,3 / 5 / 6	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
KEYSTONE DENTAL COMPATIBLE								
RESTORE								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
RENOVA								
	NP Ø 3,75	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
	RP Ø 4,5 / 4,75	30° 20° 10°	BHS04026 BHS04039 BHS04053	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
KLOCKNER COMPATIBLE								
KL (EH)								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
LEADER COMPATIBLE								
TIXOS INTERNAL HEX								
	NP Ø 3,5 / 4,0	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA		
TIXOS EXTERNAL HEX										
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098		
		20°	BHS04005							
		10°	BHS04009							
MEDENTIS COMPATIBLE										
ICX										
	Ø 3,45 / 3,7 - 4,1 / 4,8	30°	BHS04056	BHS03180	BHS03181	BHS02015	BHS02016			
		20°	BHS04057							
		10°	BHS04058							
MEGAGEN COMPATIBLE										
ANYONE INTERNAL										
	NP	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010			
		20°	BHS04057							
		10°	BHS04058							
	RP	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010			
		20°	BHS04057							
		10°	BHS04058							
	WP	30°	BHS04056	BHS03102	BHS03103	BHS02009	BHS02010			
		20°	BHS04057							
		10°	BHS04058							
ANYONE EXTERNAL										
	RP	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098		
		20°	BHS04005							
		10°	BHS04009							
ANYRIDGE										
	Small Ø3,5	30°	BHS04056	BHS03166	BHS03167	BHS02027	BHS02028			
		20°	BHS04057							
		10°	BHS04058							
	Regular Ø4,0/4,5	30°	BHS04056	BHS03166	BHS03167	BHS02027	BHS02028			
		20°	BHS04057							
		10°	BHS04058							
	Wide Ø5,0/5,5	30°	BHS04056	BHS03166	BHS03167	BHS02027	BHS02028			
		20°	BHS04057							
		10°	BHS04058							

1.2 OVERCASTING TO CrCo BASE

			PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
MICRODENT COMPATIBLE										
EKTOS (IH)										
	NP	Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130	
			20°	BHS04057						
			10°	BHS04058						
	RP	Ø 4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218	
			20°	BHS04039						
			10°	BHS04053						
SUV (EH)										
	SUV	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098	
			20°	BHS04005						
			10°	BHS04009						
UNIVERSAL (EH)										
	SU	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098	
			20°	BHS04005						
			10°	BHS04009						
MIS COMPATIBLE										
LANCE										
	STD	Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098	
			20°	BHS04005						
			10°	BHS04009						
SEVEN / M4										
	NARROW	Ø 3,3	30°	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	EM045130	
			20°	BHS04057						
			10°	BHS04058						
	STANDARD	Ø 3,75 / 4,2	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130	
			20°	BHS04057						
			10°	BHS04058						
	WIDE	Ø 5,6 / 6,0	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218	
			20°	BHS04039						
			10°	BHS04053						

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
MOZO GRAU COMPATIBLE								
INHEX								
	RP Ø 3,75 / 4,25	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
OSSEUS								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
MPI COMPATIBLE								
EXTERNA EH								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
INTERNA CM								
	RP Ø 4,0	30° 20° 10°	BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
NEOBIOTECH COMPATIBLE								
EB EXTERNAL SYSTEM								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
IS INTERNAL SYSTEM								
	Ø 3,5 / 4 / 4,5 / 5	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03102	BHS03103	BHS02009	BHS02010	

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NEODENT COMPATIBLE								
TITAMAX HE SMART								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
MINI PILAR CM								
	RP Ø 4,8	30° 20° 10°	BHS04002 BHS04006 BHS04010		BHS03008	BHS02003	BHS02004	EM045169
NOBEL COMPATIBLE								
BRÄNEMARK								
	RP Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
ACTIVE								
	NP Ø 3,5	30° 20° 10°	BHS04068 BHS04069 BHS04070	BHS03058	BHS03059	BHS02025	BHS02026	EM048148
	RP Ø 4,3 / Ø5,0	30° 20° 10°	BHS04104 BHS04105 BHS04106	BHS03046	BHS03047	BHS02023	BHS02024	EM048149
MULTIUNIT								
	RP Ø 4,8	30° 20° 10°	BHS04002 BHS04006 BHS04010		BHS03008	BHS02003	BHS02004	EM045169

1.2 OVERCASTING TO CrCo BASE

			PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
NT TRADING COMPATIBLE											
2-CONnect											
	M1,4	Ø5,0	30° 20° 10°	BHS04002 BHS04006 BHS04010			BHS03177	BHS02063	BHS02064		
	M1,6	Ø5,0	30° 20° 10°	BHS04002 BHS04006 BHS04010			BHS03177	BHS02065	BHS02066		
	M1,8	Ø5,0	30° 20° 10°	BHS04002 BHS04006 BHS04010			BHS03177	BHS02067	BHS02068		
	OSSTEM IMPLANT COMPATIBLE										
	OSSTEM TS										
		Ø 4,0 / 4,5 / 5,0 / 6,0		30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03102	BHS03103	BHS02009	BHS02010		
OSSTEM US											
	RP	Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098		
OSTEOPLUS COMPATIBLE											
OSTEOPLUS - SHe											
	RP	Ø 4,1	30° 20° 10°	BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098		
OSTEOPLUS - SHi											
	RP	Ø 3,5	30° 20° 10°	BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130		
	WP	Ø 4,5	30° 20° 10°	BHS04026 BHS04039 BHS04053	BHS03033	BHS03034	BHS02021	BHS02022	EM046218		

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OXTEIN COMPATIBLE								
HEX EXTERNA								
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
HEX INTERNA								
	NP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
	RP Ø 4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
		20°	BHS04039					
		10°	BHS04053					
PALTOP COMPATIBLE								
INTERNAL HEX CONNECTION								
	NARROW Ø 3,25	30°	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	
		20°	BHS04057					
		10°	BHS04058					
	STANDARD Ø 3,75 / 4,2 / 5,0	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
	WIDE Ø 6,0	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
		20°	BHS04039					
		10°	BHS04053					
PCA Implant System (Paltop Conical Active)								
	NARROW Ø 3,25	30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
		20°	BHS04087					
		10°	BHS04088					
	STANDARD Ø 3,75 / 4,2 / 5,0	30°	BHS04086	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
		20°	BHS04087					
		10°	BHS04088					

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHIBO COMPATIBLE								
TSH								
	S3 / S4	Ø 4,1	30° 20° 10° BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
AUREA EVO								
	NP	Ø 3,5	30° 20° 10° 0° BHS04056 BHS04057 BHS04058 BHS04055	BHS03106	BHS03107	BHS02041	BHS02042	
	RP	Ø 4,3	30° 20° 10° 0° BHS04056 BHS04057 BHS04058 BHS04055	BHS03110	BHS03111	BHS02021	BHS02022	
RADHEX IMPLANTS COMPATIBLE								
PHE / PHEA								
	RP	Ø 4,1	30° 20° 10° BHS04001 BHS04005 BHS04009	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
PHI / PHIA								
	NP	Ø 3,5	30° 20° 10° BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
	RP	Ø 4,5	30° 20° 10° BHS04026 BHS04039 BHS04053	BHS03033	BHS03034	BHS02021	BHS02022	EM046218
PCI								
	280		30° 20° 10° BHS04086 BHS04087 BHS04088	BHS03070	BHS03071	BHS02025	BHS02026	EM046371
SGS COMPATIBLE								
P1 / P7								
	RP	Ø 3,75	30° 20° 10° BHS04056 BHS04057 BHS04058	BHS03050	BHS03051	BHS02021	BHS02022	EM045130

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SIGNO VINCES COMPATIBLE								
DUO / INTEGRA								
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
SWEDEN & MARTINA COMPATIBLE								
SYRA (EH)								
	Ø 4,10 / 4,25	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
SYRA Short (EH)								
	Ø 4,10	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
TRI DENTAL IMPLANTS COMPATIBLE								
TRI - VENT								
	Ø 3,75 / 4,1 / 4,7	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
TRINON COMPATIBLE								
TRINON - Q2								
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					

1.2 OVERCASTING TO CrCo BASE

	PLATFORM Ø		CrCo CASTABLE	 CrCo BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ZIACOM COMPATIBLE								
ZINC - IH								
	NP Ø 3,2	30°	BHS04056	BHS03200	BHS03201	BHS02041	BHS02042	EM045130
		20°	BHS04057					
		10°	BHS04058					
	RP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	
		20°	BHS04057					
		10°	BHS04058					
	WP Ø4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	
		20°	BHS04039					
		10°	BHS04053					
OEX - EH								
	RP Ø 4,1	30°	BHS04001	BHS03003	BHS03004	BHS02001	BHS02002	EM045098
		20°	BHS04005					
		10°	BHS04009					
OIN - IH (M2)								
	RP Ø4,0	30°	BHS04056	BHS03050	BHS03051	BHS02011	BHS02012	EM045169
		20°	BHS04057					
		10°	BHS04058					
	WP Ø5,0	30°	BHS04026	BHS03033	BHS03034	BHS02011	BHS02012	
		20°	BHS04039					
		10°	BHS04053					
XDrive (Multi-unit)								
	RP Ø 4,8	30°	BHS04002		BHS03008	BHS02003	BHS02004	EM045169
		20°	BHS04006					
		10°	BHS04010					
ZIMMER COMPATIBLE								
INTERNAL HEX								
	NP Ø 3,5	30°	BHS04056	BHS03050	BHS03051	BHS02021	BHS02022	EM045130
		20°	BHS04057					
		10°	BHS04058					
	RP Ø4,5	30°	BHS04026	BHS03033	BHS03034	BHS02021	BHS02022	
		20°	BHS04039					
		10°	BHS04053					

1.3 CEMENTING TO Ti BASE



100 % PASSIVE FIT IN HYBRID RESTORATIONS

Continuous surface to facilitate restorations.

Conical design for multiple restorations & engaging flat wall for single ones.

Optimal emergence profile. Soft tissue care.

Optimized screw access channel design.

10, 20 & 30° angle correction.

U-Formed design, enabling smooth passage for the screw & driver up to 30°, while providing optimal height for cementation.


Retention pins & guide grip that enables axial retention.

Engaging & non-engaging design for single and multiple restorations.

1.3 CEMENTING TO Ti BASE

		PLATFORM Ø		Ti CASTABLE	Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE									
I SERIES INTERNAL HEX									
   	NARROW	Ø 3 - 3,75	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03198	BHS03199	BHS02041	BHS02042	
	STD	Ø 3,5 - 6,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
I 10C CONUS									
   	RP	Ø 4,2	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03044	BHS03045	BHS02023	BHS02024	EM048149
ACE COMPATIBLE									
INFINITY EXTERNAL HEX									
   	NP	Ø 3,3/3,75	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
	RP	Ø 4,0	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
INFINITY TRI - CAM									
   	NP	Ø 3,5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03060	BHS03061	BHS02027	BHS02028	EM046291
	RP	Ø 4,3	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03012	BHS03013	BHS02009	BHS02010	EM046292
INFINITY INTERNAL HEX									
   	NP	Ø 3,7 / 4,1	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
	RP	Ø 4,7 / 5,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03020	BHS03021	BHS02021	BHS02022	EM046218


1.3 CEMENTING TO Ti BASE

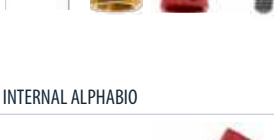	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ADIN COMPATIBLE								
SWELL								
	Ø 3,5 - 6	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
TOUAREG S / OS								
	Ø 3,5 - 6	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
TOUAREG CLOSEFIT								
	RP Ø 3,5	30°	BHS04059	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
		20°	BHS04060					
		10°	BHS04061					
	WP Ø 4,3 / 5	30°	BHS04018	BHS03044	BHS03045	BHS02023	BHS02024	EM048149
		20°	BHS04019					
		10°	BHS04020					
ALPHABIO COMPATIBLE								
EXTERNAL ALPHABIO								
	NP Ø 3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
		20°	BHS04060					
		10°	BHS04061					
	RP Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20°	BHS04019					
		10°	BHS04020					
INTERNAL ALPHABIO								
	Ø 3,3 - 6	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ALPHA DENT COMPATIBLE								
ACTIVE								
	Ø 3,75 / 4,2 / 5,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
ACTIVE PLUS								
	Ø 4,2 / 5,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
CLASSIC								
	Ø 3,3 / 3,75 / 4,2 / 5,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
ANCLADEN COMPATIBLE								
ANCLALOCK								
	Ø 3,75 - 5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060					
		10°	BHS04061					
ANTHOGRYR COMPATIBLE								
ANTHOFIT EH								
	RP Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20°	BHS04019					
		10°	BHS04020					
OSSFIT								
	RP Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
		20°	BHS04022					
		10°	BHS04023					
	SYN_RP Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375
		20°	BHS04081					
		10°	BHS04082					

1.3 CEMENTING TO Ti BASE

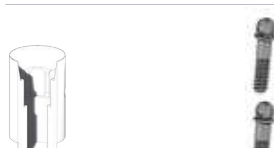	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ARDS COMPATIBLE								
SMART								
	Ø 3,75 - 4,5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
CLASSIC								
	Ø 3,3 - 6	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
PREMIUM								
	Ø 3,3 - 6	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
ASTRA COMPATIBLE								
AQUA								
	RP_3,5 / 4,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
CONO 20°								
	Ø3,5	30° 20° 10°	BHS04059 BHS04060 BHS04061		BHS03072	BHS02003	BHS02004	EM045638
CONO 45°								
	Ø3,5	30° 20° 10°	BHS04059 BHS04060 BHS04061		BHS03074	BHS02003	BHS02004	EM047786
EV								
	Ø3,6	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03152	BHS03153	BHS02053	BHS02054	
	Ø4,2	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03156	BHS03157	BHS02055	BHS02056	

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
AVINENT COMPATIBLE									
HE									
	NP	Ø 3,5	30°	BHS04059	BHS03010	BHS03011	BHS02031	BHS02032	EM045222
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
IH									
		Ø 3,5 / 4,1	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060						
		10°	BHS04061						
MULTIUNIT									
	RP	Ø 4,8	30°	BHS04004		BHS03009	BHS02003	BHS02004	EM045169
			20°	BHS04008					
			10°	BHS04012					
OCEAN CONICAL (C.C.)									
		Ø 3,5 (Implant Ø 3,5/ Ø4,0)	30°	BHS04059	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
		20°	BHS04060						
		10°	BHS04061						
		Ø 4,1 (Implant Ø 4,5/ Ø5,0)	30°	BHS04018	BHS03044	BHS03045	BHS02023	BHS02024	EM048149
		20°	BHS04019						
		10°	BHS04020						
BEGO COMPATIBLE									
IH									
		Ø 3,75 / 4,1	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20°	BHS04060						
		10°	BHS04061						

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
BIOGENESIS COMPATIBLE									
3ICON									
	RP	Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
ITICON									
	RP	Ø 4,8	30° 20° 10°	BHS04021 BHS04022 BHS04023	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
	RP	Ø 4,8	30° 20° 10°	BHS04080 BHS04081 BHS04082		BHS03019	BHS02013	BHS02014	EM046375
ATICON									
	RP	Ø 3,5 / 4,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
KISSES									
	RP	Ø 3,5 / 4,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
CONO 20°									
	RP/WP		30° 20° 10°	BHS04059 BHS04060 BHS04061		BHS03072	BHS02003	BHS02004	EM045638
CONO 45°									
	RP/WP		30° 20° 10°	BHS04059 BHS04060 BHS04061		BHS03074	BHS02003	BHS02004	EM047786
BIOHORIZONS COMPATIBLE									
INTERNAL									
	NP	Ø 3,5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
	RP	Ø 4,5	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03020	BHS03021	BHS02021	BHS02022	EM046218

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
BIOMET 3i COMPATIBLE									
EH OSSEOTITE									
	RP	Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
BIONER COMPATIBLE									
IKELT / BIKELT									
	RP	Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
HIKELT									
	NP	Ø 3,8	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
	RP	Ø 4,7	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
TOP DM									
	Ø 3,5 / 4,0 / 5,0		30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
MULTIUNIT									
	RP	Ø 4,8	30° 20° 10°	BHS04004 BHS04008 BHS04012		BHS03009	BHS02003	BHS02004	EM045169
BTI COMPATIBLE									
EH									
	RP	Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BTK COMPATIBLE								
8NECK INT								
	SR Ø 4,8	30° 20° 10°	BHS04021 BHS04022 BHS04023	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
EVO DL EXT								
	RP Ø 4,1 / 4,6	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
EVO DL INT								
	Ø 3,5 / 3,7 / 4,6 / 5,0 / 5,6	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
KLASSIC EXT								
	RP Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
KLASSIC INT								
	Ø 3,5 / 4,5 / 5,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
KONIC EXT								
	RP Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
KONIC INT								
	Ø 3,5 / 5,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
SAFETL INT								
	SR Ø 4,8	30° 20° 10°	BHS04021 BHS04022 BHS04023	BHS03016	BHS03017	BHS02011	BHS02012	EM045144

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
B&W COMPATIBLE									
HEXAGONO EXT									
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
HEXAGONO INT									
	NP	Ø 3,3	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	WP	Ø 5,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
CAMLOG COMPATIBLE									
CAMLOG									
		Ø 3,3	30°	BHS04089	BHS03128	BHS03129	BHS02049	BHS02050	
			20°	BHS04090					
			10°	BHS04091					
		Ø 3,8	30°	BHS04089	BHS03132	BHS03133	BHS02049	BHS02050	
			20°	BHS04090					
			10°	BHS04091					
		Ø 4,3	30°	BHS04095	BHS03136	BHS03137	BHS02049	BHS02050	
			20°	BHS04096					
			10°	BHS04097					
		Ø 5,0	30°	BHS04083	BHS03140	BHS03141	BHS02051	BHS02052	
			20°	BHS04084					
			10°	BHS04085					
CONEXÃO SIST. PROTESE COMPATIBLE									
MORSE INDEXADO									
		Ø 3,5 / 3,75 / 4,0 / 5,0	30°	BHS04059	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
			20°	BHS04060					
			10°	BHS04061					
HEXAGONO EXT									
		Ø 3,75 / 4,0	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
HEXAGONO INT								
	Ø 3,75 / 4,0 / 5,0	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
CORTEX COMPATIBLE								
DYNAMIX								
	Ø 3,3 / 3,8 / 4,2 / 5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
CLASSIX								
	Ø 3,3 / 3,8 / 4,2 / 5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
SATURN								
	Ø 3,8 / 4,2	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
COWELLMEDI COMPATIBLE								
EXTERNAL								
	RP Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
INTERNAL								
	RN Ø 4,8	30° 20° 10°	BHS04021 BHS04022 BHS04023	BHS03016	BHS03017	BHS02011	BHS02012	EM045144

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA						
SUBMERGED (INT)														
	Ø3,5/4/4,5/5/6	30°	BHS04059	BHS03100		BHS02009	BHS02010							
		20°	BHS04060											
		10°	BHS04061											
DENTALTECH COMPATIBLE														
IMPLAGIC GII														
	RP Ø 4,5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218						
		20°	BHS04019											
		10°	BHS04020											
IMPLASSIC HX2														
	RP Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098						
		20°	BHS04019											
		10°	BHS04020											
DENTIS COMPATIBLE														
S-CLEAN														
	RP Ø 3,4	30°	BHS04059	BHS03100		BHS02009	BHS02010							
		20°	BHS04060											
		10°	BHS04061											
I - CLEAN														
	RP Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144						
		20°	BHS04022											
		10°	BHS04023											
DENTIUM COMPATIBLE														
SIMPLE LINE II														
	RN Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144						
		20°	BHS04022											
		10°	BHS04023											

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
IMPLANTIUM									
	Ø3,6/4/4,4/4,9/6/7	30°	BHS04059	BHS03100			BHS02009	BHS02010	
		20°	BHS04060						
		10°	BHS04061						
SUPER LINE									
	Ø3,6/4/4,5/5/6	30°	BHS04059	BHS03100			BHS02009	BHS02010	
		20°	BHS04060						
		10°	BHS04061						
DENTSPLY SIRONA COMPATIBLE									
XIVES									
		D 3.4	30°	BHS04059	BHS03116	BHS03117	BHS02047	BHS02048	
			20°	BHS04060					
			10°	BHS04061					
		D 3.8	30°	BHS04059	BHS03120	BHS03121	BHS02047	BHS02048	
			20°	BHS04060					
			10°	BHS04061					
		D 4.5	30°	BHS04059	BHS03124	BHS03125	BHS02047	BHS02048	
			20°	BHS04060					
			10°	BHS04061					
DIO IMPLANTS COMPATIBLE									
EXTERNAL									
	NP Ø3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222	
		20°	BHS04060						
		10°	BHS04061						
	RP Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098	
		20°	BHS04019						
		10°	BHS04020						
UF II									
	Ø3,8/4/4,5/5/5,5	30°	BHS04059	BHS03100			BHS02009	BHS02010	
		20°	BHS04060						
		10°	BHS04061						

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø	∠	Ti CASTABLE	TI BASE ⬡	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
EDISON MEDICAL COMPATIBLE								
LAMINA (IH)								
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
FORTIS (IH)								
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
ELITE MEDICA COMPATIBLE								
EXTERNAL CONNECT								
	NP Ø 3,5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
EUROTEKNIKA COMPATIBLE								
UNEVA								
	Ø 3,6 / 4,1 / 4,8	30° 20° 10°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
NATURACTIS								
	Ø 3,5 / 4 / 4,5 / 5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
NATURALL +								
	Ø 3,5 / 4 / 4,5 / 5	30° 20° 10°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
NATEA									
	NP	Ø 3,6 / 4,1	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,8 / 6	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					
AESTHETICA									
		Ø 3,6 / 4,1 / 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
			20°	BHS04022					
			10°	BHS04023					
			GALIMPLANT COMPATIBLE						
IPXE - EXTERNA									
	RP	Ø 3,5 / 4 / 5	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
			ICI - INTERNA						
	RP	Ø 3,2 / 3,5 / 4 / 5	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					
			IPX - INTERNA						
	RP	Ø 3,5 / 4 / 4,5 / 5	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					

1.3 CEMENTING TO Ti BASE

				PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
GMI COMPATIBLE											
FRONTIER (IH)											
	RP Standard	Ø 3,3	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130		
			20°	BHS04060							
			10°	BHS04061							
	WP Wide	Ø 4,3	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218		
			20°	BHS04019							
			10°	BHS04020							
FRONTIER multi-aesthetic											
	RP	Ø 4,8	30°	BHS04004		BHS03009	BHS02003	BHS02004	EM045169		
			20°	BHS04008							
			10°	BHS04012							
			PHOENIX (EH)								
	Standard	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098		
			20°	BHS04019							
			10°	BHS04020							
			PHOENIX multi-aesthetic								
	RP	Ø 4,8	30°	BHS04004		BHS03009	BHS02003	BHS02004	EM045169		
			20°	BHS04008							
			10°	BHS04012							
			GT MEDICAL COMPATIBLE								
BEST FIT HXE											
	NP	Ø 3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222		
			20°	BHS04060							
			10°	BHS04061							
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098		
			20°	BHS04019							
			10°	BHS04020							
BEST FIT 002											
	NP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130		
			20°	BHS04060							
			10°	BHS04061							

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
BEST FIT OCT									
	RN	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
			20°	BHS04022					
			10°	BHS04023					
	SYN_RN	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375
			20°	BHS04081					
			10°	BHS04082					
HI TECH COMPATIBLE									
SELF THREAD									
	SLT	Ø 3,3 / 3,75	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	GTI	Ø 4,2 / 5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
			20°	BHS04019					
			10°	BHS04020					
LGI PLUS									
		Ø 3,5 / 4,3 / 5 / 6	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
IMPLANT DIRECT COMPATIBLE									
INTERACTIVE									
	NP	Ø 3,5	30°	BHS04059	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,3	30°	BHS04018	BHS03044	BHS03045	BHS02023	BHS02024	EM048149
			20°	BHS04019					
			10°	BHS04020					
LEGACY									
	NP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
			20°	BHS04019					
			10°	BHS04020					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø	∠	Ti CASTABLE	TI BASE ⬡	○	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
REPLUS / REPLANT									
  	NP	Ø 3,5	30°	BHS04059	BHS03060	BHS03061	BHS02027	BHS02028	EM046291
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,3	30°	BHS04018	BHS03012	BHS03013	BHS02009	BHS02010	EM046292
			20°	BHS04019					
			10°	BHS04020					
SWISHPLUS / SWISHPLANT									
 	RN	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
			20°	BHS04022					
			10°	BHS04023					
	SYN_RN	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375
			20°	BHS04081					
			10°	BHS04082					
JDENTALCARE COMPATIBLE									
JDEVOLUTION									
	Ø 3,7 / 4,3 / 5 / 6	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130	
		20°	BHS04060						
		10°	BHS04061						
KEYSTONE DENTAL COMPATIBLE									
RESTORE									
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
OCT-XT									
 	RN	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
			20°	BHS04022					
			10°	BHS04023					
	SYN_RN	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375
			20°	BHS04081					
			10°	BHS04082					
RENOVA									
 	NP	Ø 3,75	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,5 / 4,75	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
			20°	BHS04019					
			10°	BHS04020					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLOCKNER COMPATIBLE									
KL (EH)									
	NP	Ø3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
			20°	BHS04060					
			10°	BHS04061					
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
LEADER COMPATIBLE									
TIXOS INTERNAL HEX									
	NP	Ø 3,5 / 4,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
TIXOS EXTERNAL HEX									
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					
MEDENTIS COMPATIBLE									
ICX									
		Ø 3,45 / 3,7 - 4,1 / 4,8	30°	BHS04059	BHS03178	BHS03179	BHS02015	BHS02016	
			20°	BHS04060					
			10°	BHS04061					

1.3 CEMENTING TO Ti BASE

				PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
MEGAGEN COMPATIBLE											
ANYONE INTERNAL											
	NP	Ø	30°	BHS04059	BHS03100				BHS02009	BHS02010	
			20°	BHS04060							
			10°	BHS04061							
	RP		30°	BHS04059	BHS03100				BHS02009	BHS02010	
			20°	BHS04060							
			10°	BHS04061							
	WP		30°	BHS04059	BHS03100				BHS02009	BHS02010	
			20°	BHS04060							
			10°	BHS04061							
ANYONE EXTERNAL											
	NP	Ø	30°	BHS04059	BHS03010	BHS03011			BHS02005	BHS02006	EM045222
			20°	BHS04060							
			10°	BHS04061							
	RP		30°	BHS04018	BHS03005	BHS03006			BHS02001	BHS02002	EM045098
			20°	BHS04019							
			10°	BHS04020							
ANYRIDGE											
	Small	Ø3,5	30°	BHS04059	BHS03164	BHS03165			BHS02027	BHS02028	
			20°	BHS04060							
			10°	BHS04061							
	Regular	Ø4,0/4,5	30°	BHS04059	BHS03164	BHS03165			BHS02027	BHS02028	
			20°	BHS04060							
			10°	BHS04061							
	Wide	Ø5,0/5,5	30°	BHS04059	BHS03164	BHS03165			BHS02027	BHS02028	
			20°	BHS04060							
			10°	BHS04061							
OCTA Abutment											
	SYN_RN	Ø4,8	30°	BHS04080		BHS03019			BHS02013	BHS02014	EM046375
			20°	BHS04081							
			10°	BHS04082							
MICRODENT COMPATIBLE											
EKTOS (IH)											
	NP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049			BHS02021	BHS02022	EM045130
			20°	BHS04060							
			10°	BHS04061							
	RP	Ø 4,5	30°	BHS04018	BHS03020	BHS03021			BHS02021	BHS02022	EM046218
			20°	BHS04019							
			10°	BHS04020							

1.3 CEMENTING TO Ti BASE

			PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA							
SUV (EH) & UNIVERSAL (EH)																	
	SUV & SU	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098								
			20°	BHS04019													
			10°	BHS04020													
MIS COMPATIBLE																	
LANCE																	
	STD	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098								
			20°	BHS04019													
			10°	BHS04020													
SEVEN / M4																	
	Narrow	Ø 3,30	30°	BHS04059	BHS03198	BHS03199	BHS02041	BHS02042									
			20°	BHS04060													
			10°	BHS04061													
	Standard	Ø 3,75 / 4,2	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130								
			20°	BHS04060													
			10°	BHS04061													
	Wide	Ø 5,0 / 6,0	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218								
			20°	BHS04019													
			10°	BHS04020													
MOZO GRAU COMPATIBLE																	
INHEX																	
	RP	Ø 3,75 / 4,25	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371								
			20°	BHS04060													
			10°	BHS04061													
	CONO20		30°	BHS04059		BHS03072	BHS02003	BHS02004	EM045638								
			20°	BHS04060													
			10°	BHS04061													
	CONO45		30°	BHS04059		BHS03074	BHS02003	BHS02004	EM047786								
			20°	BHS04060													
			10°	BHS04061													
OSSEUS																	
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098								
			20°	BHS04019													
			10°	BHS04020													

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
MPI COMPATIBLE									
EXTERNA EH									
	RP	Ø 4,1	30°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°						
			10°						
INTERNA CM									
	RP	Ø 4,0	30°	BHS04059 BHS04060 BHS04061	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°						
			10°						
NEOBIOTECH COMPATIBLE									
EB EXTERNAL SYSTEM									
	NP	Ø 3,5	30°	BHS04059 BHS04060 BHS04061	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
			20°						
			10°						
	RP	Ø 4,1	30°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°						
			10°						
IS INTERNAL SYSTEM									
	Ø3,5/4/4,5/5		30°	BHS04059 BHS04060 BHS04061	BHS03100		BHS02009	BHS02010	
			20°						
			10°						
NEODENT COMPATIBLE									
TITAMAX EH SMART									
	RP	Ø 4,1	30°	BHS04018 BHS04019 BHS04020	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°						
			10°						
MINI PILAR CM									
	RP	Ø 4,8	30°	BHS04004 BHS04008 BHS04012		BHS03009	BHS02003	BHS02004	EM045169
			20°						
			10°						

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NOBEL COMPATIBLE								
BRÄNEMARK								
	NP Ø 3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
		20°	BHS04060					
		10°	BHS04061					
	RP Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20°	BHS04019					
		10°	BHS04020					
REPLACE								
	NP Ø 3,5	30°	BHS04059	BHS03060	BHS03061	BHS02027	BHS02028	EM046291
		20°	BHS04060					
		10°	BHS04061					
	RP Ø 4,3	30°	BHS04018	BHS03012	BHS03013	BHS02009	BHS02010	EM046292
		20°	BHS04019					
		10°	BHS04020					
ACTIVE								
	NP Ø 3,5	30°	BHS04059	BHS03056	BHS03057	BHS02025	BHS02026	EM048148
		20°	BHS04060					
		10°	BHS04061					
	RP Ø 4,3 / Ø 5,0	30°	BHS04018	BHS03044	BHS03045	BHS02023	BHS02024	EM048149
		20°	BHS04019					
		10°	BHS04020					
MULTIUNIT								
	RP Ø 4,8	30°	BHS04004		BHS03009	BHS02003	BHS02004	EM045169
		20°	BHS04008					
		10°	BHS04012					
NT TRADING COMPATIBLE								
2-Connect								
	M1,4	Ø 5,0	30°	BHS04004		BHS03176	BHS02063	BHS02064
			20°	BHS04008				
			10°	BHS04012				
	M1,6	Ø 5,0	30°	BHS04004		BHS03176	BHS02065	BHS02066
			20°	BHS04008				
			10°	BHS04012				
	M1,8	Ø 5,0	30°	BHS04004		BHS03176	BHS02067	BHS02068
			20°	BHS04008				
			10°	BHS04012				

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		TI CASTABLE	TI BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OSSTEM IMPLANT COMPATIBLE								
OSSTEM TS								
	Ø 4,0/4,5/5,0/6,0	30º	BHS04059	BHS03100		BHS02009	BHS02010	
		20º	BHS04060					
		10º	BHS04061					
OSSTEM US								
	MINI Ø 3,5	30º	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222
		20º	BHS04060					
		10º	BHS04061					
	RP Ø 4,1	30º	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20º	BHS04019					
		10º	BHS04020					
OSTEOPLUS COMPATIBLE								
OSTEOPLUS - SHe								
	RP Ø 4,1	30º	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20º	BHS04019					
		10º	BHS04020					
OSTEOPLUS - SHi								
	RP Ø 3,5	30º	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20º	BHS04060					
		10º	BHS04061					
	WP Ø 4,5	30º	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
		20º	BHS04019					
		10º	BHS04020					
OXTEIN COMPATIBLE								
HEX EXTERNA								
	RP Ø 4,1	30º	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
		20º	BHS04019					
		10º	BHS04020					
HEX INTERNA								
	NP Ø 3,5	30º	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
		20º	BHS04060					
		10º	BHS04061					
	RP Ø 4,5	30º	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
		20º	BHS04019					
		10º	BHS04020					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		TI CASTABLE	 TI BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
OCTA INTERNA									
	RN	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144
			20°	BHS04022					
			10°	BHS04023					
	RN	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375
			20°	BHS04081					
			10°	BHS04082					
PALTOP COMPATIBLE									
INTERNAL HEX CONNECTION									
	NARROW	Ø 3,25	30°	BHS04059	BHS03198	BHS03199	BHS02041	BHS02042	
			20°	BHS04060					
			10°	BHS04061					
	STANDARD	Ø 3,75 / 4,2 / 5,0	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130
			20°	BHS04060					
			10°	BHS04061					
	WIDE	Ø 6,0	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218
			20°	BHS04019					
			10°	BHS04020					
PCA Implant System (Paltop Conical Active)									
	NARROW	Ø 3,25	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					
	STANDARD	Ø 3,75 / 4,2 / 5,0	30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371
			20°	BHS04060					
			10°	BHS04061					
MULTI UNIT									
	Ø 5,00		30°	BHS04004		BHS03194	BHS02031	BHS02032	
			20°	BHS04008					
			10°	BHS04012					
SINGLE UNIT									
	Ø 4,50		30°	BHS04018	BHS03196		BHS02031	BHS02032	
			20°	BHS04019					
			10°	BHS04020					
PHIBO COMPATIBLE									
TSH									
	S3 / S4	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098
			20°	BHS04019					
			10°	BHS04020					

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA						
RADHEX IMPLANTS COMPATIBLE														
PHE / PHEA														
	NP	Ø 3,5	30°	BHS04059	BHS03010	BHS03011	BHS02005	BHS02006	EM045222					
			20°	BHS04060										
			10°	BHS04061										
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098					
			20°	BHS04019										
			10°	BHS04020										
PHI / PHIA														
	NP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130					
			20°	BHS04060										
			10°	BHS04061										
	RP	Ø 4,5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022	EM046218					
			20°	BHS04019										
			10°	BHS04020										
PCI														
	280		30°	BHS04059	BHS03068	BHS03069	BHS02025	BHS02026	EM046371					
			20°	BHS04060										
			10°	BHS04061										
SGS COMPATIBLE														
P1 / P7														
	RP	Ø 3,75	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130					
			20°	BHS04060										
			10°	BHS04061										
SIGNO VINCES COMPATIBLE														
DUO / INTEGRA														
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098					
			20°	BHS04019										
			10°	BHS04020										
STRAUMANN COMPATIBLE														
OCTA INTERN														
	RN	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144					
			20°	BHS04022										
			10°	BHS04023										
	SYN_RN	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375					
			20°	BHS04081										
			10°	BHS04082										

1.3 CEMENTING TO Ti BASE

	PLATFORM Ø		Ti CASTABLE	 Ti BASE		CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA	
BONE LEVEL									
	NC	Ø 3,3	30°	BHS04025	BHS03064	BHS03065	BHS02015	BHS02016	EM048048
			20°	BHS04060					
			10°	BHS04061					
	RC	Ø 4,1 / 4,8	30°	BHS04025	BHS03036	BHS03037	BHS02015	BHS02016	EM046904
			20°	BHS04060					
			10°	BHS04061					
SWEDEN & MARTINA COMPATIBLE									
SYRA (EH)									
	Ø 4,1 / 4,25	30° 20° 10°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098	
			BHS04019						
			BHS04020						
SYRA Short (EH)									
	Ø 4,1	30° 20° 10°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098	
			BHS04019						
			BHS04020						
TRI DENTAL IMPLANTS COMPATIBLE									
TRI - VENT									
	Ø 3,75 / 4,1 / 4,7	30° 20° 10°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130	
			BHS04060						
			BHS04061						
TRINON COMPATIBLE									
TRINON - Q2									
	RP Ø 4,1	30° 20° 10°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098	
			BHS04019						
			BHS04020						
TRINON - QK									
	RN Ø 4,8	30° 20° 10°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144	
			BHS04022						
			BHS04023						
	RN Ø 4,8	30° 20° 10°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375	
			BHS04081						
			BHS04082						

1.3 CEMENTING TO Ti BASE

				PLATFORM Ø		Ti CASTABLE		Ti BASE 	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ZIACOM COMPATIBLE											
ZINIC - IH											
	NP	Ø 3,2	30°	BHS04059	BHS03198	BHS03199	BHS02041	BHS02042	EM045130		
			20°	BHS04060							
			10°	BHS04061							
	RP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022			
			20°	BHS04060							
			10°	BHS04061							
	WP	Ø4,5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022			
			20°	BHS04019							
			10°	BHS04020							
OEX - EH											
	RP	Ø 4,1	30°	BHS04018	BHS03005	BHS03006	BHS02001	BHS02002	EM045098		
			20°	BHS04019							
			10°	BHS04020							
OIN - HI (M2)											
	RP	Ø 4,0	30°	BHS04059	BHS03048	BHS03049	BHS02011	BHS02012			
			20°	BHS04060							
			10°	BHS04061							
	WP	Ø 5,0	30°	BHS04018	BHS03020	BHS03021	BHS02011	BHS02012			
			20°	BHS04019							
			10°	BHS04020							
XDrive (Multi-unit)											
	RP	Ø 4,8	30°	BHS04004		BHS03009	BHS02003	BHS02004	EM045169		
			20°	BHS04008							
			10°	BHS04012							
ZIMMER COMPATIBLE											
INTERNAL HEX											
	NP	Ø 3,5	30°	BHS04059	BHS03048	BHS03049	BHS02021	BHS02022	EM045130		
			20°	BHS04060							
			10°	BHS04061							
	RP	Ø 4,5	30°	BHS04018	BHS03020	BHS03021	BHS02021	BHS02022			
			20°	BHS04019							
			10°	BHS04020							
SWISS - PLUS											
	RP	Ø 4,8	30°	BHS04021	BHS03016	BHS03017	BHS02011	BHS02012	EM045144		
			20°	BHS04022							
			10°	BHS04023							
	RP	Ø 4,8	30°	BHS04080		BHS03019	BHS02013	BHS02014	EM046375		
			20°	BHS04081							
			10°	BHS04082							

2.0 BHS30 SCREWDRIVER



EXCELLENT MECHANICAL PROPERTIES & SCREW HEAD PROTECTION

Universal contra angle compatible connection (ISO1797).

Hardened Steel offering excellent torsion resistance (90 N-cm).

Optimal coupling engagement ensuring screw head integrity.

TiN Coated: Surface treatment which increases the life of the driver.

Designed to correct 30° angle deviations with 100% efficient torque transmission.

SCREWDRIVER BHS



BHS01001 - L25mm
BHS01004 - L20mm
BHS01005 - L32mm

HANDLER



BHS01110 - Nob Ø 8mm / CA
BHS01111 - Str Ø 7mm / CA

3.0 PROTOCOLS


DIRECT CASTING

Products Required:


BHS30 Screw
BHS30 Screwdriver
Angled Castable Cylinder 10|20|30°

Materials:


Casting in CrCo, CrNi, Ti, Pd, Au,...

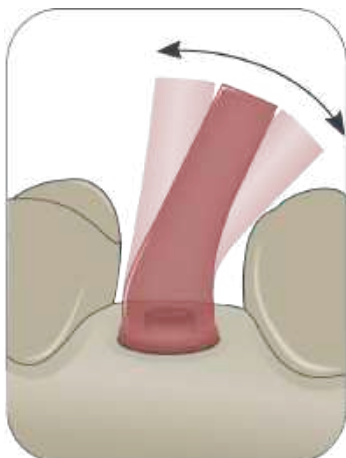


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration.

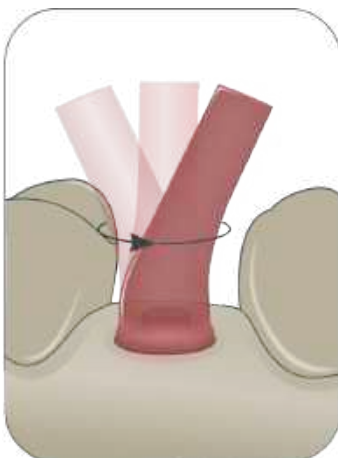


2.
LAB_ Insertion and orientation of the Angled Castable Cylinder on the analogue to obtain an angled access channel with optimum design and emergence profile. Fixation of the abutment with BHS30 Screw.

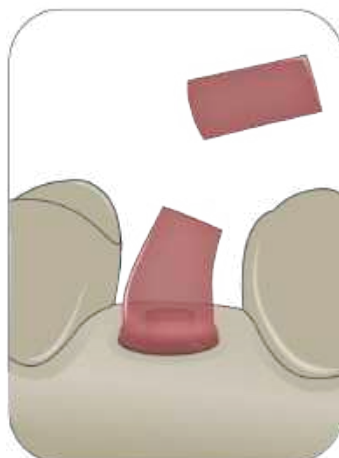


3.
LAB_ Cutting and adjustment of the Angled Castable Cylinder in accordance with the anatomy of the prosthesis and its antagonist.

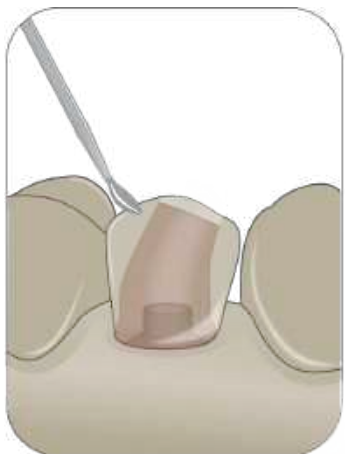

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Cylinder.

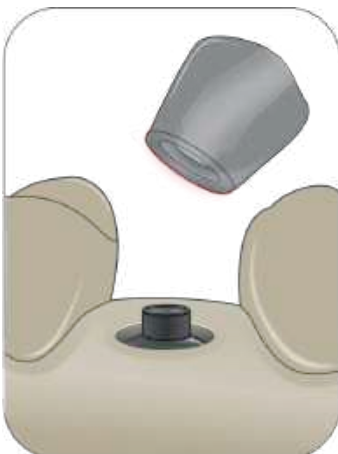


5.
LAB_ Metal casting of the structure and hand finishing of its implant connection to improve fitting on the implant analogue.

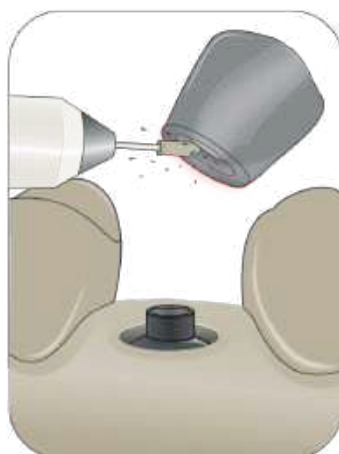

6.
LAB_ Metal casting of the structure and hand finishing of its implant connection to improve fitting on the implant analogue.

3.0 PROTOCOLS

DIRECT CASTING

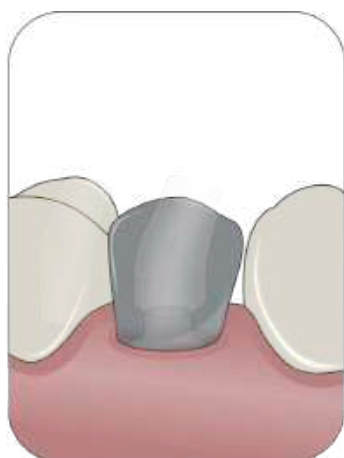

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth.



8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.



9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS


OVERCASTING ON BHS30 CrCo BASE

Products Required:


BHS30 Screw
BHS30 Screwdriver
BHS30 CrCo Base
Angled Castable Chimney 10|20|30°

Materials:


Overcasting in CrCo

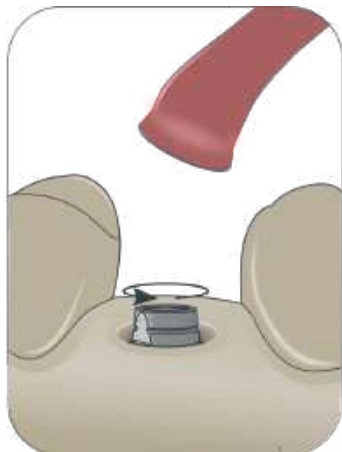

1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 CrCo Base.

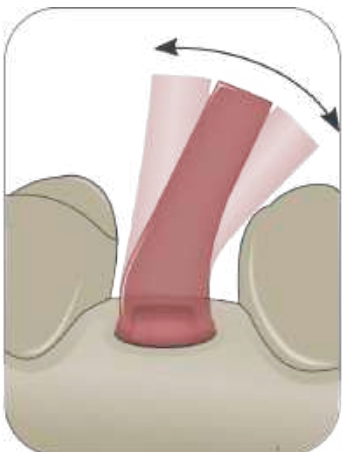

2.
LAB_ Insertion of the Base BHS30 CrCo + Chimney on the implant analogue to obtain an angled access channel with optimal design. Screwing of the BHS30 CrCo Base and final adjustment of the Chimney's orientation according to the desired position for screw access.

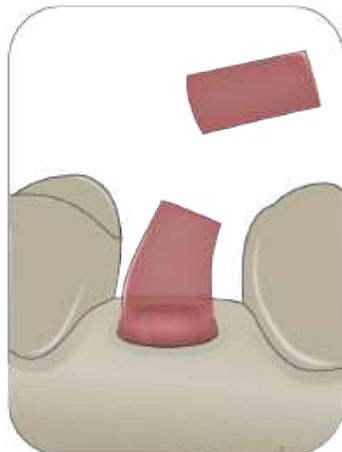


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney in accordance with the anatomy of the prosthesis and its antagonist.



4.
LAB_ Wax-up and modelling of the structure on the BHS30 CrCo BASE + Angled Castable Chimney.




5.
LAB_ Metal overcasting of the structure.

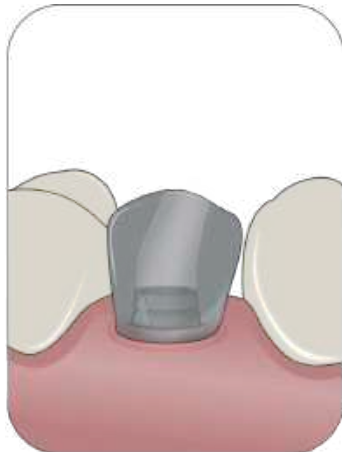


6.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth.

3.0 PROTOCOLS


OVERCASTING ON BHS30 CrCo BASE

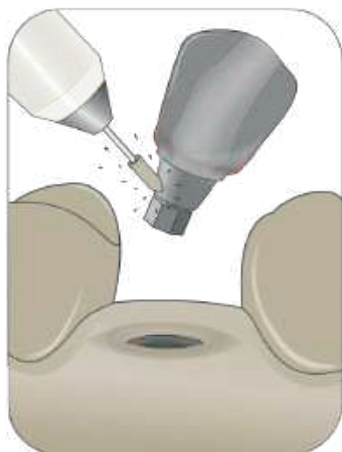

7.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.

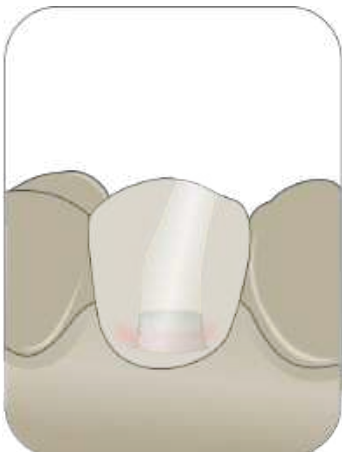

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.

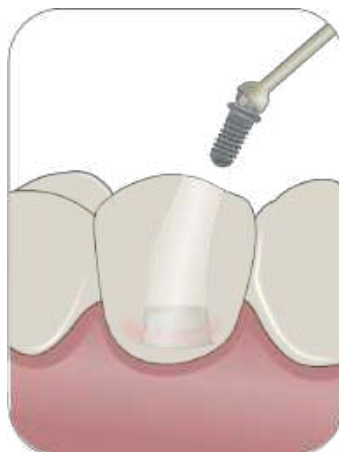

9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS


CASTING + WELDING ON BHS30 CrCo BASE

Products Required:


BHS30 Screw
BHS30 Screwdriver
BHS30 CrCo Base
Angled Castable Chimney 10|20|30°

Materials:


Casting in CrCo

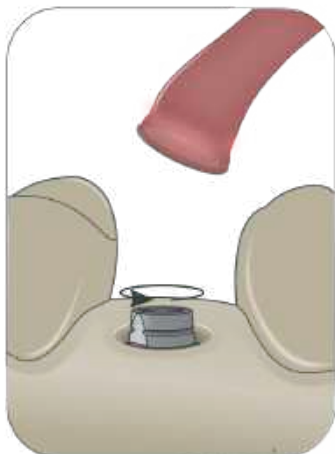


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 CrCo Base.

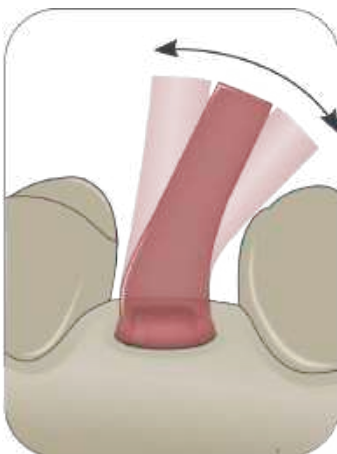

2.
LAB_ Insertion of the Base BHS30 CrCo + Chimney on the implant analogue to obtain an angled access channel with optimal design. Screwing of the BHS30 CrCo Base and final adjustment of the Chimney's orientation according to the desired position for screw access.

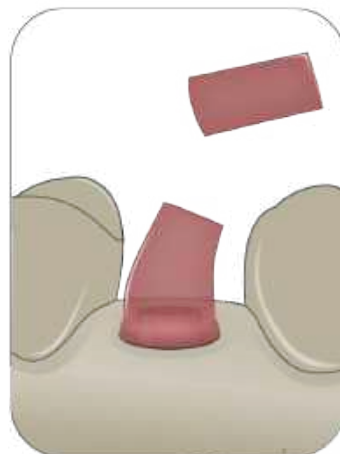


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney in accordance with the anatomy of the prosthesis and its antagonist.

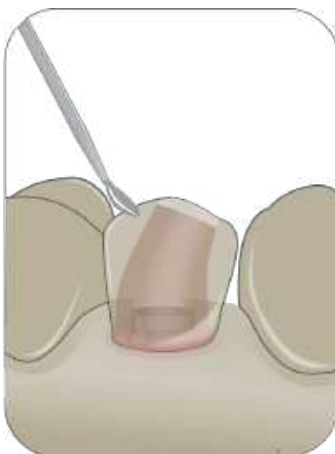

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Chimney.




5.
LAB_ Removal of the Angled Castable Chimney with the modelling of the structure for metal casting, leaving the BHS30 CrCo Base screwed onto the model.

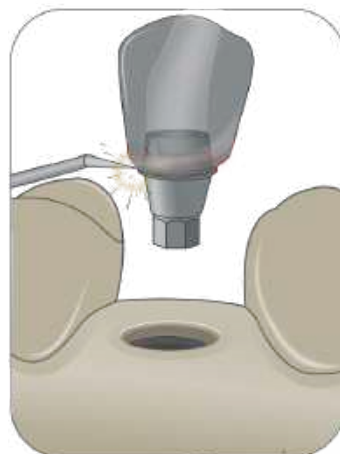


6.
A. LAB_ Adjustment of the metal cast structure on the BHS30 CrCo Base until obtaining passive fit on the model and perimetral laser welding of both parts and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS


CASTING + WELDING ON BHS30 CrCo BASE

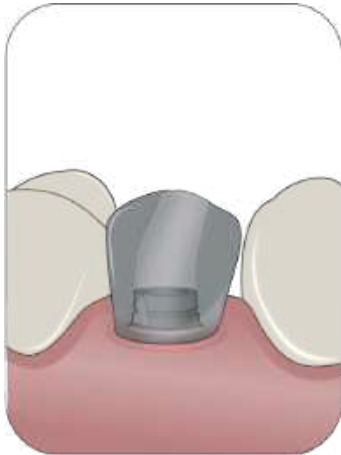

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.

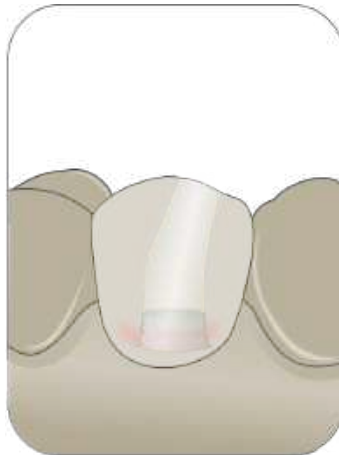

8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and perimetral laser welding of both parts. Aesthetic restoration, polishing and finishing.

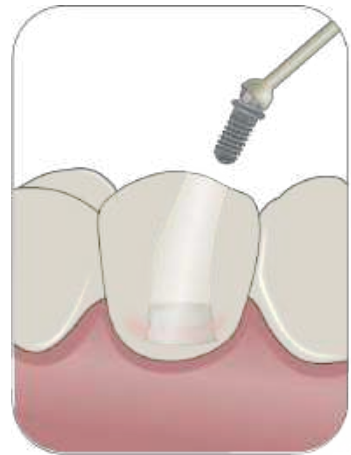

9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS


CASTING + CEMENTING ON BHS30 Ti BASE

Products Required:


BHS30 Screw
BHS30 Screwdriver
BHS30 Ti Base
Angled Castable Chimney 10|20|30°

Materials:


Casting in metal (Ti, CrCo, CrNi, noble...)

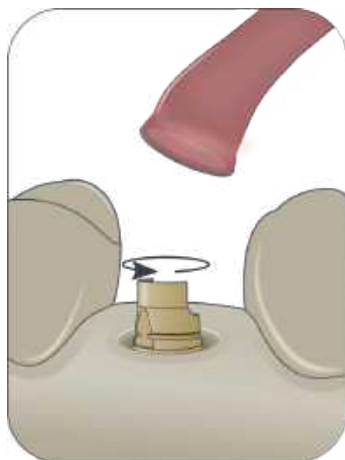


1.
LAB_ Check implant analogue position and selection of the appropriate Angled Castable Chimney Cylinder according to the undesired emergence on the restoration. Coupling of the Chimney to the axial retention guide in the BHS30 Ti Base.

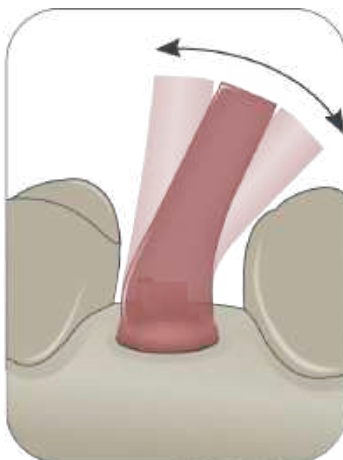


2.
LAB_ Insertion and orientation of the BHS30 Ti Base + Chimney on the implant analogue to obtain an angled access channel with optimal design and emergence profile. Screwing of the BHS30 Ti Base.

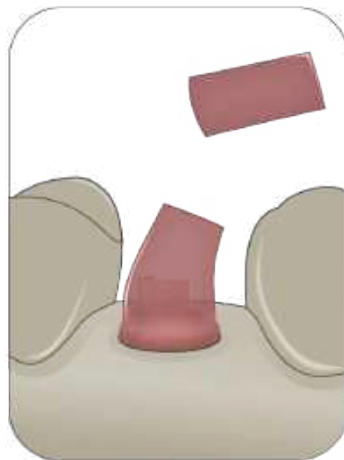


3.
LAB_ Cutting and adjustment of the Angulated Castable Chimney according to anatomy of the prosthesis and its antagonist.

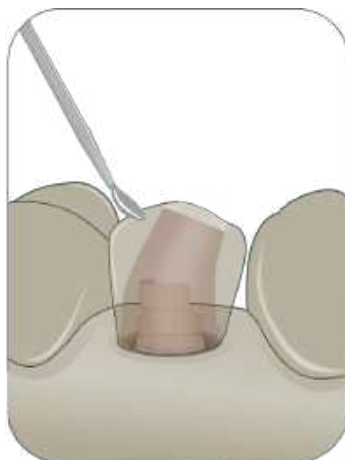

4.
LAB_ Wax-up and modelling of the structure on the Angled Castable Chimney.

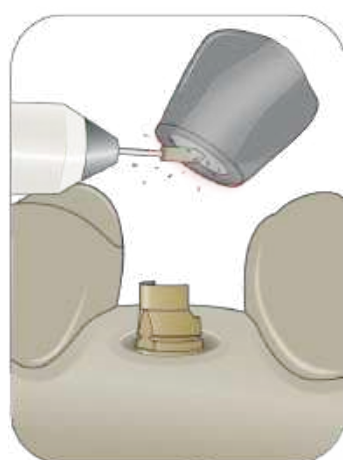


5.
LAB_ Removal of the Angled Castable Chimney with the modelling of the structure for metal casting, leaving the BHS30 Ti Base screwed onto the model.

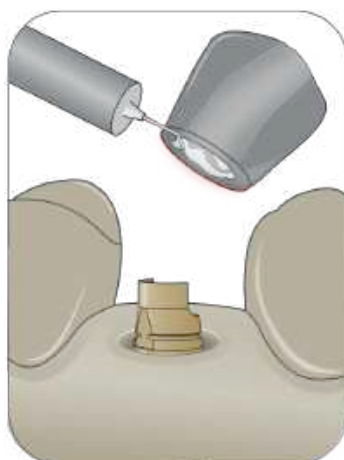

6.
A. LAB_ Adjustment of the metal cast structure on the BHS30 Ti Base until obtaining passive fit on the model and cementing of both parts, the excess of cement being expelled through the relining holes in the Chimney and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS

CASTING + CEMENTING ON BHS30 Ti BASE

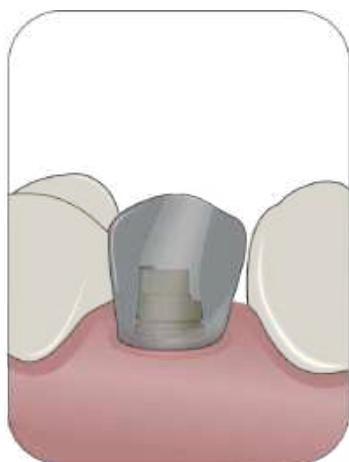

7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.



8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and cementing of both parts, the excess of cement being expelled through the relining holes in the Chimney. Aesthetic restoration, polishing and finishing.



9.
CLINIC_ Placement of the prosthesis to the patient.

4.0 COMPATIBILITY CHART

TRADEMARK	PAGES			TRADEMARK	PAGES		
	DIRECT CASTING	OVERCASTING CrCo BASE	CEMENTING Ti BASE		DIRECT CASTING	OVERCASTING CrCo BASE	CEMENTING Ti BASE
AB IMPLANTS	6	37	62	IMPLANT DIRECT	22	50	77
ACE	6	37	62	JDENTAL CARE	23	51	78
ADIN	6	37	63	KEYSTONE DENTAL	23	51	78
ALPHABIO	7	38	63	KLOCKNER	24	51	79
ALPHADENT	7	38	64	LEADER	24	51	79
ANCLADEN	8	38	64	MEDENTIS	24	52	70
ANTHOGRYR	8	39	64	MEGAGEN	25	52	80
ARDS	8	39	65	MICRODENT	25	53	80
ASTRA	9	39	65	MIS	26	53	81
AVINENT	9	40	66	MOZO GRAU	27	54	81
BEGO	10	40	66	MPI	27	54	82
BIOGENESIS	10	41	67	NEOBIOTECH	28	54	82
BIOHORIZONS	11	41	67	NEODENT	28	55	82
BIOMET 3i	12	41	68	NOBEL	28	55	83
BIONER	12	42	68	NT TRADING	29	56	83
BTI	13	42	68	OSSTEM IMPLANT	29	56	84
BTK	13	42	69	OSTEOPLUS	30	56	84
B&W	14	43	70	OXTEIN	30	57	84
CAMLOG	15	44	70	PALTOP	31	57	85
CON. SISTEMA PROTESE	15	44	70	PHIBO	31	58	85
CORTEX	16	45	71	RADHEX IMPLANTS	32	58	86
COWELLMEDI	16	45	71	SGS	32	58	86
DENTALTECH	17	45	72	SIGNO VINCES	32	59	86
DENTIS	17	46	72	STRAUMANN	33		86
DENTIUM	17	46	72	SWEDEN & MARTINA	33	59	87
DENTSPLY SIRONA	18	46	73	TRI DENTAL IMPLANTS	34	59	87
DIO IMPLANTS	18	47	73	TRINON	34	59	87
EDISON MEDICAL	18	47	74	ZIACOM	34	60	88
ELITE MEDICA	19	47	74	ZIMMER	35	60	88
EUROTEKNIKA	19	48	74				
GALIMPLANT	20	48	75				
GMI	20	49	76				
GT MEDICAL	21	49	76				
HI TECH	22	50	77				



www.ballheadtechnology.com

(+34) 973 831 391 · info@ballheadtechnology.com