

Ball Head Technology®

Pioneers in Angulation



2020 PRODUCT CATALOGUE
CAD/CAM PROSTHESIS



30° OF FREEDOM



The evolution towards customization in dental prosthetic design increasingly requires the use of angled abutments to provide versatility, simplicity and ease of use. Our Premium and Universal Solutions are pioneers in covering these needs.



Our main goal is achieving excellence to generate safety, reliability and comfort when using our products, both for the dental lab and clinical professionals, responsible for the patient's treatment.



Since our first patent in 2008 we have developed the most efficient and innovative angled solutions based on the vision and experience of the most qualified restoration and oral surgery professionals.

INDEX

1.0 CAD/CAM PROSTHESIS	5/88
1.1 CAD/CAM IMPLANT LEVEL	5/32
1.2 CAD/CAM WELDING TO CrCo BASE	33/60
1.3 CAD/CAM CEMENTING TO Ti BASE	61/88
 2.0 BHS30 SCREWDRIVER	 89
 3.0 PROTOCOLS	 90/95
 4.0 COMPATIBILITY CHART	 96

Av. de Les Garrigues 17
25001 Lleida · Spain
(+34) 973 831 391
info@ballheadtechnology.com
www.ballheadtechnology.com

All names and trademarks used in this catalogue have been registered and are propriety of their respective companies.
Images published are merely intended for commercial purposes and have no contractual nature.
The information contained in this catalogue is aimed exclusively to professionals in the dental sector.
Ball Head Technology® reserves the right to make changes to this document without prior notice.
The guarantee conditions of Ball Head Technology® will apply as long as the instructions for each reference are respected.



1.1 CAD/CAM IMPLANT LEVEL



ABSOLUTE FREEDOM AND ACCESSIBILITY

Wide flat wall for an optimum Best Fit (Position and Orientation).

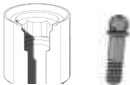
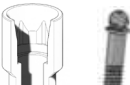
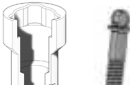
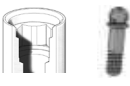
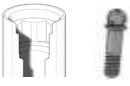
Manufactured in Titanium Grade 5 with tolerances non exceeding $\pm 5\mu\text{m}$.

Optimal surface finishing for a reliable scanning.


Exact connections for all major Implant Systems.

No need for spray.


1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE									
I SERIES INTERNAL HEX									
	Narrow Ø 3 - Ø 3,75	BHS08951	BHS30-9500 (NARROW)				BHS02041	BHS02042	
	STD Ø 3,5 - 6,0	BHS08321	BHS30-3200 (NP)				BHS02021	BHS02022	EM045130
I 10C CONUS									
	RP Ø 4,2	BHS08312	BHS30 - 3100 (RP)				BHS02023	BHS02024	EM048149
ACE COMPATIBLE									
INFINITY EXTERNAL HEX									
	NP Ø 3,3 / 3,75	BHS08611	BHS30 - 6100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,0	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INFINITY TRI - CAM									
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)				BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)				BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
INFINITY INTERNAL HEX									
	NP Ø 3,7 / 4,1	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,7 / 5,1	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
ADIN COMPATIBLE									
SWELL									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
TOUAREG S / OS									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130

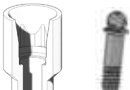1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOUAREG CLOSEFIT									
	RP Ø 3,5	BHS08311	BHS30 - 3100 (NP)				BHS02025	BHS02026	EM048148
	WP Ø 4,3 / 5	BHS08312	BHS30 - 3100 (RP)				BHS02023	BHS02024	EM048149
ALPHABIO COMPATIBLE									
EXTERNAL ALPHABIO									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
INTERNAL ALPHABIO									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
ALPHA DENT COMPATIBLE									
ACTIVE									
	Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
ACTIVE PLUS									
	Ø 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 / 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
ANCLADEN COMPATIBLE									
ANCLALOCK									
	Ø 3,75 - 5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130

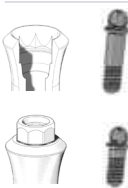1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANTHOGRYR COMPATIBLE									
ANTHOFIT EH									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
OSSFIT									
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
ARDS COMPATIBLE									
SMART									
	Ø 3,75 - 4,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
PREMIUM									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
ASTRA COMPATIBLE									
AQUA									
	RP_3.5 / 4.0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
LILAC									
	WP_4.5 / 5.0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
CONE 20°									
	Ø 3,5	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CONE 45°									
 EV	Ø 3,5	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786
 EV	Ø3,6	BHS08152	BHS30 - 1500 (Ø3,6)				BHS02053	BHS02054	
	Ø4,2	BHS08153	BHS30 - 1500 (Ø4,2)				BHS02055	BHS02056	
AVINENT COMPATIBLE									
EH									
 EH	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02031	BHS02032	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
 IH	Ø 3,5 / 4,1	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
MULTIUNIT									
 MULTIUNIT	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02003	BHS02004	EM045169
OCEAN CONICAL (C.C.)									
 OCEAN CONICAL (C.C.)	Ø 3,5 (Implant Ø 3,5/Ø 4,0)	BHS08311	BHS30 - 3100 (NP)				BHS02025	BHS02026	EM048148
	Ø 4,1 (Implant Ø 4,5/Ø 5,0)	BHS08312	BHS30 - 3100 (RP)				BHS02023	BHS02024	EM048149
BEGO COMPATIBLE									
EH									
 BEGO COMPATIBLE EH	NP_Ø 2,7 / 3,1	BHS08421	BHS30 - 4200 (NP)				BHS02013	BHS02014	EM047957


1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
IH									
	Ø 3,75 / 4,1	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
BIOGENESIS COMPATIBLE									
3ICON									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
ITICON									
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
ATICON									
	RP Ø 3,5 / 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
	WP Ø 4,5 / 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
KISSES									
	RP Ø 3,5 / 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
	WP Ø 4,5 / 5,0 / 5,5	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
CONE 20°									
	RP / WP	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638
CONE 45°									
	RP / WP	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BIOHORIZONS COMPATIBLE									
INTERNAL									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
BIOMET 3i COMPATIBLE									
EH OSSEOTITE									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
CERTAIN									
	NP Ø 3,4	BHS08511	BHS30 - 5100 (NP)				BHS02007	BHS02008	EM045166
	RP Ø 4,1	BHS08512	BHS30 - 5100 (RP)				BHS02007	BHS02008	EM045165
	WP Ø 5,0	BHS08513	BHS30 - 5100 (WP)				BHS02007	BHS02008	EM045161
BIONER COMPATIBLE									
IKELT / BIKELT									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
HIKELT									
	NP Ø 3,8	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,7	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218

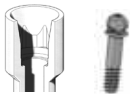
1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOP DM									
	Ø 3,5 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)				BHS02025	BHS02026	EM048148
MULTIUNIT									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02003	BHS02004	EM045169
BTI COMPATIBLE									
EH									
	NP Tiny Ø 3,5	BHS08421	BHS30 - 4200 (NP)				BHS02013	BHS02014	EM047957
	RP Ø 4,1	BHS08422	BHS30 - 4200 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,5	BHS08423	BHS30 - 4200 (WP)				BHS02013	BHS02014	EM047967
CUATRILOBULAR									
	RP Ø 4,1	BHS08231	BHS30 - 2300 (RP)				BHS02021	BHS02022	EM047743
	WP Ø 5,5	BHS08232	BHS30 - 2300 (WP)				BHS02021	BHS02022	
BTK COMPATIBLE									
8NECK INT									
	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
EVO DL EXT									
	NP Ø 3,4 / 3,7	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1 / 4,6	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0 / 5,6	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
EVO DL INT									
	Ø 3,5 / 3,7 / 4,6 / 5,0 / 5,6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLASSIC EXT									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KLASSIC INT									
	Ø 3,5/ 4,5 / 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
KONIK EXT									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KONIC INT									
	Ø 3,5/ 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
SAFETL INT									
	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
B&W COMPATIBLE									
HEXAGONO EXT									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082

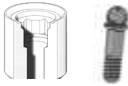1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
HEXAGONO INT									
	NP Ø 3,3	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	WP Ø 5,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
CAMLOG COMPATIBLE									
CAMLOG									
	Ø 3,3	BHS08721	BHS30-7200 (Ø3,3)				BHS02049	BHS02050	
	Ø 3,8	BHS08722	BHS30-7200 (Ø3,8)				BHS02049	BHS02050	
	Ø 4,3	BHS08723	BHS30-7200 (Ø4,3)				BHS02049	BHS02050	
	Ø 5,0	BHS08724	BHS30-7200 (Ø5,0)				BHS02051	BHS02052	
CONEXÃO SIST. PROTESE COMPATIBLE									
MORSE INDEXADO									
	Ø 3,5 / 3,75 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)				BHS02025	BHS02026	EM048148
HEXAGONO EXT									
	Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	Ø 3,75 / 4,0	BHS08742	BHS30 - 7400 (S3 y S4)				BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEXAGONO INT									
	Ø 3,75 / 4,0 / 5,0	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218


1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CORTEX COMPATIBLE									
DYNAMIX									
 	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
CLASSIX									
 	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
SATURN									
 	Ø 3,8 / 4,2	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
COWELLMEDI COMPATIBLE									
EXTERNAL									
 	RP Ø 4,1	BHS08212	BHS30 - 2600 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2600 (WP)				BHS02039	BHS02040	EM045070
INTERNAL									
 	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
SUBMERGED (INT)									
 	Ø 3,5 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)				BHS02009	BHS02010	
DENTALTECH COMPATIBLE									
IMPLOGIC GII									
 	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
IMPLASSIC HX2									
 	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DENTIS COMPATIBLE								
S-CLEAN								
 	RP Ø 3,4	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
I-CLEAN								
 	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
DENTUM COMPATIBLE								
SIMPLE LINE II								
 	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
IMPLANTUM								
 	Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
SUPER LINE								
 	Ø 3,6 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
DENTSPLY SIRONA COMPATIBLE								
XIVE S								
	D 3.4	BHS08221	BHS30-2200 (D3.4)			BHS02047	BHS02048	
	D 3.8	BHS08222	BHS30-2200 (D3.8)			BHS02047	BHS02048	
	D 4.5	BHS08223	BHS30-2200 (D4.5)			BHS02047	BHS02048	

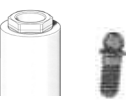1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DIO IMPLANTS COMPATIBLE									
EXTERNAL									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
UF II									
	Ø 3,8 / 4 / 4,5 / 5 / 5,5	BHS08432	BHS30 - 4600 (4.0)				BHS02009	BHS02010	
EDISON MEDICAL COMPATIBLE									
LAMINA (IH)									
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
FORTIS (IH)									
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
ELITE MEDICA COMPATIBLE									
EXTERNAL CONNECT									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
EUROTEKNIKA COMPATIBLE									
UNEVA									
	Ø 3,6 / 4,1 / 4,8	BHS08212	BH30 - 2100 (RP)				BHS02001	BHS02002	EM045098
NATURACTIS									
	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BH30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NATURALL +									
	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
NATEA									
	NP Ø 3,6 / 4,1	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
	RP Ø 4,8 / 6	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
AESTHETICA									
	Ø 3,6 / 4,1 / 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
GALIMPLANT COMPATIBLE									
IPXE - EXTERNA									
	RP Ø 3,5 / 4 / 5	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
ICI - INTERNA									
	RP Ø 3,2 / 3,5 / 4 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
IPX - INTERNA									
	RP Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
GMI COMPATIBLE									
FRONTIER (IH)									
	RP Standard Ø 3,3	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	WP Wide Ø 4,3	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
FRONTIER multi-aesthetic									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02003	BHS02004	EM045169

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHOENIX (EH)									
	Standard Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	Wide (Hex. 2,7) Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	Wide (Hex. 3,4) Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
PHOENIX multi-aesthetic									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02003	BHS02004	EM045169
GT MEDICAL COMPATIBLE									
BEST FIT HXE									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
BEST FIT HXI									
	WP 4,5 / 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
BEST FIT 002									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
BEST FIT OCT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
HI TECH COMPATIBLE									
SELF THREAD									
	SLT Ø 3,3 / 3,75	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	GTI Ø 4,2 / 5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218



1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LGI PLUS									
	Ø 3,5 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
IMPLANT DIRECT COMPATIBLE									
INTERACTIVE									
	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)				BHS02025	BHS02026	EM048148
	RP Ø 4,3	BHS08312	BHS30 - 3100 (RP)				BHS02023	BHS02024	EM048149
LEGACY									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
REPLUS / REPLANT									
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)				BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)				BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
SWISHPLUS / SWISHPLANT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
JDENTALCARE COMPATIBLE									
JDEVOLUTION									
	Ø 3,7 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130


1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KEYSTONE DENTAL COMPATIBLE									
RESTORE									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)				BHS02001	BHS02002	EM045098
OCT-XT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
RENOVA									
	NP Ø 3,75	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5 / 4,75	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
KLOCKNER COMPATIBLE									
SK2 (EH)									
	RP Ø 4,2	BHS08521	BHS30 - 5200 (RP)				BHS02031	BHS02032	EM046747
KL (EH)									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070

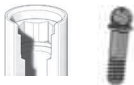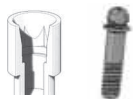1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LEADER COMPATIBLE								
TIXOS INTERNAL HEX								
	NP Ø 3,5 / 4,0	BHS08321	BHS30 - 3200 (NP)			BHS02021	BHS02022	EM045130
TIXOS EXTERNAL HEX								
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)			BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082
MEDENTIS COMPATIBLE								
ICX								
	Ø 3,45 / 3,7 - 4,1 / 4,8	BHS08751	BHS30-7700 (Ø 3,45 / 3,7 - 4,1 / 4,8)			BHS02015	BHS02016	
MEGAGEN COMPATIBLE								
ANYONE INTERNAL								
	NP	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
	RP	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
	WP	BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
ANYONE EXTERNAL								
	NP	BHS08211	BHS30 - 2100 (NP)			BHS02005	BHS02006	EM045222
	RP	BHS08212	BHS30 - 2100 (RP)			BHS02001	BHS02002	EM045098
	WP	BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082

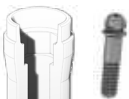1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANYRIDGE									
	Small Ø3,5	BHS08991	BHS30-9900 (4,5/5,0)				BHS02027	BHS02028	
	Regular Ø4,0/4,5	BHS08991	BHS30-9900 (4,5/5,0)				BHS02027	BHS02028	
	Wide Ø5,0/5,5	BHS08991	BHS30-9900 (4,5/5,0)				BHS02027	BHS02028	
OCTA Abutment									
	SYN_RN Ø4,8	BHS08713	BHS30-7100 (RN SYNOCTA Ø4,8)				BHS02013	BHS02014	EM046375
MICRODENT COMPATIBLE									
EKTOS (IH)									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
MINI (EH)									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
SUV (EH)									
	SUV Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	SUV Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYSTEM (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
UNIVERSAL (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	SU Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	SU Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MIS COMPATIBLE									
LANCE									
	STD Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SEVEN / M4									
	Narrow Ø 3,30	BHS08951	BHS30-9500 (NARROW)				BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	Wide Ø 5,0 / 6,0	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
MULTIUNIT MIS									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02037	BHS02038	
MOZO GRAU COMPATIBLE									
INHEX									
	RP Ø 3,75 / 4,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
	CONO 20	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638
	CONO 45	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OSSEUS									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
MPI COMPATIBLE									
EXTERNA EH									
	NP Ø 3,5	BHS08421	BHS30 - 4200 (NP)				BHS02001	BHS02002	EM047957
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INTERNA CM									
	RP Ø 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)				BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
NEOBIOTECH COMPATIBLE									
EB EXTERNAL SYSTEM									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
IS INTERNAL SYSTEM									
	Ø 3,5 / 4 / 4,5 / 5	BHS08432	BHS30 - 4600 (4.0)				BHS02009	BHS02010	



1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NEODENT COMPATIBLE								
TITAMAX HE SMART								
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)			BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082
MINI PILAR CM								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS02003	BHS02004	EM045169
NOBEL COMPATIBLE								
BRÄNEMARK								
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)			BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)			BHS02039	BHS02040	EM045070
REPLACE								
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)			BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)			BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)			BHS02009	BHS02010	EM046293
ACTIVE								
	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)			BHS02025	BHS02026	EM048148
	RP Ø 4,3 / Ø 5,0	BHS08312	BHS30 - 3100 (RP)			BHS02023	BHS02024	EM048149
MULTIUNIT								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS02003	BHS02004	EM045169

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NT TRADING COMPATIBLE									
2-CONnect									
	M1,4	Ø5,0	BHS08771	BHS30-7700 (M1,4)			BHS02063	BHS02064	
	M1,6	Ø5,0	BHS08771	BHS30-7700 (M1,6)			BHS02065	BHS02066	
	M1,8	Ø5,0	BHS08771	BHS30-7700 (M1,8)			BHS02067	BHS02068	
OSSTEM IMPLANT COMPATIBLE									
OSSTEM TS									
	Ø 4 / 4,5 / 5 / 6		BHS08432	BHS30 - 4600 (4.0)			BHS02009	BHS02010	
OSSTEM US									
	MINI Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
OSTEOPLUS COMPATIBLE									
ARGENTUM									
	RP Ø 4,2	BHS08521	BHS30 - 5200 (RP)				BHS02031	BHS02032	EM046747
OSTEOPLUS - SHe									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
OSTEOPLUS - SHi									
	RP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	WP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218

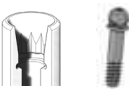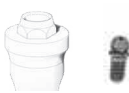1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OXTEIN COMPATIBLE									
HEX EXTERNA									
	MINI Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)				BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEX INTERNA									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
OCTA INTERNA									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
PALTOP COMPATIBLE									
INTERNAL HEX CONNECTION									
	Narrow Ø 3,25	BHS08951	BHS30-9500 (NARROW)				BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30-3200 (NP)				BHS02021	BHS02022	EM045130
	Wide Ø 6,0	BHS08322	BHS30-3200 (RP)				BHS02021	BHS02022	EM046218
PCA Implant System (Paltop Conical Active)									
	Narrow Ø 3,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)				BHS02025	BHS02026	EM046371
	Standard Ø 3,75 / 4,2 / 5,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)				BHS02025	BHS02026	EM046371
MULTI UNIT									
	Ø 5,0	BHS08271	BHS30-2700 (MU)				BHS02031	BHS02032	
SINGLE UNIT									
	Ø 4,5	BHS08251	BHS30-2500 (SINGLE)				BHS02031	BHS02032	

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHIBO COMPATIBLE									
TSH									
	S2 Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	S3 / S4 Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)				BHS02001	BHS02002	EM045098
AUREA EVO									
	NP Ø 3,5						BHS02041	BHS02042	
	RP Ø 4,3						BHS02021	BHS02022	
RADHEX IMPLANTS COMPATIBLE									
PHE / PHEA									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
PHI / PHIA									
	NP Ø3,5	BHS08321	BHS30-3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø4,5	BHS08322	BHS30-3200 (RP)				BHS02021	BHS02022	EM046218
PCI									
	280	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)				BHS02025	BHS02026	EM046371
	350	BHS08123	BHS30 - 1200 (WP Ø 4.5 /5.0)				BHS02009	BHS02010	EM046370
SGS COMPATIBLE									
P1 / P7									
	RP Ø 3,75	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SIGNO VINCES COMPATIBLE								
DUO / INTEGRA								
	RP Ø 4,1	BHS08212	BHS30-2100 (RP)			BHS02001	BHS02002	EM045098
DUOCON / CM								
	WP_4.5 / 5.0	BHS08123	BHS30-1200 (WP Ø 4.5 / 5.0)			BHS02009	BHS02010	EM046370
STRAUMANN COMPATIBLE								
OCTA INTERN								
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)			BHS02013	BHS02014	EM046375
BONE LEVEL								
	NC Ø 3,3	BHS08811	BHS30 - 8100 (NC)			BHS02015	BHS02016	EM048048
	RC Ø 4,1/ Ø 4,8	BHS08812	BHS30 - 8100 (RC)			BHS02015	BHS02016	EM046904
BL MULTIBASE								
	RC Ø 4,8					BHS02017	BHS02018	EM047780

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SWEDEN & MARTINA COMPATIBLE									
PREMIUM KOHNO									
	NP Ø 3,3	BHS08961	BHS30 - 9600 (NP Ø 3,3)				BHS02021	BHS02022	
	RP Ø 3,8	BHS08962	BHS30 - 9600 (RP Ø 3,8)				BHS02021	BHS02022	
	WP Ø 4,25	BHS08963	BHS30 - 9600 (WP Ø 4,25)				BHS02011	BHS02012	
SYRA (EH)									
	Ø 4,10 / 4,25	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYRA Short (EH)									
	Ø 4,10	BHS08612	BHS30 - 6100 (RP)				BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
TRI DENTAL IMPLANTS COMPATIBLE									
TRI - VENT									
	Ø 3,75 / 4,1 / 4,7	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
TRINON COMPATIBLE									
TRINON - Q2									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
TRINON - QK									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCCTA Ø 4,8)				BHS02013	BHS02014	EM046375

1.1 CAD/CAM IMPLANT LEVEL

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ZIACOM COMPATIBLE									
ZINIC - IH									
	NP Ø 3,2	BHS08951	BHS30-9500 (NARROW)				BHS02041	BHS02042	
	RP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	WP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
OEX - EH									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)				BHS02001	BHS02002	EM045098
OIN - IH (M2)									
	RP Ø 4,0	BHS08321	BHS30-3200 (NP)				BHS02011	BHS02012	
	WP Ø 5,0	BHS08322	BHS30-3200 (RP)				BHS02011	BHS02012	
XDrive (Multi-unit)									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)				BHS02003	BHS02004	EM045169
ZIMMER COMPATIBLE									
INTERNAL HEX									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)				BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)				BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
SWISS - PLUS									
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375

1.2 CAD/CAM WELDING TO CrCo BASE



AFFORDABLE STRUCTURES WITH PERFECT MECHANIZED SEATING

CrCo Base:

Optimum height allowing 30° angulation for BHS screwdriver.

Conical design for multiple restorations and engaging flat wall for single ones.

Scanbody:

Wide flat wall for an optimum Best Fit (Position and Orientation).

Manufactured in Titanium Grade 5 with tolerances non exceeding $\pm 5\mu\text{m}$.

Optimal surface finishing for a reliable scanning.

Exact connections for all major Implant Systems.

No need for spray.


1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE									
I SERIES INTERNAL HEX									
	Narrow Ø 3 - Ø 3,75	BHS08951	BHS30-9500 (NARROW)	 	BHS03200 BHS03201		BHS02041	BHS02042	
	STD Ø 3,5 - 6,0	BHS08321	BHS30-3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
I 10C CONUS									
	RP Ø 4,2	BHS08312	BHS30 - 3100 (RP)	 	BHS03046 BHS03047		BHS02023	BHS02024	EM048149
ACE COMPATIBLE									
INFINITY EXTERNAL HEX									
	NP Ø 3,3 / 3,75	BHS08611	BHS30 - 6100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,0	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INFINITY TRI - CAM									
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)				BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)				BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
INFINITY INTERNAL HEX									
	NP Ø 3,7 / 4,1	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,7 / 5,1	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
ADIN COMPATIBLE									
SWELL									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
TOUAREG S / OS									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOUAREG CLOSEFIT									
	RP Ø 3,5	BHS08311	BHS30 - 3100 (NP)	 	BHS03058 BHS03059		BHS02025	BHS02026	EM048148
	WP Ø 4,3 / 5	BHS08312	BHS30 - 3100 (RP)	 	BHS03046 BHS03047		BHS02023	BHS02024	EM048149
ALPHABIO COMPATIBLE									
EXTERNAL ALPHABIO									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
INTERNAL ALPHABIO									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
ALPHA DENT COMPATIBLE									
ACTIVE									
	Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
ACTIVE PLUS									
	Ø 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 / 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
ANCLADEN COMPATIBLE									
ANCLALOCK									
	Ø 3,75 - 5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130

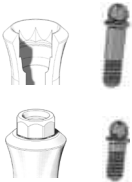1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANTHOGRYR COMPATIBLE									
ANTHOFIT EH									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
OSSFIT									
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
ARDS COMPATIBLE									
SMART									
	Ø 3,75 - 4,5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
PREMIUM									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
ASTRA COMPATIBLE									
AQUA									
	RP_3.5 / 4.0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
LILAC									
	WP_4.5 / 5.0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
CONE 20°									
	Ø 3,5	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638


1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CONE 45°									
	Ø 3,5	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786
EV									
	Ø3,6	BHS08152	BHS30 - 1500 (Ø3,6)	 	BHS03154 BHS03155		BHS02053	BHS02054	
	Ø4,2	BHS08153	BHS30 - 1500 (Ø4,2)	 	BHS03158 BHS03159		BHS02055	BHS02056	
AVINENT COMPATIBLE									
EH									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02031	BHS02032	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
IH									
	Ø 3,5 / 4,1	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
MULTIUNIT									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)		BHS03008		BHS02003	BHS02004	EM045169
OCEAN CONICAL (C.C.)									
	Ø 3,5 (Implant Ø 3,5/Ø 4,0)	BHS08311	BHS30 - 3100 (NP)	 	BHS03058 BHS03059		BHS02025	BHS02026	EM048148
	Ø 4,1 (Implant Ø 4,5/Ø 5,0)	BHS08312	BHS30 - 3100 (RP)	 	BHS03046 BHS03047		BHS02023	BHS02024	EM048149
BEGO COMPATIBLE									
EH									
	NP_Ø 2,7 / 3,1	BHS08421	BHS30 - 4200 (NP)				BHS02013	BHS02014	EM047957

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
IH									
  	Ø 3,75 / 4,1	BHS08321	BHS30 - 3200 (NP)	 BHS03050  BHS03051			BHS02021	BHS02022	EM045130
BIOGENESIS COMPATIBLE									
3ICON									
  	NP Ø 3,5 RP Ø 4,1 WP Ø 5,0	BHS08611 BHS08612 BHS08613	BHS30 - 6100 (NP) BHS30 - 6100 (RP) BHS30 - 6100 (WP)	 BHS03003  BHS03004			BHS02001 BHS02001 BHS02001	BHS02002 BHS02002 BHS02002	EM045468 EM045098 EM045082
ITICON									
   	RP Ø 4,8 RP Ø 4,8	BHS08712 BHS08713	BHS30 - 7100 (RN Ø 4,8) BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02011 BHS02013	BHS02012 BHS02014	EM045144 EM046375
ATICON									
  	RP Ø 3,5 / 4,0 WP Ø 4,5 / 5,0	BHS08122 BHS08123	BHS30 - 1200 (RP Ø 3.5 / 4.0) BHS30 - 1200 (WP Ø 4.5 / 5.0)	 BHS03070  BHS03071			BHS02025 BHS02009	BHS02026 BHS02010	EM046371 EM046370
KISSES									
  	RP Ø 3,5 / 4,0 WP Ø 4,5 / 5,0 / 5,5	BHS08122 BHS08123	BHS30 - 1200 (RP Ø 3.5 / 4.0) BHS30 - 1200 (WP Ø 4.5 / 5.0)	 BHS03070  BHS03071			BHS02025 BHS02009	BHS02026 BHS02010	EM046371 EM046370
CONE 20°									
 	RP / WP	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638
CONE 45°									
 	RP / WP	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786

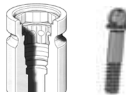1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BIOHORIZONS COMPATIBLE									
INTERNAL									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
BIOMET 3i COMPATIBLE									
EH OSSEOTITE									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
CERTAIN									
	NP Ø 3,4	BHS08511	BHS30 - 5100 (NP)				BHS02007	BHS02008	EM045166
	RP Ø 4,1	BHS08512	BHS30 - 5100 (RP)				BHS02007	BHS02008	EM045165
	WP Ø 5,0	BHS08513	BHS30 - 5100 (WP)				BHS02007	BHS02008	EM045161
BIONER COMPATIBLE									
IKELT / BIKELT									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
HIKELT									
	NP Ø 3,8	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,7	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOP DM								
	Ø 3,5 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)	 BHS03058  BHS03059		BHS02025	BHS02026	EM048148
MULTIUNIT								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 BHS03008		BHS02003	BHS02004	EM045169
BTI COMPATIBLE								
EH								
	NP Tiny Ø 3,5	BHS08421	BHS30 - 4200 (NP)			BHS02013	BHS02014	EM047957
	RP Ø 4,1	BHS08422	BHS30 - 4200 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,5	BHS08423	BHS30 - 4200 (WP)			BHS02013	BHS02014	EM047967
CUATRILOBULAR								
	RP Ø 4,1	BHS08231	BHS30 - 2300 (RP)			BHS02021	BHS02022	EM047743
	WP Ø 5,5	BHS08232	BHS30 - 2300 (WP)			BHS02021	BHS02022	
BTK COMPATIBLE								
8NECK INT								
	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
EVO DL EXT								
	NP Ø 3,4 / 3,7	BHS08611	BHS30 - 6100 (NP)			BHS02035	BHS02036	
	RP Ø 4,1 / 4,6	BHS08612	BHS30 - 6100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0 / 5,6	BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082
EVO DL INT								
	Ø 3,5 / 3,7 / 4,6 / 5,0 / 5,6	BHS08321	BHS30 - 3200 (NP)	 BHS03050  BHS03051		BHS02021	BHS02022	EM045130

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	⬡ ○	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLASSIC EXT									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	⬡ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KLASSIC INT									
	Ø 3,5/ 4,5 / 5,0	BHS08321	BHS30 - 3200 (NP)	⬡ ○	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
KONIK EXT									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	⬡ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KONIK INT									
	Ø 3,5/ 5,0	BHS08321	BHS30 - 3200 (NP)	⬡ ○	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
SAFETL INT									
	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
B&W COMPATIBLE									
HEXAGONO EXT									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	⬡ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
HEXAGONO INT									
	NP Ø 3,3	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	WP Ø 5,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
CAMLOG COMPATIBLE									
CAMLOG									
	Ø 3,3	BHS08721	BHS30-7200 (Ø3,3)	 	BHS03130 BHS03131		BHS02049	BHS02050	
	Ø 3,8	BHS08722	BHS30-7200 (Ø3,8)	 	BHS03134 BHS03135		BHS02049	BHS02050	
	Ø 4,3	BHS08723	BHS30-7200 (Ø4,3)	 	BHS03138 BHS03139		BHS02049	BHS02050	
	Ø 5,0	BHS08724	BHS30-7200 (Ø5,0)	 	BHS03142 BHS03143		BHS02051	BHS02052	
CONEXÃO SIST. PROTESE COMPATIBLE									
MORSE INDEXADO									
	Ø 3,5 / 3,75 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)	 	BHS03058 BHS03059		BHS02025	BHS02026	EM048148
HEXAGONO EXT									
	Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	Ø 3,75 / 4,0	BHS08742	BHS30 - 7400 (S3 y S4)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEXAGONO INT									
	Ø 3,75 / 4,0 / 5,0	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CORTEX COMPATIBLE									
DYNAMIX									
  	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
CLASSIX									
  	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
SATURN									
  	Ø 3,8 / 4,2	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
COWELLMEDI COMPATIBLE									
EXTERNAL									
  	RP Ø 4,1 WP Ø 5,1	BHS08212 BHS08213	BHS30 - 2600 (RP) BHS30 - 2600 (WP)	 	BHS03003 BHS03004		BHS02001 BHS02039	BHS02002 BHS02040	EM045098 EM045070
INTERNAL									
 	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
SUBMERGED (INT)									
  	Ø 3,5 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103		BHS02009	BHS02010	
DENTALTECH COMPATIBLE									
IMPLAGIC GII									
  	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
IMPLASSIC HX2									
  	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DENTIS COMPATIBLE								
S-CLEAN								
  	RP Ø 3,4	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103	BHS02009	BHS02010	
I-CLEAN								
 	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
DENTUM COMPATIBLE								
SIMPLE LINE II								
 	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
IMPLANTUM								
  	Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103	BHS02009	BHS02010	
SUPER LINE								
  	Ø 3,6 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103	BHS02009	BHS02010	
DENTSPLY SIRONA COMPATIBLE								
XIVE S								
	D 3.4	BHS08221	BHS30-2200 (D3.4)	 	BHS03118 BHS03119	BHS02047	BHS02048	
	D 3.8	BHS08222	BHS30-2200 (D3.8)	 	BHS03122 BHS03123	BHS02047	BHS02048	
	D 4.5	BHS08223	BHS30-2200 (D4.5)	 	BHS03126 BHS03127	BHS02047	BHS02048	

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DIO IMPLANTS COMPATIBLE									
EXTERNAL									
  	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
UF II									
  	Ø 3,8 / 4 / 4,5 / 5 / 5,5	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103		BHS02009	BHS02010	
EDISON MEDICAL COMPATIBLE									
LAMINA (IH)									
  	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
FORTIS (IH)									
  	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
ELITE MEDICA COMPATIBLE									
EXTERNAL CONNECT									
  	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
EUROTEKNIKA COMPATIBLE									
UNEVA									
  	Ø 3,6 / 4,1 / 4,8	BHS08212	BH30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
NATURACTIS									
  	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BH30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NATURALL +									
  	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
NATEA									
  	NP Ø 3,6 / 4,1	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
	RP Ø 4,8 / 6	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
AESTHETICA									
 	Ø 3,6 / 4,1 / 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
GALIMPLANT COMPATIBLE									
IPXE - EXTERNA									
  	RP Ø 3,5 / 4 / 5	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
ICI - INTERNA									
  	RP Ø 3,2 / 3,5 / 4 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
IPX - INTERNA									
  	RP Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
GMI COMPATIBLE									
FRONTIER (IH)									
  	RP Standard Ø 3,3	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	WP Wide Ø 4,3	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
FRONTIER multi-aesthetic									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)		BHS03008		BHS02003	BHS02004	EM045169


1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHOENIX (EH)									
  	Standard Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	Wide (Hex. 2,7) Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	Wide (Hex. 3,4) Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
PHOENIX multi-aesthetic									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 	BHS03008		BHS02003	BHS02004	EM045169
GT MEDICAL COMPATIBLE									
BEST FIT HXE									
  	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
BEST FIT HXI									
 	WP 4,5 / 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
BEST FIT 002									
  	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
BEST FIT OCT									
   	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
HI TECH COMPATIBLE									
SELF THREAD									
  	SLT Ø 3,3 / 3,75	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	GTI Ø 4,2 / 5	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LGI PLUS									
  	Ø 3,5 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
IMPLANT DIRECT COMPATIBLE									
INTERACTIVE									
  	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)	 	BHS03058 BHS03059		BHS02025	BHS02026	EM048148
	RP Ø 4,3	BHS08312	BHS30 - 3100 (RP)	 	BHS03046 BHS03047		BHS02023	BHS02024	EM048149
  	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
REPLUS / REPLANT									
 	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)				BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)				BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
SWISHPLUS / SWISHPLANT									
   	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
JDENTALCARE COMPATIBLE									
JDEVOLUTION									
  	Ø 3,7 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130

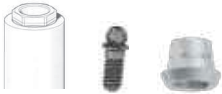
1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KEYSTONE DENTAL COMPATIBLE									
RESTORE									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
OCT-XT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
RENOVA									
	NP Ø 3,75	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5 / 4,75	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
KLOCKNER COMPATIBLE									
SK2 (EH)									
	RP Ø 4,2	BHS08521	BHS30 - 5200 (RP)				BHS02031	BHS02032	EM046747
KL (EH)									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LEADER COMPATIBLE									
TIXOS INTERNAL HEX									
	NP Ø 3,5 / 4,0	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
TIXOS EXTERNAL HEX									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MEDENTIS COMPATIBLE									
ICX									
	Ø 3,45 / 3,7 - 4,1 / 4,8	BHS08751	BHS30-7700 (Ø 3,45 / 3,7 - 4,1 / 4,8)	 	BHS03180 BHS03181		BHS02015	BHS02016	
MEGAGEN COMPATIBLE									
ANYONE INTERNAL									
	NP	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103		BHS02009	BHS02010	
	RP	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103		BHS02009	BHS02010	
	WP	BHS08432	BHS30 - 4600 (4.0)	 	BHS03102 BHS03103		BHS02009	BHS02010	
ANYONE EXTERNAL									
	NP	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082

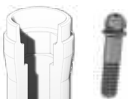1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANYRIDGE									
	Small Ø3,5	BHS08991	BHS30-9900 (4,5/5,0)		BHS03166 BHS03167		BHS02027	BHS02028	
	Regular Ø4,0/4,5	BHS08991	BHS30-9900 (4,5/5,0)		BHS03166 BHS03167		BHS02027	BHS02028	
	Wide Ø5,0/5,5	BHS08991	BHS30-9900 (4,5/5,0)		BHS03166 BHS03167		BHS02027	BHS02028	
OCTA Abutment									
	SYN_RN Ø4,8	BHS08713	BHS30-7100 (RN SYN OCTA Ø4,8)				BHS02013	BHS02014	EM046375
MICRODENT COMPATIBLE									
EKTOS (IH)									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)		BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)		BHS03033 BHS03034		BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
MINI (EH)									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
SUV (EH)									
	SUV Ø 4,1	BHS08612	BHS30 - 6100 (RP)		BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	SUV Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYSTEM (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	


1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
UNIVERSAL (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	SU Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	SU Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MIS COMPATIBLE									
LANCE									
	STD Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SEVEN / M4									
	Narrow Ø 3,30	BHS08951	BHS30-9500 (NARROW)	 	BHS03200 BHS03201		BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	Wide Ø 5,0 / 6,0	BHS08322	BHS30 - 3200 (RP)	 	BHS03033		BHS02021	BHS02022	EM046218
MULTIUNIT MIS									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)		BHS03008		BHS02037	BHS02038	
MOZO GRAU COMPATIBLE									
INHEX									
	RP Ø 3,75 / 4,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
	CONO 20	BHS08124	BHS30 - 1200 (20°)				BHS02003	BHS02004	EM045638
	CONO 45	BHS08125	BHS30 - 1200 (45°)				BHS02003	BHS02004	EM047786

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	○ ○	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OSSEUS									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	○ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
MPI COMPATIBLE									
EXTERNA EH									
	NP Ø 3,5	BHS08421	BHS30 - 4200 (NP)				BHS02001	BHS02002	EM047957
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	○ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INTERNA CM									
	RP Ø 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	○ ○	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
NEOBIOTECH COMPATIBLE									
EB EXTERNAL SYSTEM									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	○ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
IS INTERNAL SYSTEM									
	Ø 3,5 / 4 / 4,5 / 5	BHS08432	BHS30 - 4600 (4.0)	○ ○	BHS03102 BHS03103		BHS02009	BHS02010	

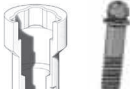
1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NEODENT COMPATIBLE								
TITAMAX HE SMART								
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082
MINI PILAR CM								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 BHS03008		BHS02003	BHS02004	EM045169
NOBEL COMPATIBLE								
BRÄNEMARK								
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)			BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)			BHS02039	BHS02040	EM045070
REPLACE								
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)			BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)			BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)			BHS02009	BHS02010	EM046293
ACTIVE								
	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)	 BHS03058  BHS03059		BHS02025	BHS02026	EM048148
	RP Ø 4,3 / Ø 5,0	BHS08312	BHS30 - 3100 (RP)	 BHS03046  BHS03047		BHS02023	BHS02024	EM048149
MULTIUNIT								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 BHS03008		BHS02003	BHS02004	EM045169

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NT TRADING COMPATIBLE									
2-CONnect									
	M1,4	Ø5,0	BHS08771	BHS30-7700 (M1,4)	 BHS03177		BHS02063	BHS02064	
	M1,6	Ø5,0	BHS08771	BHS30-7700 (M1,6)	 BHS03177		BHS02065	BHS02066	
	M1,8	Ø5,0	BHS08771	BHS30-7700 (M1,8)	 BHS03177		BHS02067	BHS02068	
OSSTEM IMPLANT COMPATIBLE									
OSSTEM TS									
	Ø 4 / 4,5 / 5 / 6		BHS08432	BHS30 - 4600 (4.0)	 BHS03102  BHS03103		BHS02009	BHS02010	
OSSTEM US									
	MINI Ø 3,5		BHS08211	BHS30 - 2100 (NP)			BHS02005	BHS02006	EM045222
	RP Ø 4,1		BHS08212	BHS30 - 2100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,1		BHS08213	BHS30 - 2100 (WP)			BHS02039	BHS02040	EM045070
OSTEOPLUS COMPATIBLE									
ARGENTUM									
	RP Ø 4,2		BHS08521	BHS30 - 5200 (RP)			BHS02031	BHS02032	EM046747
OSTEOPLUS - SHe									
	RP Ø 4,1		BHS08612	BHS30 - 6100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0		BHS08613	BHS30 - 6100 (WP)			BHS02001	BHS02002	EM045082
OSTEOPLUS - SHi									
	RP Ø 3,5		BHS08321	BHS30 - 3200 (NP)	 BHS03050  BHS03051		BHS02021	BHS02022	EM045130
	WP Ø 4,5		BHS08322	BHS30 - 3200 (RP)	 BHS03033  BHS03034		BHS02021	BHS02022	EM046218

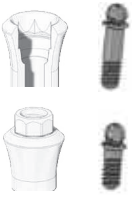1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OXTEIN COMPATIBLE									
HEX EXTERNA									
	MINI Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEX INTERNA									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
OCTA INTERNA									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
PALTOP COMPATIBLE									
INTERNAL HEX CONNECTION									
	Narrow Ø 3,25	BHS08951	BHS30-9500 (NARROW)	 	BHS03200 BHS03201		BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30-3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	Wide Ø 6,0	BHS08322	BHS30-3200 (RP)	 	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
PCA Implant System (Paltop Conical Active)									
	Narrow Ø 3,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
	Standard Ø 3,75 / 4,2 / 5,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)	 	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
MULTI UNIT									
	Ø 5,0	BHS08271	BHS30-2700 (MU)				BHS02031	BHS02032	
SINGLE UNIT									
	Ø 4,5	BHS08251	BHS30-2500 (SINGLE)				BHS02031	BHS02032	

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	⬡ ○	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHIBO COMPATIBLE									
TSH									
	S2 Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	S3 / S4 Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	⬡ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
AUREA EVO									
	NP Ø 3,5			⬡ ○	BHS03106 BHS03107		BHS02041	BHS02042	
	RP Ø 4,3			⬡ ○	BHS03110 BHS03111		BHS02021	BHS02022	
RADHEX IMPLANTS COMPATIBLE									
PHE / PHEA									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)				BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	⬡ ○	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
PHI / PHIA									
	NP Ø3,5	BHS08321	BHS30-3200 (NP)	⬡ ○	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
	RP Ø4,5	BHS08322	BHS30-3200 (RP)	⬡ ○	BHS03033 BHS03034		BHS02021	BHS02022	EM046218
PCI									
	280	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○	BHS03070 BHS03071		BHS02025	BHS02026	EM046371
	350	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
SGS COMPATIBLE									
P1 / P7									
	RP Ø 3,75	BHS08321	BHS30 - 3200 (NP)	⬡ ○	BHS03050 BHS03051		BHS02021	BHS02022	EM045130

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SIGNO VINCES COMPATIBLE									
DUO / INTEGRA									
	RP Ø 4,1	BHS08212	BHS30-2100 (RP)	 BHS03003  BHS03004			BHS02001	BHS02002	EM045098
DUOCON / CM									
	WP_4.5 / 5.0	BHS08123	BHS30-1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
STRAUMANN COMPATIBLE									
OCTA INTERN									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375
BONE LEVEL									
	NC Ø 3,3	BHS08811	BHS30 - 8100 (NC)				BHS02015	BHS02016	EM048048
	RC Ø 4,1/ Ø 4,8	BHS08812	BHS30 - 8100 (RC)				BHS02015	BHS02016	EM046904
BL MULTIBASE									
	RC Ø 4,8						BHS02017	BHS02018	EM047780

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SWEDEN & MARTINA COMPATIBLE									
PREMIUM KOHNO									
	NP Ø 3,3	BHS08961	BHS30 - 9600 (NP Ø 3,3)				BHS02021	BHS02022	
	RP Ø 3,8	BHS08962	BHS30 - 9600 (RP Ø 3,8)				BHS02021	BHS02022	
	WP Ø 4,25	BHS08963	BHS30 - 9600 (WP Ø 4,25)				BHS02011	BHS02012	
SYRA (EH)									
	Ø 4,10 / 4,25	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYRA Short (EH)									
	Ø 4,10	BHS08612	BHS30 - 6100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
TRI DENTAL IMPLANTS COMPATIBLE									
TRI - VENT									
	Ø 3,75 / 4,1 / 4,7	BHS08321	BHS30 - 3200 (NP)	 	BHS03050 BHS03051		BHS02021	BHS02022	EM045130
TRINON COMPATIBLE									
TRINON - Q2									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 	BHS03003 BHS03004		BHS02001	BHS02002	EM045098
TRINON - QK									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)				BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)				BHS02013	BHS02014	EM046375

1.2 CAD/CAM WELDING TO CrCo BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ZIACOM COMPATIBLE								
ZINIC - IH								
	NP Ø 3,2	BHS08951	BHS30-9500 (NARROW)	 BHS03200  BHS03201		BHS02041	BHS02042	
	RP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 BHS03050  BHS03051		BHS02021	BHS02022	EM045130
	WP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 BHS03033  BHS03034		BHS02021	BHS02022	EM046218
OEX - EH								
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)			BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 BHS03003  BHS03004		BHS02001	BHS02002	EM045098
OIN - IH (M2)								
	RP Ø 4,0	BHS08321	BHS30-3200 (NP)	 BHS03050  BHS03051		BHS02011	BHS02012	
	WP Ø 5,0	BHS08322	BHS30-3200 (RP)	 BHS03033  BHS03034		BHS02011	BHS02012	
XDrive (Multi-unit)								
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 BHS03008		BHS02003	BHS02004	EM045169
ZIMMER COMPATIBLE								
INTERNAL HEX								
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 BHS03050  BHS03051		BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 BHS03033  BHS03034		BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)			BHS02021	BHS02022	
SWISS - PLUS								
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS02011	BHS02012	EM045144
	RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNCTA Ø 4,8)			BHS02013	BHS02014	EM046375

1.3 CAD/CAM CEMENTING TO Ti BASE



PROVIDING AESTHETIC AND FUNCTIONAL RESULTS

Ti Base:

U-shape cut for angled screwdriver access up to 30° and provide an optimum vestibular height for a safe cementation.

Indented body design to improve axial retention.

Optimum emergence profile for soft tissue care.

Scanbody:

Wide flat wall for an optimum Best Fit (Position and Orientation).

Manufactured in Titanium Grade 5 with tolerances non exceeding $\pm 5\mu\text{m}$.

Optimal surface finishing for a reliable scanning.

Exact connections for all major Implant Systems.

No need for spray.


1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
AB IMPLANTS COMPATIBLE									
I SERIES INTERNAL HEX									
	Narrow Ø 3 - Ø 3,75	BHS08951	BHS30-9500 (NARROW)	 		BHS03198 BHS03199	BHS02041	BHS02042	
	STD Ø 3,5 - 6,0	BHS08321	BHS30-3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
I 10C CONUS									
	RP Ø 4,2	BHS08312	BHS30 - 3100 (RP)	 		BHS03044 BHS03045	BHS02023	BHS02024	EM048149
ACE COMPATIBLE									
INFINITY EXTERNAL HEX									
	NP Ø 3,3 / 3,75	BHS08611	BHS30 - 6100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,0	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INFINITY TRI - CAM									
	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)	 		BHS03060 BHS03061	BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)	 		BHS03012 BHS03013	BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
INFINITY INTERNAL HEX									
	NP Ø 3,7 / 4,1	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,7 / 5,1	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
ADIN COMPATIBLE									
SWELL									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
TOUAREG S / OS									
	Ø 3,5 - 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOUAREG CLOSEFIT									
	RP Ø 3,5	BHS08311	BHS30 - 3100 (NP)			BHS03056 BHS03057	BHS02025	BHS02026	EM048148
	WP Ø 4,3 / 5	BHS08312	BHS30 - 3100 (RP)			BHS03044 BHS03045	BHS02023	BHS02024	EM048149
ALPHABIO COMPATIBLE									
EXTERNAL ALPHABIO									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)			BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
INTERNAL ALPHABIO									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
ALPHA DENT COMPATIBLE									
ACTIVE									
	Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
ACTIVE PLUS									
	Ø 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 / 3,75 / 4,2 / 5,0	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
ANCLADEN COMPATIBLE									
ANCLALOCK									
	Ø 3,75 - 5	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130

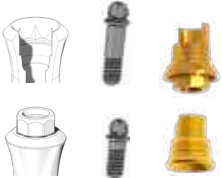1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANTHOGRYR COMPATIBLE									
ANTHOFIT EH									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
OSSFIT									
	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	SYN_RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOC TA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375
ARDS COMPATIBLE									
SMART									
	Ø 3,75 - 4,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
CLASSIC									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
PREMIUM									
	Ø 3,3 - 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
ASTRA COMPATIBLE									
AQUA									
	RP_3.5 / 4.0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
LILAC									
	WP_4.5 / 5.0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
CONE 20°									
	Ø 3,5	BHS08124	BHS30 - 1200 (20°)			BHS03072	BHS02003	BHS02004	EM045638

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CONE 45°									
	Ø 3,5	BHS08125	BHS30 - 1200 (45°)			BHS03074	BHS02003	BHS02004	EM047786
EV									
	Ø3,6	BHS08152	BHS30 - 1500 (Ø3,6)	 		BHS03152 BHS03153	BHS02053	BHS02054	
	Ø4,2	BHS08153	BHS30 - 1500 (Ø4,2)	 		BHS03156 BHS03157	BHS02055	BHS02056	
AVINENT COMPATIBLE									
EH									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02031	BHS02032	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
IH									
	Ø 3,5 / 4,1	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
MULTIUNIT									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02003	BHS02004	EM045169
OCEAN CONICAL (C.C.)									
	Ø 3,5 (Implant Ø 3,5/Ø 4,0)	BHS08311	BHS30 - 3100 (NP)	 		BHS03056 BHS03057	BHS02025	BHS02026	EM048148
	Ø 4,1 (Implant Ø 4,5/Ø 5,0)	BHS08312	BHS30 - 3100 (RP)	 		BHS03044 BHS03045	BHS02023	BHS02024	EM048149
BEGO COMPATIBLE									
EH									
	NP_Ø 2,7 / 3,1	BHS08421	BHS30 - 4200 (NP)				BHS02013	BHS02014	EM047957

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
IH									
  	Ø 3,75 / 4,1	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
BIOGENESIS COMPATIBLE									
3ICON									
  	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
  	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
ITICON									
  	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
  	RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375
ATICON									
  	RP Ø 3,5 / 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	WP Ø 4,5 / 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
KISSES									
  	RP Ø 3,5 / 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	WP Ø 4,5 / 5,0 / 5,5	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
CONE 20°									
  	RP / WP	BHS08124	BHS30 - 1200 (20°)			BHS03072	BHS02003	BHS02004	EM045638
CONE 45°									
  	RP / WP	BHS08125	BHS30 - 1200 (45°)			BHS03074	BHS02003	BHS02004	EM047786

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
BIOHORIZONS COMPATIBLE									
INTERNAL									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)			BHS03020 BHS03021	BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
BIOMET 3i COMPATIBLE									
EH OSSEOTITE									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
CERTAIN									
	NP Ø 3,4	BHS08511	BHS30 - 5100 (NP)				BHS02007	BHS02008	EM045166
	RP Ø 4,1	BHS08512	BHS30 - 5100 (RP)				BHS02007	BHS02008	EM045165
	WP Ø 5,0	BHS08513	BHS30 - 5100 (WP)				BHS02007	BHS02008	EM045161
BIONER COMPATIBLE									
IKELT / BIKELT									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
HIKELT									
	NP Ø 3,8	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,7	BHS08322	BHS30 - 3200 (RP)			BHS03020 BHS03021	BHS02021	BHS02022	EM046218

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
TOP DM									
  	Ø 3,5 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)	 		BHS03056 BHS03057	BHS02025	BHS02026	EM048148
MULTIUNIT									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02003	BHS02004	EM045169
BTI COMPATIBLE									
EH									
  	NP Tiny Ø 3,5	BHS08421	BHS30 - 4200 (NP)				BHS02013	BHS02014	EM047957
	RP Ø 4,1	BHS08422	BHS30 - 4200 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,5	BHS08423	BHS30 - 4200 (WP)				BHS02013	BHS02014	EM047967
CUATRILOBULAR									
 	RP Ø 4,1	BHS08231	BHS30 - 2300 (RP)				BHS02021	BHS02022	EM047743
	WP Ø 5,5	BHS08232	BHS30 - 2300 (WP)				BHS02021	BHS02022	
BTK COMPATIBLE									
8NECK INT									
  	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
EVO DL EXT									
  	NP Ø 3,4 / 3,7	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1 / 4,6	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0 / 5,6	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
EVO DL INT									
  	Ø 3,5 / 3,7 / 4,6 / 5,0 / 5,6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KLASSIC EXT									
  	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KLASSIC INT									
  	Ø 3,5/ 4,5 / 5,0	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
KONIK EXT									
  	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02035	BHS02036	
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
KONIK INT									
  	Ø 3,5/ 5,0	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
SAFETL INT									
  	SR Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
B&W COMPATIBLE									
HEXAGONO EXT									
  	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
HEXAGONO INT									
	NP Ø 3,3	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,0	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	WP Ø 5,0	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
CAMLOG COMPATIBLE									
CAMLOG									
	Ø 3,3	BHS08721	BHS30-7200 (Ø3,3)	 		BHS03128 BHS03129	BHS02049	BHS02050	
	Ø 3,8	BHS08722	BHS30-7200 (Ø3,8)	 		BHS03132 BHS03133	BHS02049	BHS02050	
	Ø 4,3	BHS08723	BHS30-7200 (Ø4,3)	 		BHS03136 BHS03137	BHS02049	BHS02050	
	Ø 5,0	BHS08724	BHS30-7200 (Ø5,0)	 		BHS03140 BHS03141	BHS02051	BHS02052	
CONEXÃO SIST. PROTESE COMPATIBLE									
MORSE INDEXADO									
	Ø 3,5 / 3,75 / 4,0 / 5,0	BHS08311	BHS30 - 3100 (NP)	 		BHS03056 BHS03057	BHS02025	BHS02026	EM048148
HEXAGONO EXT									
	Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	Ø 3,75 / 4,0	BHS08742	BHS30 - 7400 (S3 y S4)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEXAGONO INT									
	Ø 3,75 / 4,0 / 5,0	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
CORTEX COMPATIBLE									
DYNAMIX									
  	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
CLASSIX									
  	Ø 3,3 / 3,8 / 4,2 / 5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
SATURN									
  	Ø 3,8 / 4,2	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
COWELLMEDI COMPATIBLE									
EXTERNAL									
  	RP Ø 4,1	BHS08212	BHS30 - 2600 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2600 (WP)				BHS02039	BHS02040	EM045070
INTERNAL									
  	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
SUBMERGED (INT)									
  	Ø 3,5 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
DENTALTECH COMPATIBLE									
IMPLAGIC GII									
  	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
IMPLASSIC HX2									
  	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DENTIS COMPATIBLE									
S-CLEAN									
  	RP Ø 3,4	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
I-CLEAN									
  	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
DENTIUM COMPATIBLE									
SIMPLE LINE II									
  	RN Ø 4,8	BHS08712	BHS30 - 7600 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
IMPLANTUM									
  	Ø 3,6 / 4 / 4,4 / 4,9 / 6 / 7	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
SUPER LINE									
  	Ø 3,6 / 4 / 4,5 / 5 / 6	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
DENTSPLY SIRONA COMPATIBLE									
XIVE S									
	D 3.4	BHS08221	BHS30-2200 (D3.4)	 		BHS03116 BHS03117	BHS02047	BHS02048	
	D 3.8	BHS08222	BHS30-2200 (D3.8)	 		BHS03120 BHS03121	BHS02047	BHS02048	
	D 4.5	BHS08223	BHS30-2200 (D4.5)	 		BHS03124 BHS03125	BHS02047	BHS02048	

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
DIO IMPLANTS COMPATIBLE									
EXTERNAL									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)			BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
UF II									
	Ø 3,8 / 4 / 4,5 / 5 / 5,5	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
EDISON MEDICAL COMPATIBLE									
LAMINA (IH)									
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
FORTIS (IH)									
	Ø 3,3 / 3,75 / 4,2 / 5,0 / 6,0	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
ELITE MEDICA COMPATIBLE									
EXTERNAL CONNECT									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)			BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
EUROTEKNIKA COMPATIBLE									
UNEVA									
	Ø 3,6 / 4,1 / 4,8	BHS08212	BH30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
NATURACTIS									
	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BH30 - 1200 (RP Ø 3.5 / 4.0)			BHS03068 BHS03069	BHS02025	BHS02026	EM046371

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	⬡ ○	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NATURALL +									
  	Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
NATEA									
  	NP Ø 3,6 / 4,1	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	RP Ø 4,8 / 6	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
AESTHETICA									
  	Ø 3,6 / 4,1 / 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	⬡ ○		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
GALIMPLANT COMPATIBLE									
IPXE - EXTERNA									
  	RP Ø 3,5 / 4 / 5	BHS08212	BHS30 - 2100 (RP)	⬡ ○		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
ICI - INTERNA									
  	RP Ø 3,2 / 3,5 / 4 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
IPX - INTERNA									
  	RP Ø 3,5 / 4 / 4,5 / 5	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	⬡ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
GMI COMPATIBLE									
FRONTIER (IH)									
  	RP Standard Ø 3,3	BHS08321	BHS30 - 3200 (NP)	⬡ ○		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	WP Wide Ø 4,3	BHS08322	BHS30 - 3200 (RP)	⬡ ○		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
FRONTIER multi-aesthetic									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	○		BHS03009	BHS02003	BHS02004	EM045169

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHOENIX (EH)									
	Standard Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	Wide (Hex. 2,7) Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
	Wide (Hex. 3,4) Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
PHOENIX multi-aesthetic									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02003	BHS02004	EM045169
GT MEDICAL COMPATIBLE									
BEST FIT HXE									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
BEST FIT HXI									
	WP 4,5 / 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
BEST FIT 002									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
BEST FIT OCT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCCTA Ø 4,8)				BHS02013	BHS02014	EM046375
						BHS03019			
HI TECH COMPATIBLE									
SELF THREAD									
	SLT Ø 3,3 / 3,75	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	GTI Ø 4,2 / 5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LGI PLUS									
  	Ø 3,5 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
IMPLANT DIRECT COMPATIBLE									
INTERACTIVE									
  	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)	 		BHS03056 BHS03057	BHS02025	BHS02026	EM048148
	RP Ø 4,3	BHS08312	BHS30 - 3100 (RP)	 		BHS03044 BHS03045	BHS02023	BHS02024	EM048149
LEGACY									
  	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
REPLUS / REPLANT									
  	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)	 		BHS03060 BHS03061	BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)	 		BHS03012 BHS03013	BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
SWISHPLUS / SWISHPLANT									
     	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYN OCTA Ø 4,8)	 		BHS03019	BHS02013	BHS02014	EM046375
JDEENTALCARE COMPATIBLE									
JDEVOLUTION									
  	Ø 3,7 / 4,3 / 5 / 6	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
KEYSTONE DENTAL COMPATIBLE									
RESTORE									
	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
OCT-XT									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCITA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375
RENOVA									
	NP Ø 3,75	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,5 / 4,75	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
KLOCKNER COMPATIBLE									
SK2 (EH)									
	RP Ø 4,2	BHS08521	BHS30 - 5200 (RP)	 			BHS02031	BHS02032	EM046747
KL (EH)									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
LEADER COMPATIBLE									
TIXOS INTERNAL HEX									
  	NP Ø 3,5 / 4,0	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
TIXOS EXTERNAL HEX									
  	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MEDENTIS COMPATIBLE									
ICX									
  	Ø 3,45 / 3,7 - 4,1 / 4,8	BHS08751	BHS30-7700 (Ø 3,45 / 3,7 - 4,1 / 4,8)	 		BHS03178 BHS03179	BHS02015	BHS02016	
MEGAGEN COMPATIBLE									
ANYONE INTERNAL									
  	NP	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
	RP	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
	WP	BHS08432	BHS30 - 4600 (4.0)			BHS03100	BHS02009	BHS02010	
ANYONE EXTERNAL									
  	NP	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ANYRIDGE									
	Small Ø3,5	BHS08991	BHS30-9900 (4,5/5,0)			BHS03164 BHS02027 BHS03165		BHS02028	
	Regular Ø4,0/4,5	BHS08991	BHS30-9900 (4,5/5,0)			BHS03164 BHS02027 BHS03165		BHS02028	
	Wide Ø5,0/5,5	BHS08991	BHS30-9900 (4,5/5,0)			BHS03164 BHS02027 BHS03165		BHS02028	
OCTA Abutment									
	SYN_RN Ø4,8	BHS08713	BHS30-7100 (RN SYNOCTA Ø4,8)			BHS03019 BHS02013		BHS02014	EM046375
MICRODENT COMPATIBLE									
EKTOS (IH)									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS02021 BHS03049		BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)			BHS03020 BHS02021 BHS03021		BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)			BHS02021		BHS02022	
MINI (EH)									
	NP Ø 3,5	BHS08611	BHS30 - 6100 (NP)			BHS02031		BHS02032	
SUV (EH)									
	SUV Ø 4,1	BHS08612	BHS30 - 6100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	SUV Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYSTEM (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)			BHS02031		BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)			BHS02031		BHS02032	

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
UNIVERSAL (EH)									
	MST Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	MS Ø 3,5	BHS08611	BHS30 - 6100 (NP)				BHS02031	BHS02032	
	SU Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	SU Ø 5,1	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MIS COMPATIBLE									
LANCE									
	STD Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SEVEN / M4									
	Narrow Ø 3,30	BHS08951	BHS30-9500 (NARROW)	 		BHS03198 BHS03199	BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	Wide Ø 5,0 / 6,0	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
MULTIUNIT MIS									
	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02037	BHS02038	
MOZO GRAU COMPATIBLE									
INHEX									
	RP Ø 3,75 / 4,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
	CONO 20	BHS08124	BHS30 - 1200 (20°)			BHS03072	BHS02003	BHS02004	EM045638
	CONO 45	BHS08125	BHS30 - 1200 (45°)			BHS03074	BHS02003	BHS02004	EM047786

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	○ ○	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OSSEUS									
	NP Ø 3,4	BHS08611	BHS30 - 6100 (NP)				BHS02001	BHS02002	EM045468
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	○ ○		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
MPI COMPATIBLE									
EXTERNA EH									
	NP Ø 3,5	BHS08421	BHS30 - 4200 (NP)				BHS02001	BHS02002	EM047957
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	○ ○		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
INTERNA CM									
	RP Ø 4,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	○ ○		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	WP Ø 5,0	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
NEOBIOTECH COMPATIBLE									
EB EXTERNAL SYSTEM									
	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	○ ○		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	○ ○		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
IS INTERNAL SYSTEM									
	Ø 3,5 / 4 / 4,5 / 5	BHS08432	BHS30 - 4600 (4.0)	○		BHS03100	BHS02009	BHS02010	

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NEODENT COMPATIBLE									
TITAMAX HE SMART									
  	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
MINI PILAR CM									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02003	BHS02004	EM045169
NOBEL COMPATIBLE									
BRÄNEMARK									
  	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
REPLACE									
  	NP Ø 3,5	BHS08111	BHS30 - 1100 (NP)	 		BHS03060 BHS03061	BHS02027	BHS02028	EM046291
	RP Ø 4,3	BHS08112	BHS30 - 1100 (RP)	 		BHS03012 BHS03013	BHS02009	BHS02010	EM046292
	WP Ø 5,0	BHS08113	BHS30 - 1100 (WP)				BHS02009	BHS02010	EM046293
ACTIVE									
  	NP Ø 3,5	BHS08311	BHS30 - 3100 (NP)	 		BHS03056 BHS03057	BHS02025	BHS02026	EM048148
	RP Ø 4,3 / Ø 5,0	BHS08312	BHS30 - 3100 (RP)	 		BHS03044 BHS03045	BHS02023	BHS02024	EM048149
MULTIUNIT									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)			BHS03009	BHS02003	BHS02004	EM045169

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
NT TRADING COMPATIBLE									
2-CONnect									
	M1,4	Ø5,0	BHS08771	BHS30-7700 (M1,4)	○	BHS03176	BHS02063	BHS02064	
	M1,6	Ø5,0	BHS08771	BHS30-7700 (M1,6)	○	BHS03176	BHS02065	BHS02066	
	M1,8	Ø5,0	BHS08771	BHS30-7700 (M1,8)	○	BHS03176	BHS02067	BHS02068	
OSSTEM IMPLANT COMPATIBLE									
OSSTEM TS									
	Ø 4 / 4,5 / 5 / 6		BHS08432	BHS30 - 4600 (4.0)		BHS03100	BHS02009	BHS02010	
OSSTEM US									
	MINI Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,1	BHS08213	BHS30 - 2100 (WP)				BHS02039	BHS02040	EM045070
OSTEOPLUS COMPATIBLE									
ARGENTUM									
	RP Ø 4,2	BHS08521	BHS30 - 5200 (RP)				BHS02031	BHS02032	EM046747
OSTEOPLUS - SHe									
	RP Ø 4,1	BHS08612	BHS30 - 6100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
OSTEOPLUS - SHi									
	RP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	WP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
OXTEIN COMPATIBLE									
HEX EXTERNA									
	MINI Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	WP Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
HEX INTERNA									
	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)			BHS03020 BHS03021	BHS02021	BHS02022	EM046218
OCTA INTERNA									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375
PALTOP COMPATIBLE									
INTERNAL HEX CONNECTION									
	Narrow Ø 3,25	BHS08951	BHS30-9500 (NARROW)			BHS03198 BHS03199	BHS02041	BHS02042	
	Standard Ø 3,75 / 4,2 / 5,0	BHS08321	BHS30-3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	Wide Ø 6,0	BHS08322	BHS30-3200 (RP)			BHS03020 BHS03021	BHS02021	BHS02022	EM046218
PCA Implant System (Paltop Conical Active)									
	Narrow Ø 3,25	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)			BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	Standard Ø 3,75 / 4,2 / 5,0	BHS08122	BHS30 - 1200 (RP Ø 3.5 /4.0)			BHS03068 BHS03069	BHS02025	BHS02026	EM046371
MULTI UNIT									
	Ø 5,0	BHS08271	BHS30-2700 (MU)			BHS03194	BHS02031	BHS02032	
SINGLE UNIT									
	Ø 4,5	BHS08251	BHS30-2500 (SINGLE)			BHS03196	BHS02031	BHS02032	

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
PHIBO COMPATIBLE									
TSH									
  	S2 Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	S3 / S4 Ø 4,1	BHS08742	BHS30 - 7400 (S3 y S4)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
AUREA EVO									
 	NP Ø 3,5						BHS02041	BHS02042	
	RP Ø 4,3						BHS02021	BHS02022	
RADHEX IMPLANTS COMPATIBLE									
PHE / PHEA									
  	NP Ø 3,5	BHS08211	BHS30 - 2100 (NP)	 		BHS03010 BHS03011	BHS02005	BHS02006	EM045222
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
PHI / PHIA									
  	NP Ø3,5	BHS08321	BHS30-3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø4,5	BHS08322	BHS30-3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
PCI									
  	280	BHS08122	BHS30 - 1200 (RP Ø 3.5 / 4.0)	 		BHS03068 BHS03069	BHS02025	BHS02026	EM046371
	350	BHS08123	BHS30 - 1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
SGS COMPATIBLE									
P1 / P7									
  	RP Ø 3,75	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SIGNO VINCES COMPATIBLE									
DUO / INTEGRA									
	RP Ø 4,1	BHS08212	BHS30-2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
DUOCON / CM									
	WP_4.5 / 5.0	BHS08123	BHS30-1200 (WP Ø 4.5 / 5.0)				BHS02009	BHS02010	EM046370
STRAUMANN COMPATIBLE									
OCTA INTERN									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	SYN_RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375
BONE LEVEL									
	NC Ø 3,3	BHS08811	BHS30 - 8100 (NC)	 		BHS03064 BHS03065	BHS02015	BHS02016	EM048048
	RC Ø 4,1/ Ø 4,8	BHS08812	BHS30 - 8100 (RC)	 		BHS03036 BHS03037	BHS02015	BHS02016	EM046904
BL MULTIBASE									
	RC Ø 4,8						BHS02017	BHS02018	EM047780

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM		CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
SWEDEN & MARTINA COMPATIBLE									
PREMIUM KOHNO									
	NP Ø 3,3	BHS08961	BHS30 - 9600 (NP Ø 3,3)				BHS02021	BHS02022	
	RP Ø 3,8	BHS08962	BHS30 - 9600 (RP Ø 3,8)				BHS02021	BHS02022	
	WP Ø 4,25	BHS08963	BHS30 - 9600 (WP Ø 4,25)				BHS02011	BHS02012	
SYRA (EH)									
	Ø 4,10 / 4,25	BHS08612	BHS30 - 6100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
SYRA Short (EH)									
	Ø 4,10	BHS08612	BHS30 - 6100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
	Ø 5,0	BHS08613	BHS30 - 6100 (WP)				BHS02001	BHS02002	EM045082
TRI DENTAL IMPLANTS COMPATIBLE									
TRI - VENT									
	Ø 3,75 / 4,1 / 4,7	BHS08321	BHS30 - 3200 (NP)			BHS03048 BHS03049	BHS02021	BHS02022	EM045130
TRINON COMPATIBLE									
TRINON - Q2									
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)			BHS03005 BHS03006	BHS02001	BHS02002	EM045098
TRINON - QK									
	RN Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)			BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	RN Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOCTA Ø 4,8)			BHS03019	BHS02013	BHS02014	EM046375

1.3 CAD/CAM CEMENTING TO Ti BASE

	PLATFORM Ø	SCANBODY	LIBRARY PLATFORM	 	CrCo BASE	Ti BASE	CLINICAL SCREW	LAB SCREW	IMPLANT REPLICA
ZIACOM COMPATIBLE									
ZINIC - IH									
  	NP Ø 3,2	BHS08951	BHS30-9500 (NARROW)	 		BHS03198 BHS03199	BHS02041	BHS02042	
	RP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	WP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
OEX - EH									
  	NP Ø 3,3	BHS08741	BHS30 - 7400 (S2)				BHS02035	BHS02036	EM048022
	RP Ø 4,1	BHS08212	BHS30 - 2100 (RP)	 		BHS03005 BHS03006	BHS02001	BHS02002	EM045098
OIN - IH (M2)									
  	RP Ø 4,0	BHS08321	BHS30-3200 (NP)	 		BHS03048 BHS03049	BHS02011	BHS02012	
	WP Ø 5,0	BHS08322	BHS30-3200 (RP)	 		BHS03020 BHS03021	BHS02011	BHS02012	
XDrive (Multi-unit)									
  	RP Ø 4,8	BHS08411	BHS30 - 4100 (RP)	 		BHS03009	BHS02003	BHS02004	EM045169
ZIMMER COMPATIBLE									
INTERNAL HEX									
  	NP Ø 3,5	BHS08321	BHS30 - 3200 (NP)	 		BHS03048 BHS03049	BHS02021	BHS02022	EM045130
	RP Ø 4,5	BHS08322	BHS30 - 3200 (RP)	 		BHS03020 BHS03021	BHS02021	BHS02022	EM046218
	WP Ø 5,7	BHS08323	BHS30 - 3200 (WP)				BHS02021	BHS02022	
SWISS - PLUS									
  	RP Ø 4,8	BHS08712	BHS30 - 7100 (RN Ø 4,8)	 		BHS03016 BHS03017	BHS02011	BHS02012	EM045144
	RP Ø 4,8	BHS08713	BHS30 - 7100 (RN SYNOC TA Ø 4,8)	 		BHS03019	BHS02013	BHS02014	EM046375

2.0 BHS30 SCREWDRIVER



EXCELLENT MECHANICAL PROPERTIES & SCREW HEAD PROTECTION

Universal contra angle compatible connection (ISO1797).

Hardened Steel offering excellent torsion resistance (90 N-cm).

Optimal coupling engagement ensuring screw head integrity.

TiN Coated: Surface treatment which increases the life of the driver.

Designed to correct 30° angle deviations with 100% efficient torque transmission.

SCREWDRIVER BHS



BHS01001 - L25mm
BHS01004 - L20mm
BHS01005 - L32mm

HANDLER



BHS01110 - Nob Ø 8mm / CA
BHS01111 - Str Ø 7mm / CA

3.0 PROTOCOLS

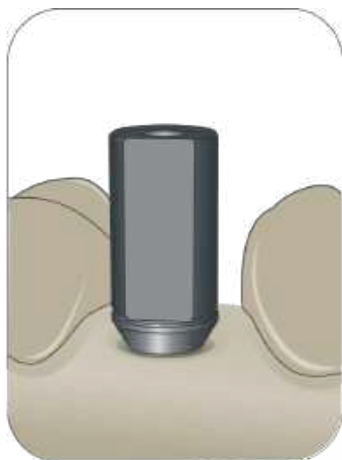
CAD/CAM IMPLANT LEVEL

Products required:


BHS30 Screw
BHS30 Screwdriver
BHS30 Scanbody and associated BHS30 Cad-Cam Library
BHS30 T-Tool

Materials:


Metal Sintering or Milling



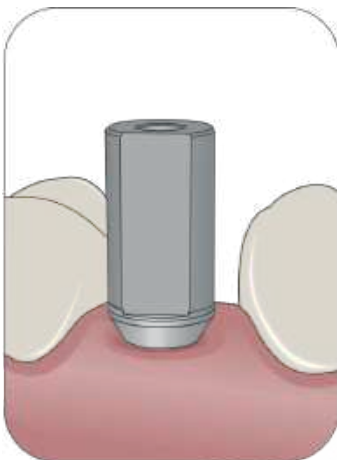1.
CLINIC_ Assembly of the BHS30 IO Scanbody on the patient's implants (flat surface aligned to vestibular) and Intraoral scanning. Delivery of the resulting .stl file to the laboratory.



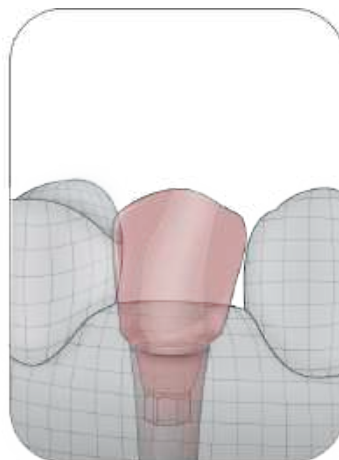
2.
LAB_ Assembly of the BHS EO Scanbody on the implant analogues (flat surface aligned to vestibular) and model scanning in accordance with the software protocols to obtain the .stl file.



3.
LAB_ CAD: Design of the Implant Level structure from .stl file with the corresponding BHS30 Cad Library and selection of the desired Angulation -Free Angulation- for the BHS Screw access channels, according to the undesired emergence on the restoration.



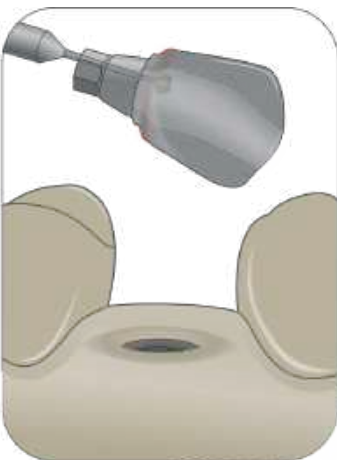4.
LAB_ CAM: Production of the structure as per the software protocols. For a perfect fit and precise finishing of the BHS30 Screw seating platform, the use of the BHS30 T-Tool is recommended. Please refer to our Milling Support service to check on software and T-Tool compatibilities.




5.
LAB_ Adjustment of the structure on the analogues until obtaining passive fit on the model.




6.
LAB_ Adjustment of the structure on the analogues until obtaining passive fit on the model.

3.0 PROTOCOLS

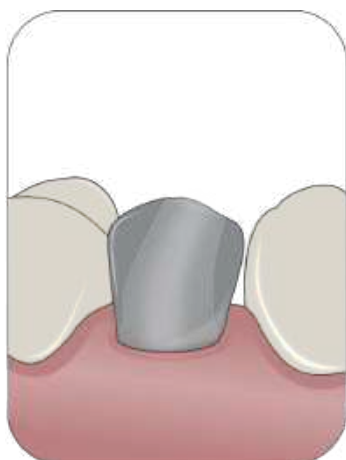CAD/CAM IMPLANT LEVEL



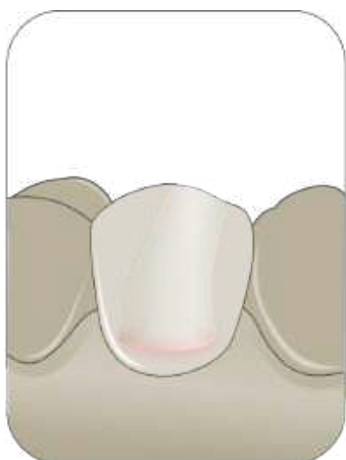7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth.



8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing.



9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS

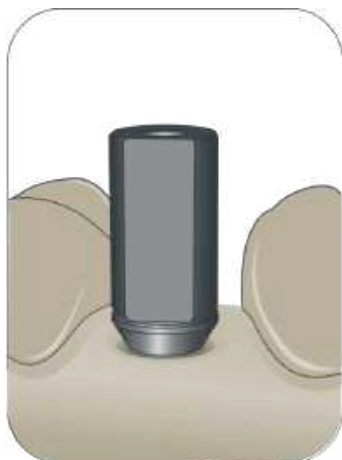
CAD/CAM SINTERING TO CrCo BASE

Products required:


BHS30 Screw
BHS30 Screwdriver
BHS30 Scanbody and associated BHS30 Cad-Cam Library
BHS30 CrCo Base

Materials:


CrCo Sintering or Milling



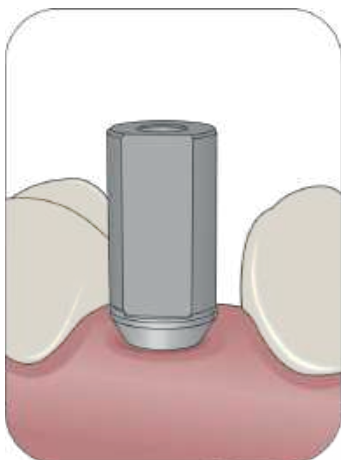1.
CLINIC_ Assembly of the BHS30 IO Scanbody on the patient's implants (flat surface aligned to vestibular) and Intraoral scanning. Delivery of the resulting .stl file to the laboratory.



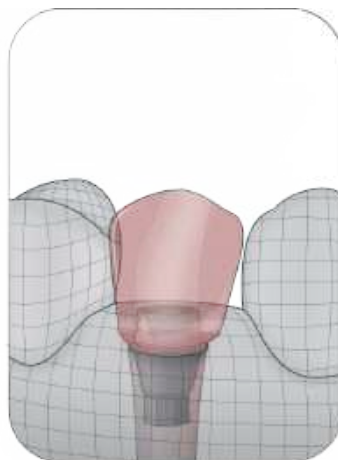
2.
LAB_ Assembly of the BHS EO Scanbody on the implant analogues (flat surface aligned to vestibular) and model scanning in accordance with the software protocols to obtain the .stl file.



3.
LAB_ CAD: Design of the structure onto BHS30 CrCo Base from .stl file with the corresponding BHS30 Cad Library and selection of the desired Angulation -Predefined Angulation- for the BHS30 Screw access channels, according to the undesired emergence on the restoration.



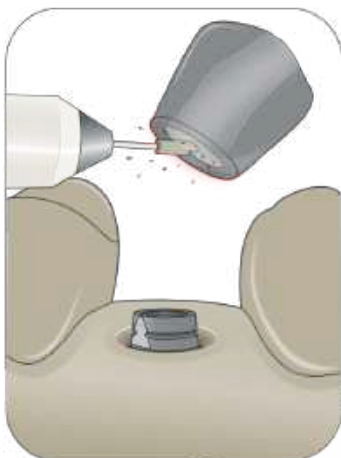4.
LAB_ CAM: Production of the structure onto BHS30 CrCo Base as per the software protocols. Please refer to our Milling Support service to check on software compatibilities.




5.
LAB_ Adjustment of the metal cast structure on the BHS30 CrCo Base until obtaining passive fit on the model and perimetral laser welding of both parts.



6.
LAB_ Adjustment of the metal cast structure on the BHS30 CrCo Base until obtaining passive fit on the model and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS

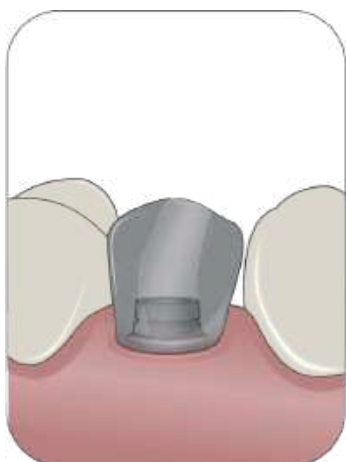CAD/CAM SINTERING TO CrCo BASE



7.
CLINIC_ Metal test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.



8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and perimetral laser welding of both parts. Aesthetic restoration, polishing and finishing.



9.
CLINIC_ Placement of the prosthesis to the patient.

3.0 PROTOCOLS

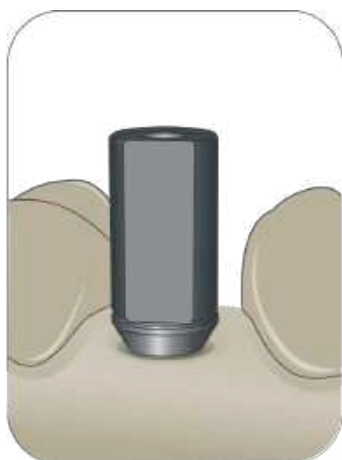
CAD/CAM MILLING TO BHS30 Ti BASE

Products required:


BHS30 Screw
BHS30 Screwdriver
BHS30 Scanbody and associated BHS30 Cad-Cam Library
BHS30 CrCo Base

Materials:


Metal Sintering or Milling



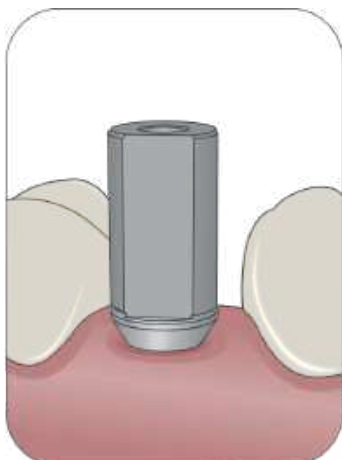1.
CLINIC_ Assembly of the BHS30 IO Scanbody on the patient's implants (flat surface aligned to vestibular) and Intraoral scanning. Delivery of the resulting .stl file to the laboratory.



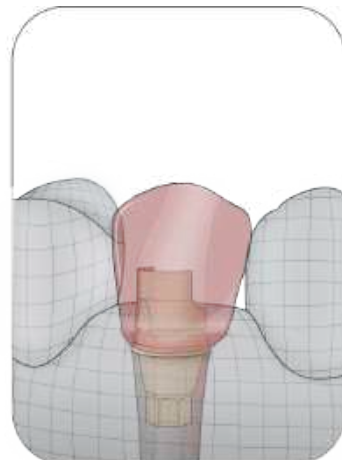
2.
LAB_ Assembly of the BHS EO Scanbody on the implant analogues (flat surface aligned to vestibular) and model scanning in accordance with the software protocols to obtain the .stl file.




3.
LAB_ CAD: Design of the structure onto BHS30 Ti Base from .stl file with the corresponding BHS30 Cad Library and selection of the desired Angulation -Predefined Angulation- for the BHS30 Screw access channels, according to the undesired emergence on the restoration.



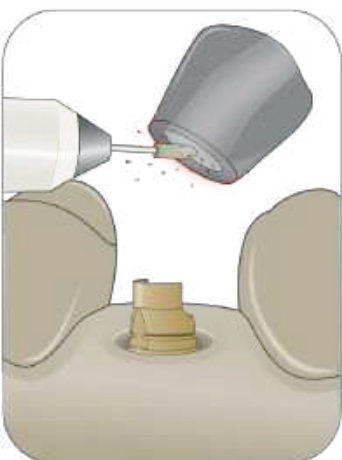4.
LAB_ CAM: Production of the structure onto BHS30 CrCo Base as per the software protocols. Please refer to our Milling Support service to check on software compatibilities.



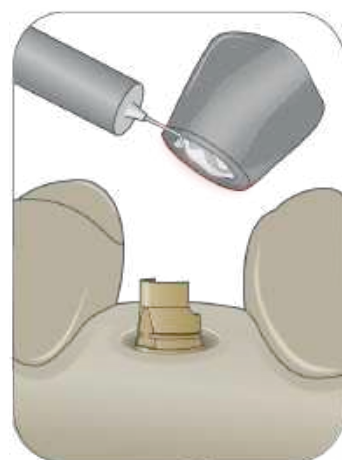
5.
LAB_ Adjustment of the metal cast structure on the BHS30 Ti Base until obtaining passive fit on the model and cementing of both parts.



6.
LAB_ Adjustment of the metal cast structure on the BHS30 CrCo Base until obtaining passive fit on the model and sending the Base and the structure to the Clinic separately for metal testing according to the actual relation in the mouth.

3.0 PROTOCOLS

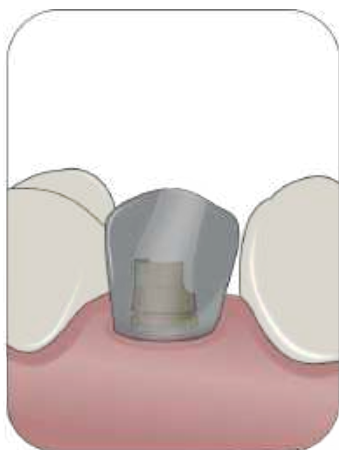CAD/CAM MILLING TO Ti BASE



7.
CLINIC_ Restoration test, especially in case of multiple structures, to confirm adjustment on patient's implants and their positioning in the mouth. Provisional fixation with passive fit in the mouth and delivery back to the lab.



8.
LAB_ Final adjustments in the structure, where appropriate, aesthetic restoration, polishing and finishing. Reassembly of the replicas according to the actual relation in the mouth and cementing of both parts. Aesthetic restoration, polishing and finishing.



9.
CLINIC_ Placement of the prosthesis to the patient.

4.0 COMPATIBILITY CHART

TRADEMARK	PAGES			TRADEMARK	PAGES		
	CAD/CAM IMPLANT LEVEL	CAD/CAM SINTERING CrCo BASE	CAD/CAM MILLING Ti BASE		CAD/CAM IMPLANT LEVEL	CAD/CAM SINTERING CrCo BASE	CAD/CAM MILLING Ti BASE
AB IMPLANTS	6	34	62	IMPLANT DIRECT	20	48	76
ACE	6	34	62	JDENTAL CARE	20	48	76
ADIN	6	34	62	KEYSTONE DENTAL	21	49	77
ALPHABIO	7	35	63	KLOCKNER	21	49	77
ALPHADENT	7	35	63	LEADER	22	50	78
ANCLADEN	7	35	63	MEDENTIS	22	50	78
ANTHOGYR	8	36	64	MEGAGEN	22	50	78
ARDS	8	36	64	MICRODENT	23	51	79
ASTRA	8	36	64	MIS	24	52	80
AVINENT	9	37	65	MOZO GRAU	24	52	80
BEGO	9	37	65	MPI	25	53	81
BIOGENESIS	10	38	66	NEOBIOTECH	25	53	81
BIOHORIZONS	11	39	67	NEODENT	26	54	82
BIOMET 3i	11	39	67	NOBEL	26	54	82
BIONER	11	39	67	NT TRADING	27	55	83
BTI	12	40	68	OSSTEM IMPLANT	27	55	83
BTK	12	40	68	OSTEOPLUS	27	55	83
B&W	13	41	69	OXTEIN	28	56	84
CAMLOG	14	42	70	PALTOP	28	56	84
CON. SISTEMA PROTESE	14	42	70	PHIBO	29	57	85
CORTEX	15	43	71	RADHEX IMPLANTS	29	57	85
COWELLMEDI	15	43	71	SGS	29	57	85
DENTALTECH	15	43	71	SIGNO VINCES	30	58	86
DENTIS	16	44	72	STRAUMANN	30	58	86
DENTIUM	16	44	72	SWEDEN & MARTINA	31	59	87
DENTSPLY SIRONA	16	44	72	TRI DENTAL IMPLANTS	31	59	87
DIO IMPLANTS	17	45	73	TRINON	31	59	87
EDISON MEDICAL	17	45	73	ZIACOM (OSSEOLIFE)	32	60	88
ELITE MEDICA	17	45	73	ZIMMER	32	60	88
EUROTEKNIKA	17	45	73				
GALIMPLANT	18	46	74				
GMI	18	46	74				
GT MEDICAL	19	47	75				
HI TECH	19	47	75				



www.ballheadtechnology.com

(+34) 973 831 391 · info@ballheadtechnology.com